



Breas Vivo 45 LS  
User Manual  
007300 F-2 en-EU

Vivo 45 LS User Manual



Vivo 45 LS, User's  
  
ENGLISH (EU)



Breas Medical AB, Företagsvägen 1, SE-435 33 Mölnlycke, Sweden  
Phone +46 31 86 88 00  
Order +46 31 86 88 20 , Technical Support +46 31 86 88 60  
Fax +46 31 86 88 10,, [www.breas.com](http://www.breas.com), Date: 2023-12-14 | 25425

CE 2797 UK CA 0086

7 321820 073009



# Table of Contents

|       |  |    |
|-------|--|----|
| 1     | Introduction .....                               | 7  |
| 1.1   | What is the Vivo 45 LS? .....                    | 7  |
| 1.2   | Intended Use .....                               | 9  |
| 1.3   | Contraindications .....                          | 9  |
| 1.4   | Intended Environment .....                       | 9  |
| 1.5   | Operation by Lay Users .....                     | 10 |
| 1.6   | About this Manual .....                          | 10 |
| 1.6.1 | Audience .....                                   | 10 |
| 1.6.2 | Icons in this Manual .....                       | 10 |
| 1.7   | Manufacturer Information .....                   | 11 |
| 2     | Safety Information .....                         | 12 |
| 2.1   | General User Precautions .....                   | 12 |
| 2.1.1 | Requirements for Life Supporting Treatment ..... | 14 |
| 2.2   | Electrical Safety .....                          | 15 |
| 2.3   | Environmental Conditions .....                   | 15 |
| 2.4   | Usage of Patient Circuit .....                   | 16 |
| 2.5   | Usage of Filters .....                           | 18 |
| 2.6   | Humidification .....                             | 19 |
| 2.7   | Cleaning and Maintenance .....                   | 20 |
| 2.8   | Usage of Oxygen .....                            | 21 |
| 3     | Product Description .....                        | 23 |
| 3.1   | Main Components .....                            | 23 |
| 3.2   | Front View .....                                 | 25 |
| 3.3   | Side Views .....                                 | 26 |
| 3.3.1 | Detaching and Reattaching the Side Panels .....  | 27 |
| 3.4   | Equipment Designation .....                      | 29 |
| 3.4.1 | Additional Symbols .....                         | 30 |
| 4     | Preparing the Vivo 45 LS for Use .....           | 32 |
| 4.1   | Checking the Vivo 45 LS before First Use .....   | 32 |
| 4.2   | Placing the Vivo 45 LS .....                     | 32 |
| 4.3   | Connecting the Vivo 45 LS to Mains .....         | 34 |
| 4.4   | Connecting the Patient Circuit .....             | 35 |
| 4.5   | Inspecting the Vivo 45 LS before Use .....       | 36 |
| 4.6   | Adjusting the Vivo 45 LS Patient Settings .....  | 37 |
| 4.7   | Performing the Pre-use Test .....                | 37 |
| 4.7.1 | Actions At Pre-Use Test Failure .....            | 38 |
| 5     | How to Use the Vivo 45 LS .....                  | 40 |
| 5.1   | Switch the Vivo 45 LS On and Off .....           | 40 |
| 5.2   | Using the Menu .....                             | 41 |
| 5.2.1 | Navigating with the Buttons .....                | 41 |

|        |   |    |
|--------|---|----|
| 5.2.2  | Symbols Used in the Menu .....                    | 42 |
| 5.2.3  | Menu Overview .....                               | 44 |
| 5.2.4  | Display Overview .....                            | 45 |
| 5.2.5  | The Main Display, Home Mode .....                 | 46 |
| 5.2.6  | The Setup Section, Home Mode .....                | 47 |
| 5.2.7  | The Alarm Section .....                           | 48 |
| 5.2.8  | The Monitor Section .....                         | 49 |
| 5.2.9  | The Others Section .....                          | 49 |
| 5.3    | Monitored Values in the Vivo 45 LS .....          | 50 |
| 5.3.1  | Ppeak .....                                       | 51 |
| 5.3.2  | PEEP .....  | 51 |
| 5.3.3  | Pmean .....                                       | 51 |
| 5.3.4  | Leakage .....                                     | 51 |
| 5.3.5  | MVe .....   | 52 |
| 5.3.6  | MVi .....   | 52 |
| 5.3.7  | Vte .....   | 52 |
| 5.3.8  | Vti .....   | 52 |
| 5.3.9  | FiO <sub>2</sub> .....                            | 52 |
| 5.3.10 | % in TgV .....                                    | 52 |
| 5.3.11 | Total Rate .....                                  | 52 |
| 5.3.12 | Spont Rate .....                                  | 52 |
| 5.3.13 | % Spont .....                                     | 52 |
| 5.3.14 | SpO <sub>2</sub> (Oxygen Saturation) .....        | 53 |
| 5.3.15 | Pulse Rate .....                                  | 53 |
| 5.3.16 | I:E .....   | 53 |
| 5.3.17 | Insp. Time .....                                  | 53 |
| 5.3.18 | Rise Time .....                                   | 53 |
| 5.3.19 | EtCO <sub>2</sub> .....                           | 53 |
| 5.3.20 | InspCO <sub>2</sub> .....                         | 53 |
| 5.3.21 | PtcCO <sub>2</sub> .....                          | 53 |
| 5.4    | Modes in the Vivo 45 LS .....                     | 53 |
| 5.4.1  | Standby and Operating Mode .....                  | 54 |
| 5.4.2  | Device Mode .....                                 | 54 |
| 5.4.3  | Patient Mode .....                                | 54 |
| 5.4.4  | Ventilation and Breath Modes .....                | 54 |
| 5.5    | Functions and Parameters in the Vivo 45 LS .....  | 59 |
| 5.5.1  | Parameters' Availability in Different Modes ..... | 59 |
| 5.5.2  | Insp. Pressure .....                              | 62 |
| 5.5.3  | PEEP .....  | 62 |
| 5.5.4  | Breath Rate .....                                 | 62 |
| 5.5.5  | Backup Rate .....                                 | 62 |

|        |  |    |
|--------|--|----|
| 5.5.6  | SIMV Rate .....  | 62 |
| 5.5.7  | Insp. Time (Inspiratory Time) .....                          | 63 |
| 5.5.8  | Backup Insp. Time (Backup Inspiratory Time) .....            | 63 |
| 5.5.9  | Sigh Parameters .....  | 63 |
| 5.5.10 | Rise Time .....  | 64 |
| 5.5.11 | Insp. Trigger (Inspiratory Trigger) .....                    | 64 |
| 5.5.12 | Sup. Pressure (SIMV mode) .....                              | 65 |
| 5.5.13 | Exp. Trigger (Expiratory Trigger) .....                      | 65 |
| 5.5.14 | Max Insp. Time (Maximum Inspiratory Time) .....              | 65 |
| 5.5.15 | Min Insp. Time (Minimum Inspiratory Time) .....              | 65 |
| 5.5.16 | Target Volume .....  | 66 |
| 5.5.17 | Tidal Volume .....   | 66 |
| 5.5.18 | Flow Pattern .....   | 67 |
| 5.5.19 | CPAP .....   | 67 |
| 5.5.20 | Humidifier .....   | 67 |
| 5.5.21 | Humidifier Setting .....                                     | 67 |
| 5.5.22 | Heated Circuit Temp .....                                    | 67 |
| 5.5.23 | Circuit Heating .....  | 67 |
| 5.6    | Using Batteries .....  | 67 |
| 5.6.1  | Power Source Priority .....                                  | 68 |
| 5.6.2  | Power Switchover Tests .....                                 | 68 |
| 5.6.3  | Charging the Batteries .....                                 | 69 |
| 5.6.4  | Battery Icons .....  | 70 |
| 5.6.5  | Internal Battery .....                                       | 70 |
| 5.6.6  | Click-in Battery .....                                       | 71 |
| 5.6.7  | Battery Operating Time (Internal and Click-in) .....         | 72 |
| 5.6.8  | Storing the Internal Battery and the Click-in Battery .....  | 72 |
| 5.6.9  | External DC .....  | 72 |
| 5.7    | Using Accessories .....                                      | 74 |
| 5.7.1  | Connecting and Disconnecting the Cables .....                | 74 |
| 5.7.2  | Using the ventilator with a Nurse Call System .....          | 74 |
| 5.7.3  | Using the ventilator with the FiO <sub>2</sub> Sensor .....  | 75 |
| 5.7.4  | Using the ventilator with the Remote Alarm .....             | 76 |
| 5.7.5  | Using the ventilator with the EtCO <sub>2</sub> Sensor ..... | 77 |
| 5.7.6  | Using the ventilator with the PtcCO <sub>2</sub> Cable ..... | 81 |
| 5.7.7  | Using the Ventilator with the SpO <sub>2</sub> module .....  | 81 |
| 5.7.8  | Using the ventilator with the Remote Start/Stop .....        | 81 |
| 5.7.9  | Using the ventilator with the Protective Cover .....         | 81 |
| 5.7.10 | Using the Vivo 45 LS with the Trolley .....                  | 82 |
| 5.7.11 | Using the Click-in Humidifier .....                          | 83 |
| 5.7.12 | Using the Patient Circuit with Heated Wire .....             | 88 |

|        |  |     |
|--------|--|-----|
| 6      | Alarms .....   | 90  |
| 6.1    | Alarm Function.....  | 90  |
| 6.1.1  | Alarm Indication .....   | 90  |
| 6.1.2  | Audible Signal Pause .....                                     | 91  |
| 6.1.3  | Audible Signal Presilence.....                                 | 91  |
| 6.1.4  | Alarm Reset .....  | 92  |
| 6.2    | Operator's Position.....                                       | 92  |
| 6.3    | Physiological Alarms .....                                     | 92  |
| 6.3.1  | High Pressure Alarm .....                                      | 93  |
| 6.3.2  | Low Pressure Alarm .....                                       | 94  |
| 6.3.3  | High PEEP Alarm.....   | 95  |
| 6.3.4  | Low PEEP Alarm .....   | 95  |
| 6.3.5  | High $V_t$ (High Inspired Tidal Volume Alarm) .....            | 96  |
| 6.3.6  | Low $V_t$ Alarm (Low Inspired Tidal Volume).....               | 96  |
| 6.3.7  | High $MV_i$ Alarm (High Inspired Minute Volume Alarm).....     | 97  |
| 6.3.8  | Low $MV_i$ (Low Inspired Minute Volume Alarm) .....            | 97  |
| 6.3.9  | High $V_{t_e}$ Alarm (High Expired Tidal Volume) .....         | 98  |
| 6.3.10 | Low $V_{t_e}$ Alarm (Low Expired Tidal Volume).....            | 98  |
| 6.3.11 | High $MV_e$ (High Expired Minute Volume Alarm).....            | 99  |
| 6.3.12 | Low $MV_e$ Alarm (Low Expired Minute Volume) .....             | 99  |
| 6.3.13 | High Breath Rate Alarm.....                                    | 100 |
| 6.3.14 | Low Breath Rate Alarm .....                                    | 100 |
| 6.3.15 | Apnea Alarm.....   | 101 |
| 6.3.16 | Disconnection Alarm .....                                      | 101 |
| 6.3.17 | Rebreathing Alarm .....  | 103 |
| 6.3.18 | Obstruction Alarm .....  | 103 |
| 6.3.19 | High $FiO_2$ Alarm .....                                       | 104 |
| 6.3.20 | Low $FiO_2$ Alarm.....   | 104 |
| 6.3.21 | High $SpO_2$ Alarm .....                                       | 105 |
| 6.3.22 | Low $SpO_2$ Alarm.....   | 105 |
| 6.3.23 | High $EtCO_2$ Alarm .....                                      | 106 |
| 6.3.24 | Low $EtCO_2$ Alarm.....  | 106 |
| 6.3.25 | High $InspCO_2$ Alarm (High Inspired $CO_2$ ).....             | 107 |
| 6.3.26 | High Pulse Rate Alarm.....                                     | 107 |
| 6.3.27 | Low Pulse Rate Alarm .....                                     | 108 |
| 6.3.28 | $PtcCO_2$ Alarm.....   | 108 |
| 6.4    | Technical Alarms .....   | 109 |
| 6.4.1  | Power Fail Alarm .....   | 109 |
| 6.4.2  | High Patient Air Temp. (High Patient Air Temperature) .....    | 109 |
| 6.4.3  | Low Patient Air Temp. (Low Patient Air Temperature Alarm)..... | 109 |
| 6.4.4  | Low Last Power Source Alarm .....                              | 110 |

|        |   |     |
|--------|---|-----|
| 6.4.5  | Crit. Low Last Power Source Alarm .....   | 110 |
| 6.4.6  | Lost Mains Alarm .....  | 110 |
| 6.4.7  | Exhalation Valve Control Error Alarm .....  | 111 |
| 6.4.8  | SpO <sub>2</sub> Disconnected (SpO <sub>2</sub> Sensor Failure/Disconnection Alarm) ..... | 111 |
| 6.4.9  | SpO <sub>2</sub> Signal Lost Alarm .....  | 111 |
| 6.4.10 | Poor SpO <sub>2</sub> Signal .....  | 112 |
| 6.4.11 | CO <sub>2</sub> Disconnected (CO <sub>2</sub> Sensor Failure/Disconnection Alarm) .....   | 112 |
| 6.4.12 | CO <sub>2</sub> Accuracy Error Alarm .....  | 112 |
| 6.4.13 | Check CO <sub>2</sub> Adapter Alarm .....   | 113 |
| 6.4.14 | CO <sub>2</sub> Sensor Error Alarm .....  | 113 |
| 6.4.15 | FiO <sub>2</sub> Disconnected (FiO <sub>2</sub> Sensor Failure/Disconnection Alarm) ..... | 113 |
| 6.4.16 | Ambient Pressure Compensation Lost Alarm .....  | 114 |
| 6.4.17 | Temperature Comp. Lost (Ambient Temperature Compensation Lost Alarm) .....                | 114 |
| 6.4.18 | Humidity Comp. Lost (Humidity Compensation Lost Alarm) .....                              | 114 |
| 6.4.19 | LED Failure Alarm .....   | 115 |
| 6.4.20 | Low Alarm Battery Alarm .....   | 115 |
| 6.4.21 | Alarm Battery Error Alarm .....   | 115 |
| 6.4.22 | Internal/Click-In Battery Hot Alarm .....   | 116 |
| 6.4.23 | Heated Circuit Temp. Alarm .....  | 116 |
| 6.4.24 | High Humidifier Temp. Alarm .....   | 116 |
| 6.4.25 | Humidifier Fault Alarm .....  | 117 |
| 6.4.26 | Heated Circuit Fault Alarm .....  | 117 |
| 6.4.27 | Internal Function Failure .....   | 118 |
| 6.4.28 | Air Temp. Sensor Fail Alarm .....   | 118 |
| 6.4.29 | Internal Error Alarm .....  | 118 |
| 6.4.30 | Database Integrity Fail Alarm .....   | 119 |
| 6.4.31 | Cooling Fan Error Alarm .....   | 119 |
| 6.4.32 | Clock Failure Alarm .....   | 119 |
| 6.4.33 | Internal Temp High Alarm .....  | 119 |
| 6.4.34 | Humidifier/Bypass Loose Alarm .....   | 120 |
| 6.5    | Alarm Test .....  | 120 |
| 6.5.1  | Alarm Signal Test .....   | 120 |
| 6.5.2  | Mandatory Alarm Tests .....   | 120 |
| 6.5.3  | Optional Alarm Tests .....  | 123 |
| 7      | Cleaning and Maintenance .....  | 124 |
| 7.1    | Cleaning the Vivo 45 LS .....   | 124 |
| 7.1.1  | Main Unit .....   | 124 |
| 7.1.2  | Air Pathway Disinfection .....  | 125 |

|     |       |   |     |
|-----|-------|---|-----|
|     | 7.1.3 | Patient Circuit.....  | 125 |
| 7.2 |       | Cleaning and Replacing the Filters.....   | 126 |
|     | 7.2.1 | Washing a coarse filter.....  | 127 |
| 7.3 |       | Change of Patients.....   | 127 |
| 7.4 |       | Regular Maintenance.....  | 127 |
| 7.5 |       | Service and Repair.....   | 127 |
| 7.6 |       | Storage.....  | 128 |
| 7.7 |       | Disposal.....   | 128 |
| 8   |       | Technical Specifications.....   | 129 |
| 8.1 |       | System Description.....   | 129 |
|     | 8.1.1 | Pneumatic Diagram for the ventilator.....   | 130 |
| 8.2 |       | Data.....   | 130 |
|     | 8.2.1 | Worst Case Accuracy.....  | 130 |
|     | 8.2.2 | Modes Specifications.....   | 130 |
|     | 8.2.3 | Parameter Specifications.....   | 131 |
|     | 8.2.4 | Monitored Values Specifications.....  | 135 |
|     | 8.2.5 | Power Supply.....   | 136 |
|     | 8.2.6 | Environmental Conditions.....   | 136 |
|     | 8.2.7 | Other.....  | 137 |
| 8.3 |       | Emission and Immunity Declaration.....  | 138 |
|     | 8.3.1 | Vivo 45 LS Essential Performance.....   | 138 |
|     | 8.3.2 | Guidance and Manufacturer's Declaration – Electromagnetic Immunity.....   | 139 |
|     | 8.3.3 | Guidance and Manufacturer's Declaration – Electromagnetic Emission.....   | 140 |
|     | 8.3.4 | Frequencies of portable and mobile transmitters for which the recommended separation distance is 30 cm (12 inches)..... | 141 |
|     | 8.3.5 | Recommended separation distances between external power conductors and the ventilator.....                              | 141 |
|     | 8.3.6 | Proximity fields from RF wireless communications equipment.....   | 142 |
| 8.4 |       | Delivery Settings.....  | 142 |
| 9   |       | Accessories and Parts.....  | 144 |
|     | 9.1   | Patient Circuits and Air Delivery Accessories.....  | 144 |
|     | 9.2   | Power Accessories.....  | 148 |
|     | 9.3   | Monitoring Accessories.....   | 150 |
|     | 9.4   | Ventilator Filters and Detachable Parts.....  | 154 |
|     | 9.5   | Other Accessories.....  | 155 |
| 10  |       | Patient Settings.....   | 159 |
| 11  |       | FAA Compliance.....   | 160 |
|     |       | Index.....  | 161 |



# Introduction

## WARNING!

### Risk of Personal Injury

The Vivo 45 LS must only be used:

- For the intended treatment in accordance with this manual and with the instructions given by the responsible clinical personnel.
- In accordance with the operating conditions specified in this manual.
- In original and unmodified shape and only with compatible accessories.

Every other use may lead to risk of personal injury!



## CAUTION!

Read this manual thoroughly so that you completely understand how the Vivo 45 LS is operated and maintained before taking it into use, to ensure correct usage, maximum performance and serviceability. Non-professional caregivers (e.g. family members) should consult the medical equipment provider's respiratory therapist if they have any questions about the function, proper use, operation, service or maintenance of the Vivo 45 LS.

Breas Medical reserves the right to make changes to this product without any prior notification.

U.S. Federal law restricts this device to sale by or on order of a physician.



## 1.1

## What is the Vivo 45 LS?

The Vivo 45 LS is a pressure and volume ventilator capable of delivering continuous or intermittent ventilatory support for patients who require invasive or non-invasive mechanical ventilation. The Vivo 45 LS is capable of running 24 hours/day.

The Vivo 45 LS can be operated in the following ventilation modes:

- PSV (Pressure Support Ventilation). May be combined with:
  - TgV (Target Volume)
  - AE (Auto-EPAP)
- PCV (Pressure Controlled Ventilation). May be combined with:
  - A (Assist function — inspiration trigger)
  - TgV (Target Volume)
  - AE (Auto-EPAP)
  - MPV (Mouthpiece Ventilation)
  - SIMV (Synchronized Intermittent Mandatory Ventilation)
- VCV (Volume Controlled Ventilation). May be combined with:
  - A (Assist function — inspiration trigger)
  - MPV (Mouthpiece Ventilation)
  - SIMV (Synchronized Intermittent Mandatory Ventilation)
- CPAP (Continuous Positive Airway Pressure)

## Compatible Patient Circuits

The Vivo 45 LS can be used with a leakage circuit, an MPV circuit or a circuit with active exhalation valve. See 9 *Accessories and Parts*, page 144 for detailed information about compatible patient circuits.

For leakage circuits: The patient circuit shall comply to ISO 17510. The leakage should be at least 12 l/min at 4 cmH<sub>2</sub>O, to prevent rebreathing of exhaled air. The recommended leakage is 20 to 50 l/min at 10 cmH<sub>2</sub>O pressure.

## Compatible Patient Interfaces

For invasive use, the patient interface may be a tracheostomy tube (cuffed or uncuffed).

For non-invasive use it may be a mask, mouthpiece or pillow interface. See the patient interface's instructions for use when selecting the interface to use.

## Data Log

The Vivo 45 LS has an internal memory with a data log that holds the following data:

- Running hours
- Technical alarms
- Settings
- Asset data
- Treatment hours
- Treatment settings
- Device serial number
- Physiological alarms
- Detailed log, containing at least 24 h data of clinical data (monitored values)
- Breath log, containing at least 30 day data of (monitored values)
- Usage log (containing at least 1 year data of non-clinical events, alarms and settings)

The log data is maintained when device is powered off or in case of power failure. When logs are filled up, the oldest data will be discarded. The data can be transferred to a computer, printed out, and analyzed via Breas software products, see .



For more information about Breas software products, please contact your Breas representative.

## Multiple Use

This is a multiple patient multiple use ventilator. If it should be used by multiple patients, see the cleaning instructions in 7.3 *Change of Patients*, page 127 before assigning it to a new patient. Note that accessories to the ventilator might be for single patient use and should in that case be replaced at change of patient.

## Expected Service Life

The expected service life of the Vivo 45 LS is 8 years.

## 1.2 Intended Use

The Vivo 45 LS ventilator (with or without the SpO<sub>2</sub> and CO<sub>2</sub> sensor) is intended to provide continuous or intermittent ventilatory support for the care of individuals who require mechanical ventilation. Specifically, the ventilator is applicable for paediatric through adult patients weighing more than 5 kg (11 lbs.)

The Vivo 45 LS with the SpO<sub>2</sub> is intended to measure functional oxygen saturation of arterial hemoglobin (%SpO<sub>2</sub>) and pulse rate.

The Vivo 45 LS with the CO<sub>2</sub> sensor is intended to measure CO<sub>2</sub> in the inspiratory and expiratory gas.

The device is intended to be used in home, institution, hospitals and portable applications such as wheelchairs and gurneys. It may be used for both invasive and non-invasive ventilation. The Vivo 45 LS is not intended to be used as an emergency transport or critical care ventilator.

## 1.3 Contraindications

- The use of the Vivo 45 LS is contraindicated for patients who need to be ventilated with oxygen concentrations (FiO<sub>2</sub>) higher than achievable when combining inlet from a low pressure oxygen source at 30 l/min with actual ventilator settings.
- Generally, after surgery, the surgeon should be consulted to avoid organ damage and help determine ventilator parameters that do not adversely affect hemodynamics or have a negative impact on the patient's health status.
- In case of facial surgery, make sure to choose a suitable patient interface, in order to avoid discomfort and injury.
- The Vivo 45 LS is not intended to be used as an emergency transport ventilator or critical care ventilator.

## Undesirable Side Effects

If the patient experiences chest discomfort, pain, severe headache or shortness of breath while using the Vivo 45 LS, a physician or responsible clinician should be contacted immediately.

The following side effects may occur during the course of therapy with the Vivo 45 LS, patients are advised to report any new or changing adverse effects to their physician:

- Nasal, mouth or throat dryness
- Nosebleeds
- Abdominal bloating
- Ear or sinus discomfort
- Eye irritation
- Skin rashes

## 1.4 Intended Environment

The Vivo 45 LS is intended to be used in clinical settings (e.g. hospitals, sub-acute care institutions), public spaces and home environments as well as during portable applications such as wheelchairs, personal family vehicles, ground ambulances and civil aircraft, excluding helicopters.

It is not intended for use during emergency transportation.

# 1.5 Operation by Lay Users

Day-to-day caregivers, patients, relatives and other non-professional users may operate the Vivo 45 LS with the *Home mode* activated, after it has been set up according to the prescribed treatment. In Home mode, some settings and controls are locked or hidden.

The lay operator may also perform basic maintenance that doesn't require special equipment or a service environment.

## Training

The lay operator shall be trained to basic knowledge of the Vivo 45 LS and in the specific operations they are assigned to perform. The training shall be based on this user manual and the responsible clinical personnel shall assess the level of training required for each lay operator. This manual shall be available for training and as reference when operating the Vivo 45 LS.

# 1.6 About this Manual

## CAUTION!

Always read this manual before setting up and using the Vivo 45 LS or performing maintenance on the machine, to ensure correct usage, maximum performance and serviceability.






## 1.6.1 Audience





This manual is intended for patients and other lay users operating the Vivo 45 LS.

- Care providers, clinical personnel, physicians and others who require a working knowledge of the Vivo 45 LS will find additional information on settings and functions in the Clinician's Manual. The Clinician's manual should be of the same revision as the User's Manual.
- Service personnel may order the Service Manual that contains detailed technical information for maintenance, service, repair and disposal procedure. The Service manual's revision is independent of the User's Manual revision.

## 1.6.2 Icons in this Manual

In this manual, icons are used to highlight specific information. The meaning of each icon is explained in the table below.

| Icon  | Explanation  |
|---|--|
|  | <b>Warning!</b><br>Risk of death or personal injury. |
|  | <b>Warning!</b><br>Risk of Cross-contamination.      |
|  | <b>Warning!</b><br>Risk of electric shock.           |
|  | <b>Warning!</b><br>Hot surface, risk of burns.       |
|  | <b>Warning!</b><br>Flammable material, risk of fire. |

| Icon  | Explanation  |
|---|--|
|  | <b>Caution!</b><br>Risk of equipment damage, loss of data, extra work, or unexpected results.              |
|  | MR Unsafe. The device should not enter a magnetic resonance (MR) environment, such as an MRI scanner room. |
|  | <b>Note</b><br>Information that may be valuable but is not of critical importance, tips.                   |
|  | <b>Reference</b><br>Reference to other manuals with additional information on a specific topic.            |

## 1.7

### Manufacturer Information

#### Legal Manufacturer



#### Postal Address

Breas Medical AB  
Företagsvägen 1  
SE-435 33 Mölnlycke Sweden

#### Web Address

[www.breas.com](http://www.breas.com)

#### Email address

[breas@breas.com](mailto:breas@breas.com)

#### Phone

+46 31 868800, Order: +46 31 868820, Technical support: +46 31 868860

#### Fax

+46 31 868810

#### UK Sales and Support Information

**UK Responsible Person:** Breas Medical Ltd

#### Phone UK Head Office

+44 (0)1789293460

#### Fax UK Head Office

+44 (0)1789262470

#### General Enquiries and Ordering

[orders@nippyventilator.com](mailto:orders@nippyventilator.com)

#### Tech Support

[techsupport@nippyventilator.com](mailto:techsupport@nippyventilator.com)

#### Local Representatives Information

[www.breas.com/contact-us/](http://www.breas.com/contact-us/)

## General User Precautions

**WARNING!****Risk of Personal Injury**

The Vivo 45 LS must only be used:

- For the intended treatment in accordance with this manual and with the instructions given by the responsible clinical personnel.
- In accordance with the operating conditions specified in this manual.
- In original and unmodified shape and only with compatible accessories.

Every other use may lead to risk of personal injury!

**Risk of Insufficient Ventilation**

Usage outside the specified operating conditions may cause reduced performance.

The Vivo 45 LS must only be used in accordance with the operating conditions specified in this manual.

**Risk of Reduced Safety and Performance**

Accessories that have not been verified to be compatible with the Vivo 45 LS might affect safety features and performance negatively.

Only use the Vivo with accessories that are compatible with the ventilator. Use of incompatible parts to connect the ventilator to the patient can result in degraded performance and change of pressure gradient.

Breas Medical has verified the compatibility between the Vivo 45 LS and the accessories listed in the Clinician's manual.

The responsible organization must ensure the compatibility of the ventilator with all parts used to connect to the patient before starting the intended treatment. If incompatible accessories are used, Breas Medical has no responsibility for the safe and effective use of the Vivo 45 LS.

Changes to the patient circuit, like adding or removing accessories or changing type or length of breathing tube, may affect both circuit compliance and alarm triggering conditions.

It's recommended to perform a pre-use test and re-test the alarm function after making changes to the patient circuit.

When a patient is treated, there must be a supervising person present during the treatment in order to take care of alarms and conditions that the patient cannot solve on their own.

Always have immediate access to an alternative means of ventilation, which is ready for use, to avoid patient death or serious injury. Failure to have an alternate means of ventilation can result in serious injury or patient death if ventilator fails.

To prevent death or serious injury, monitor the patient and the ventilator regularly in order to determine the need to provide emergency ventilation when an alarm sounds or the ventilator malfunctions.

If using the ventilator in a transport case, only use the protective cover specified in the accessory list to prevent adverse ventilator performance, which can consequently result in patient death.

Measured values for volume and expired CO<sub>2</sub> may differ from the actual patient values due to unintentional leaks.

Do not obstruct the gas intake port.





To prevent disconnection of the patient circuit during use, only tubes in compliance with ISO 5367 or ISO 80601-2-74 should be used.



The ventilator may not work properly if any part has been dropped, damaged or submerged in water.

## **WARNING!**



### **Risk of Burns**

Covering breathing tubes with a blanket or heating them with an overhead heater can affect the quality of the therapy or injure the patient.

## **WARNING!**



### **Risk of Faulty Treatment**

If the patient is admitted to a hospital or is prescribed any other form of medical treatment, always inform the medical staff that the patient is on mechanical ventilation treatment.

### **Risk of Faulty Treatment**

Do not use the Vivo 45 LS in the event of:

- Suspected damage to the device, including the occurrence of Internal Functional Failure alarms.
- Unexpected patient symptoms during treatment.
- Unexplainable or sudden changes of pressure, performance or sound during operation.
- Delivered air being abnormally hot or emitting an odor.

Contact your responsible care provider for an inspection.



### **Risk of Faulty Treatment**

The responsible organization should periodically reassess the settings of the therapy for effectiveness.



The ventilator therapy settings must always be based on medical advice and must be carried out by authorised clinical personnel only. When changing treatment settings or changing to another device a clinical assessment must be performed to determine if blood gas measurement needs to be performed.



Before starting treatment, always perform the procedure 4.5 *Inspecting the Vivo 45 LS before Use*, page 36.

### **Risk of Unnoticed Critical Conditions**

- All the physiological alarms of the Vivo 45 LS must be set at safe levels that will effectively warn the user of any risk.  
The alarm levels should be assessed considering the patient's treatment settings.
- Any change of treatment settings or change of components in the ventilation system may require readjustment of the alarm levels.
- The alarm sound level should be set to a clearly audible level. Setting the alarm sound level below that of the ambient sound level can impede recognition of alarm conditions.



## **CAUTION!**

Clinical personnel must read the Clinician's manual thoroughly and understand the ventilator operation before setting up and using the ventilator.



- Handle the ventilator with care.
- Do not use the ventilator while in the carry bag.
- Do not use the ventilator with nitric oxide, helium or helium mixtures.



Contact Injuries: Skin irritation may occur due to prolonged exposure to either a mask (if used) or the SpO<sub>2</sub> module.



Ensure that the cooling air intake vents are not blocked. If the vents are blocked, especially in hot use environments, the surface temperature of the patient circuit may rise above 41°C (106°F). In a 40°C (104°F) environment and with the cooling air intake vents blocked, surface temperatures as high as 50°C (122°F) can occur. Before an unsafe temperature is reached, the “High Patient Air Temp” alarm will occur. If this alarm occurs, assure that the ventilator cooling air intake path is free of obstruction and that the patient circuit surface is not heating the patient’s skin.

## **WARNING!**

### **Risk of Electric Shock**

Modifying or using the ventilator with accessories that are not compatible may cause cardiac arrhythmia.

The Vivo 45 LS must only be used in original and unmodified shape and only with compatible accessories.

Inadequate use of device or accessories may cause loss of treatment or decreased performance.



## **CAUTION!**

If you suspect that the device has been mistreated, perform a functional check before taking it to use. A basic functional check can be performed as described in 4.5 *Inspecting the Vivo 45 LS before Use*, page 36. A complete functional check can be performed by an authorized service technician.



## **NOTE**

Any serious incident that has occurred in relation to this device should be reported to the competent authority and the manufacturer.

## **2.1.1 Requirements for Life Supporting Treatment**

Life supporting- treatment requires that:

- An emergency equipment (e.g. resuscitation bag) is available.
- One of the following means of surveillance is used:
  - Using the EtCO<sub>2</sub> sensor accessory, or an external EtCO<sub>2</sub> monitor (capnometer). This surveillance method can be used for ventilation with active exhalation valve circuits as well as leakage port circuits.  
The CO<sub>2</sub> sensor must be connected between the patient and the exhalation valve or leakage port to be able to measure exhaled gases. If using an external CO<sub>2</sub> monitor, it shall fulfil the ISO 80601-2-55 standard (Medical electrical equipment – Particular requirements for the basic safety and essential performance of respiratory gas monitors).
  - Supervising the ventilator’s monitoring of exhaled volume. This surveillance method can be used for ventilation with leakage circuits only.



## 2.2

## Electrical Safety

### WARNING!

#### Risk of Electric Shock

High voltage contact may cause cardiac arrhythmia.

- Do not operate the Vivo 45 LS if it has a damaged power cord, power supply or casing.
- To avoid electrical shock, only clean the Vivo 45 LS according to instructions in this manual. Do not soak or immerse the Vivo 45 LS into any fluids.
- Use the approved power supply units only.  
Use of unapproved power supply units may compromise the electrical isolation and lead to risk of electric shock.
- Do not use more than one multiple portable socket-outlet or extension cord.  
If a multiple portable socket-outlet is used, it must not be placed on the floor.
- The operator must not touch accessible contacts of connectors and the patient simultaneously.
- Nurse Call must only be connected to a safety extra low voltage system with an isolation from AC power (Mains) voltage which complies with the requirements of IEC 60601-1.



### WARNING!

#### Risk of Faulty Treatment

Electromagnetic Interference may cause electrical equipment to malfunction.

- The aspects of electromagnetic compatibility must be considered.
  - The Vivo 45 LS should not be used adjacent to or stacked with other equipment; if adjacent or stacked use is necessary, the Vivo 45 LS should be observed to verify normal operation in that configuration.
  - Mobile or transportable radio transmitters may interfere with the Vivo 45 LS.
  - Further guidance for safe installation of the ventilator can be found in the chapter about emission and immunity declaration.
- If a portable AC power supply is used, make sure that the voltage variations are within the operating limits of the Vivo 45 LS.
- Portable RF communications equipment (including peripherals such as antenna cables and external antennas) should be used no closer than 30 cm (12 inches) to any part of the Vivo 45 LS, including cables specified. Otherwise, degradation of the performance of this equipment could result.



### WARNING!

Avoid touching the contacts within the ventilator click-in battery compartment. Under certain circumstances touch current limits per IEC 60601-1 may be exceeded.



## 2.3

## Environmental Conditions

### WARNING!

#### Risk of Intoxication

Do not use the Vivo 45 LS in a toxic environment.



## **WARNING!**



### **Risk of Fire**

Do not use the Vivo 45 LS in environments where explosive gases or flammable anesthetic agents present.

## **WARNING!**



The delivered patient air can be as much as 4°C (7°F) higher than ambient temperature. Caution should be exercised if the room temperature is greater than 36°C (97°F).



### **Risk of Faulty Treatment**

If a room humidifier is used, place it at least 2 meters away from the Vivo 45 LS.

### **Risk of faulty Treatment**


The performance of the Vivo 45 LS may deteriorate at altitudes or ambient temperatures outside the operation conditions specified in the chapter *Technical Specifications*.



- Do not use the ventilator while positioned in a warm place, such as direct sunlight or close to a radiator as this might lead to temperature outside the specifications.
- Do not use the ventilator in an hyperbaric chamber, as this would cause an ambient pressure outside the specifications.
- Do not use the ventilator immediately after storage or transport outside the recommended operating conditions.

### **Risk of Faulty Treatment**



 Do not use or store the Vivo 45 LS in a magnetic resonance (MR) environment. Use of the Vivo 45 LS in an MR environment may result in malfunction of the Vivo 45 LS and pose unacceptable risk to the patient, medical staff or other persons.



Unsteady indicated values for delivered volumes or pressures and the occurrence of alarm conditions without apparent cause may be an indication of loss of performance due to electromagnetic disturbances. Follow the instructions above and the guidance provided in 8.3 *Emission and Immunity Declaration*, page 138 to mitigate the effects of electromagnetic disturbances.

## **CAUTION!**



The ventilator, any accessories and all replaced parts, must be disposed of in accordance with the local environmental regulations regarding the disposal of used equipment and waste.

## **2.4 Usage of Patient Circuit**

### **WARNING!**



The ventilator supports leakage circuits, circuits with an active exhalation valve and circuits with mouthpiece interface. Failure to use a mask or accessory that minimizes rebreathing of carbon dioxide or permits spontaneous breathing can cause asphyxiation.



For the ventilator to deliver treatment according to settings, it is important that the selection of the patient circuit type is correctly set.

### **Risk of Insufficient Ventilation**

Insufficient ventilation may cause transient hypoxia.

Before use:

- Make sure that the patient circuit and joined parts are undamaged and correctly connected, in order to avoid unwanted leakage.
- For leakage circuits: The leakage port of the patient circuit or patient interface prevents rebreathing by flushing the exhaled air. It should be located as near the patient interface as possible (this is even more important for treatments with low pressure). Make sure that it is not blocked or obstructed.
- For active exhalation valve circuits: Check the function of the exhalation valve and that it is not blocked or obstructed.
- The Vivo 45 LS should be turned on and the function of the leakage port should be checked before use:

### **Risk of Suffocation**

If the patient needs assistance to remove the patient interface, the patient shall not be left alone. This is to avoid the risk of re-breathing of CO<sub>2</sub> in case of accidental ventilator failure.

Do not breathe through the connected patient circuit unless the ventilator is turned on and operating properly.

### **WARNING!**

#### **Risk of Electric Shock**

Do not use antistatic or electrically conductive hoses or tubing with the ventilator breathing system. This could result in electrical shock.

### **WARNING!**

Patient connected parts and all filters must be replaced regularly to ensure correct function of the ventilator. All replaced parts must be disposed of according to local environmental regulations regarding the disposal of used equipment and parts.

Patient circuits used with the ventilator shall have the following characteristics:

- Length: Max 6 feet (2 meters).
- Connector: 22 mm.
- Resistance: Max 2 mbar at 40 l/min.

If an active exhalation valve is used:

- Inner diameter of the exhalation valve control pressure tube: 3 mm.
- The exhalation valve must be of the type that is open (letting the exhaled patient air out) when unpressurised by the control pressure

By conducting a pre-use test (see 4.7 *Performing the Pre-use Test*, page 37) the compatibility of the complete patient circuit configuration with the ventilator can be verified. If a pre-use test is successfully performed the circuit configuration meets the required characteristics.

Always perform a pre-use test when the patient circuit is replaced or modified. Also check the alarm settings as changes to the patient circuit may affect the alarm triggering.

### **Risk of Suffocation**

Periodically check for moisture in the patient circuit.

When present, remove the moisture. Before attempting to dry the circuit, disconnect it from the Vivo 45 LS to ensure no water flows back into the Vivo 45 LS.

The frequency at which these checks must be performed will depend on the patient's condition and the device used. The responsible caregiver should assess this on an individual basis in accordance with the patient's needs.

When using the ventilator invasively, the low volume alarm and the low breath rate alarm must be carefully set, to ensure safe use.

### **Risk of Insufficient Ventilation**

Insufficient ventilation may cause transient hypoxia.

The use of equipment such as endotracheal tubes, oral/nasal tubes, adaptors etc. with small inner diameters or high resistance filters, humidifiers etc. increase the resistance in the patient circuit which may interfere with the operation of the patient disconnect function. It may also interfere with the device trigger function.

This impact can be reduced by conducting a pre-use test.

Make sure that the exhalation valve or the leakage port never is blocked or obstructed.

### **Risk of Constriction**

Entanglement with cables or tubing constricting airways may cause asphyxiation.

Do not leave long lengths of air tubing or cables around the top of the bed. They could twist around the patient's head or neck while sleeping.

The ventilator is equipped with a rebreathing alarm. The alarm is not a substitute for operator vigilance in ensuring that the leakage port or exhalation valve remains clear at all times. Periodically check the patient circuit during therapy.

In general, as pressure decreases, the potential of rebreathing increases. Lower pressures produce less flow through the leakage port which may not clear all CO<sub>2</sub> from the circuit to prevent rebreathing.

### **WARNING!**

#### **Risk of Cross-Contamination**

Patient circuits might get contaminated by exhaled gases. To avoid cross-contamination, always use a properly cleaned or a new patient circuit when the Vivo 45 LS is to be used by a new patient.

### **NOTE**

For masks and accessories, always follow the manufacturer's instructions.

## **2.5 Usage of Filters**

### **WARNING!**

Always use the ventilator with patient air inlet filters installed. Only use the ventilator with accessories recommended by Breas Medical.



### **Risk of Overheating**

Replace or clean the air inlet filters as specified in the *Maintenance* chapter.

Using old or clogged filters may cause the Vivo 45 LS to operate at higher temperatures than intended.

When operating the Vivo 45 LS, make sure that the air inlet and filters are not obstructed or occluded.



### **Risk of Insufficient Ventilation**

Insufficient ventilation may cause transient hypoxia.

Do not use high resistance bacteria filter at the air outlet of the Vivo 45 LS. High resistance bacteria filters placed between the air outlet and the patient interface may interfere with the operation of the patient disconnect function. It may also interfere with the device trigger function.



### **WARNING!**

#### **Risk of Cross-Contamination**

Deep tissue or mucosal contact with infectious agents may cause infections.

If the Vivo 45 LS is used by several patients, a low resistance bacteria filter shall be used between the air outlet and the patient circuit, for preventing patient cross-contamination. Reuse of bacteria filter, patient circuit or mask may expose patients to contagious agents.

## **2.6**

### **Humidification**



### **WARNING!**

When adding or removing an HME (Heat and Moisture Exchanger, artificial nose) or HCH (Hygrosopic Condenser Humidifier), always reassess the settings, including alarm settings, and perform a pre-use test.



#### **Risk of Suffocation**

When the attachable humidifier is installed, the Vivo 45 LS must be located below the patient and on a flat surface. This is to prevent personal injury from accidental spillage or from excess water or condensation flowing down the patient tube and into the patient's mask. Extra cautions should be taken for patients who are unable to guard their airways or cannot pull off the mask.



When using a humidifier or a nebuliser any patient air filter will need more frequent replacement to prevent increased resistance or blockage.



The ventilator accuracy can be adversely affected by the gas added by the use of a pneumatic nebuliser.



### **WARNING!**

#### **Risk of electric shock**

If using the protective cover or the carry bag, first remove the attachable water chamber. Water spillage may cause electric shocks.



### **WARNING!**

The use of a heated wire patient circuit decreases condensation in the patient circuit.



In case of invasive application, the use of an appropriate external heated humidifier or HME is recommended.



If the condensation in the patient circuit is excessive, the use of a heated humidifier may require the installation of a water trap in the circuit. The water trap prevents any condensed water in the patient circuit from running into the patient airways and causing personal injury.



Any external humidifier connected to the ventilator must comply with ISO 8185 or 80601-2-74.



Any HME connected to the ventilator must comply with ISO 9360.



Do not add any attachments or accessories to the humidifier that are not listed in the instruction for use of the humidifier or the humidifier might not function correctly affecting the quality of the therapy or injuring the patient.



The use of an HME or an external humidifier may require readjustment of the ventilator low-pressure alarm.



Certain HMEs and HCHs are sufficient to provide humidification when the ventilator is used invasively. Check specific suppliers' recommended use.



#### NOTE

The ventilator has been tested and validated with the Fisher & Paykel MR850 heated humidifier.

## 2.7

### Cleaning and Maintenance



#### WARNING!

##### **Risk of Electric Shock**

Cleaning with excessive water or opening the device's casing without certified training may cause electric shocks.

The Vivo 45 LS should be regularly cleaned and maintained in accordance with this operating manual.



#### WARNING!

##### **Risk of Faulty Treatment**

Service and Maintenance of the Vivo 45 LS shall not be performed when the Vivo 45 LS is in use.



#### WARNING!

##### **Risk of Electric Shock**

High voltage contact may cause cardiac arrhythmia.

Repairs and modifications must be carried out by authorized technicians only and in accordance with instructions from Breas Medical

- The Vivo 45 LS must not be opened, repaired or modified by unauthorized personnel or interconnected with incompatible equipment. If subjected to unauthorized modifications or operations, Breas Medical is no longer responsible for the device's performance and safety and all warranties will become invalid.



#### CAUTION!

Do not attempt to autoclave or sterilize the Vivo 45 LS.

## 2.8 Usage of Oxygen

When using the Vivo 45 LS with oxygen, always follow the oxygen provider's instructions and use only medical grade oxygen complying with local regulations.

### **WARNING!**



As this medical device uses an alternative small-bore connector design different from those specified in the ISO80369 series, there is a possibility that a misconnection can occur between this medical device and a medical device using a different alternative small-bore connector, which can result in a hazardous situation causing harm to the patient. Special measures need to be taken by the user to mitigate these reasonable foreseeable risks.



It is the responsibility of the responsible organization to ensure that the oxygen source is compatible with the rated range of pressure, flowrate and oxygen concentration as indicated in the instructions for use as this can affect the performance of the ventilator that can consequently result in patient death or serious deterioration of health.

### **WARNING!**



#### **Risk of fire**

The presence of oxygen can speed up combustion of inflammable materials.



#### **Risk of Fire**

When oxygen is used with the Vivo 45 LS, the oxygen flow must be turned off when the Vivo 45 LS is not operating. Oxygen delivered into the patient tubing may accumulate within the machine enclosure. Oxygen accumulated in the machine enclosure increases the risk of fire.

### **WARNING!**



Do not use a humidifier between the oxygen source and the ventilator, in order to humidify the oxygen flow.  
If humidification is required, use the attachable humidifier or an external humidifier after the patient air outlet.

### **WARNING!**



#### **Risk of Fire**

Ventilate the room adequately. Do not smoke in a room where oxygen is being used.



#### **Risk of Fire**

Naked light bulbs and other sources of ignition must be kept a minimum of 2 meters (6 feet) away from the oxygen cylinder, the patient circuit or any other oxygen carrying parts.



#### **Risk of Fire**

Do not use aerosols or solvents close to the oxygen supply, even when the oxygen supply is shut off.

### **WARNING!**



Supplemental oxygen with a flow up to 30 l/min can be added by an oxygen source with rotameter such as oxygen cylinder, central oxygen supply system or an oxygen concentrator.



### **Risk of faulty Treatment**

At a fixed flow rate of supplemental oxygen flow, the inhaled oxygen concentration will vary, depending on the pressure delivered, the patient's breathing pattern, the patient interface and the leak rate.

To monitor the oxygen concentration, use the FiO<sub>2</sub> sensor accessory.



Supplemental oxygen flow and pressure must not exceed 30 l/min and 100 kPa.



### **CAUTION!**

Supplemental oxygen is added before the volume measurement sensor and thereby included in the measurements. However, the oxygen concentration still has influence on the volume measurement of the delivered air.

This measurement is based on a normal oxygen concentration of 21%. If the oxygen concentration is higher, the actual inspired volume will deviate from the monitored volume as follows:

- 40% oxygen concentration: -2.5% deviation
- 60% oxygen concentration: -5% deviation
- 80% oxygen concentration: -7.5% deviation



## 3 Product Description

### 3.1 Main Components

This section describes the components of the Vivo 45 LS medical electric equipment.

#### NOTE

- There might be local variations of the main components configuration.
- The standard Vivo 45 LS and its packaging do not contain any natural rubber latex.

#### Carry bag

**Function:** Storage for transportation

**Part No:** 006014



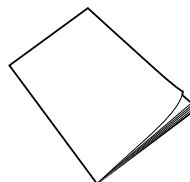
#### Manual

**Function:** Product and usage information

**Part No:**

User's manual: 007300


Clinician's manual: 007298



#### Patient Circuit

**Function:** Delivers air to the patient (applied part)

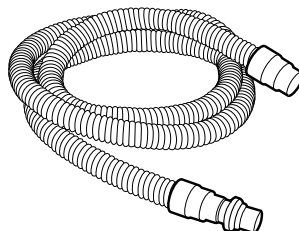Delivered patient circuit depends on the sales configuration. See 9 *Accessories and Parts*, page 144 for compatible patient circuits.



### Patient air inlet filter, fine, white, disposable

**Function:** Fine inlet air filtration.

**Material:** AS 100

**NaCl Penetration:** (0.65  $\mu\text{m}$  NaCl @ 95 l/min) = <7.35%

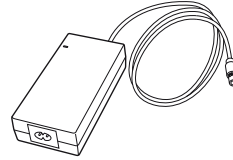
**Part No:** 007103 (5pcs)



### Power Supply

**Function:** Deliver power to the ventilator

**Part No:** 006396



### Power cord

**Function:** Deliver power to the AC power supply

**Part No:**

GB: 003521

CN: 005304

EU: 003520

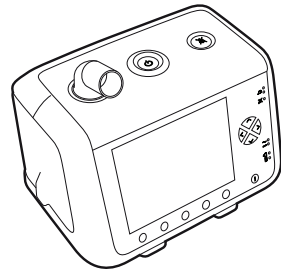
JP: 004834

US: 003522




### Vivo 45 LS Main Unit

Main Unit



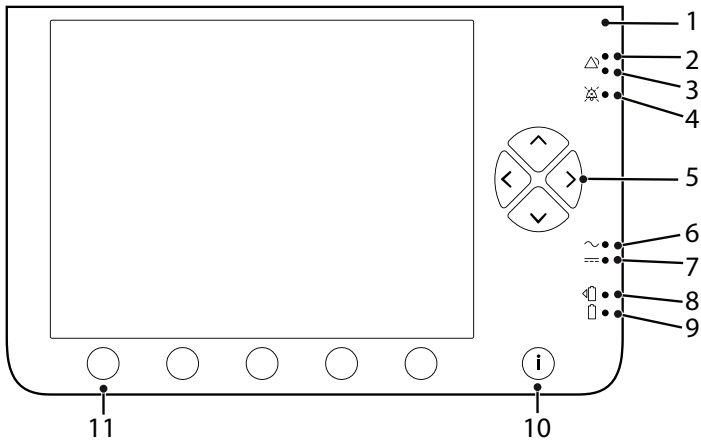
3.2 Front View



| No   | Item                         | Function   |
|------|------------------------------|--|
| 1    | Sensor                       | Ambient light sensor   |
| 2, 3 | Alarm (red & yellow) LEDs    | Alarm indication:<br>Red = High priority<br>Yellow = Medium priority |
| 4    | Audio pause LED              | Paused alarm sound indication  |
| 5    | Navigation /Setting buttons  | Navigate on the display /Define settings                             |
| 6    | Mains LED                    | Power source indication: Mains                                       |
| 7    | External DC LED              | Power source indication: External DC                                 |
| 8    | Click-in battery LED         | Power source indication: Click-in battery                            |
| 9    | Internal battery LED         | Power source indication: Internal battery                            |
| 10   | Information button           | Show / Hide information texts  |
| 11   | Navigation /Function buttons | Functions according to the display                                   |


3.3 Side Views

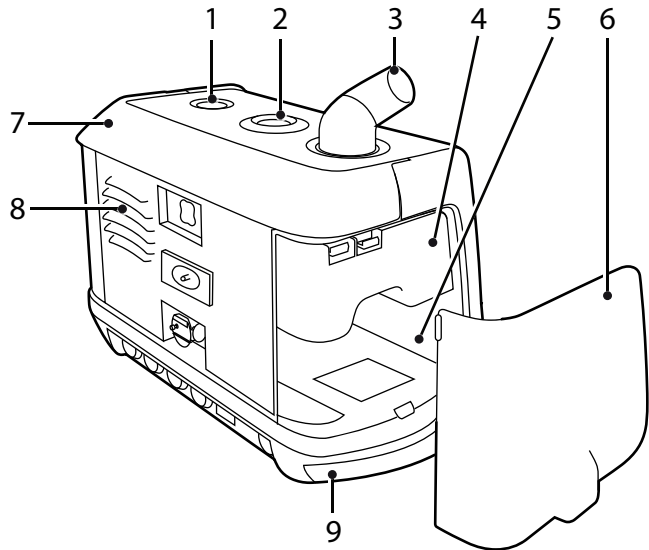Click-in compartment side



| No | Item                 | Function  |
|----|----------------------|---|
| 1  | Audio pause          | Pause the alarm sound   |
| 2  | Start/Stop           | Start/Stop ventilation treatment  |
| 3  | Patient air outlet   | Connection for patient circuit  |
| 4  | Air bypass unit      | Click-in airway/silencer for use without the click-in humidifier. (If the click-in humidifier is used, it replaces the air bypass unit) |
| 5  | Click-in compartment | Compartment for either of the accessories click-in humidifier or click-in battery.  |
| 6  | Side panel           | Cover   |
| 7  | Carrying handle      | Handle for lifting the ventilator   |
| 8  | Cooling air outlet   | Outlet internal cooling   |
| 9  | Cooling air inlet    | Inlet internal cooling  |

### Filter Side

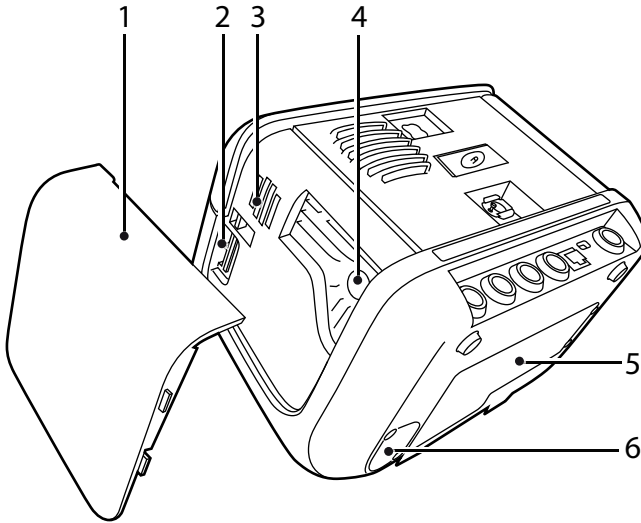


| No | Item                          | Function   |
|----|-------------------------------|--|
| 1  | Side panel                    | Cover  |
| 2  | Memory card slot (SD card)    | Memory download                                      |
| 3  | Alarm beeper                  | Alarm Sounds Output                                  |
| 4  | Patient air inlet             | Air bypass unit in, replaceable filters              |
| 5  | Internal battery              | Compartment for the internal battery                 |
| 6  | FiO <sub>2</sub> sensor hatch | Compartment for the optional FiO <sub>2</sub> sensor |

## 3.3.1 Detaching and Reattaching the Side Panels

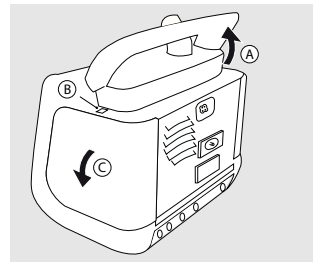
### Detaching the Filter Side Panel

- 1 Lift the handle to access the release button (A).
- 2 Looking from behind, to dismount the filter side panel press the button above the panel (B). The panel is released.
- 3 Remove the panel. (C)




### Reattaching the Filter Side Panel

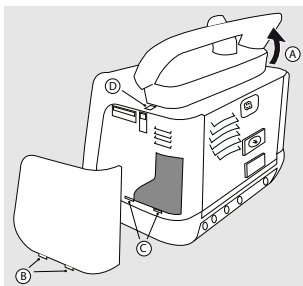- 1 Lift the handle to access the release button (A).
- 2 To mount the filter side panel, insert the tabs (B) on the lower side of the panel into the holes (C).
- 3 Press the side panel into the casing until it clicks in place at the button (D).




### Detaching the Click-in Compartment Side Panel

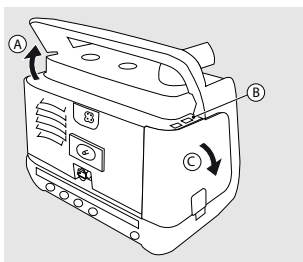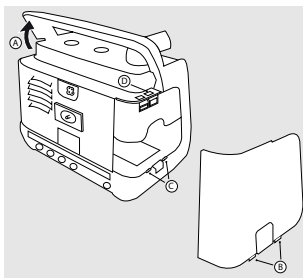- 1 Lift the handle to access the release button (A).
- 2 Press the button marked "1". (B).  
⇒The panel is released.
- 3 Remove the panel (C).

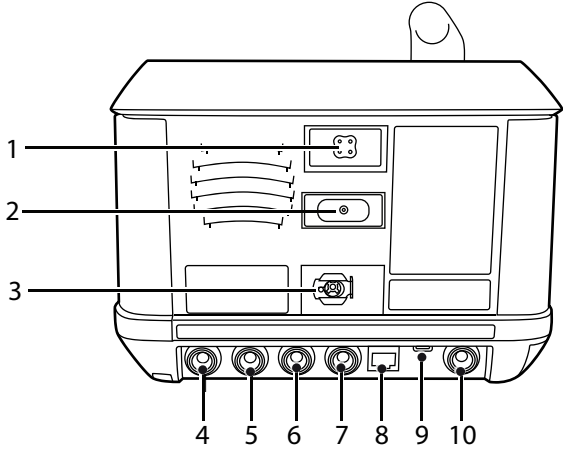



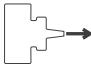











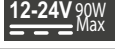

### Reattaching the Click-in Compartment Side Panel

- 1 Lift the handle to access the release button (A).
- 2 To mount the click-in compartment side panel, insert the tabs (B) on the lower side of the panel into the holes (C).
- 3 Press the side panel into the casing until it clicks in place at the button (D).

















3.4 Equipment Designation









| No | Item/Symbol   | Description  | Colour  |
|----|---|--|---|
| 1  |    | Electrical connector, for power to the heated patient circuit. |   |
| 2  |    | Exhalation valve port  |   |
| 3  |    | Connection for low pressure/bleed-in oxygen source             |   |
| 4  |    | CO2 interface port   |    |
| 5  |  | SpO2 interface port  |  |
| 6  |  | Remote start/stop, Audio pause, and effort belt interface port |  |
| 7  |  | Remote alarm and Nurse call interface port                     |  |
| 8  |  | Network connection port  |   |
| 9  |  | USB data connection port                                       |   |
| 10 |  | Mains/External DC inlet  |  |

### 3.4.1 Additional Symbols

This section describes symbols and markings that might appear on the parts or packaging of the Vivo 45 LS.

| Symbol  | Description  |
|---|--|
|    | Internal battery   |
|    | Product number   |
|    | Read user instructions.  |
|    | Attention: Read the intended use in the manual. Read the Safety chapter in the manual.   |
|    | This product must not be exposed to open fire.   |
|    | This product should be recycled.   |
|    | Read 7.7 <i>Disposal</i> , page 128 for information about recycling and disposal.  |
| IP22  | Degree of protection provided by enclosure: IP22. See 8.2.6 <i>Environmental Conditions</i> , page 136 for detailed information.                   |
|    | Manufacturer   |
|    | Serial number  |
|    | This product is a Medical Device.  |
|   | Date of Manufacture  |
|  | IEC protection Class II: Double insulated equipment.   |
|  | Indication of applied parts (IEC 60601-1 Type BF <sub>I</sub> Isolated Applied Part)   |
| Rx Only   | (Symbol only applicable in U.S.) Caution: U.S. Federal law restricts this device to sale by or on the order of a licensed healthcare practitioner. |
|  | Meets all requirements for CE marking according to relevant European health, safety and environmental protection legislation                       |



| Symbol  | Description   |
|---|---|
|  | Meets all requirements for UKCA marking according to applicable United Kingdom health, safety and environmental protection legislation. |
|  | Do not obstruct air inlets or outlet.   |
|  | Single patient multiple use.  |
|  | Single patient multiple use.  |
|  | Hot Surface. Do not touch. (Heating plate in click-in compartment.)   |
|  | MR Unsafe. The device should not enter a magnetic resonance (MR) environment, such as an MRI scanner room.                              |

## 4 Preparing the Vivo 45 LS for Use



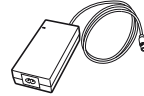
### WARNING!

Read 2 *Safety Information*, page 12 before setting up the Vivo 45 LS.

### 4.1 Checking the Vivo 45 LS before First Use

When using the Vivo 45 LS for the first time, follow the instructions below:


1 Check that all main components and ordered accessories have been delivered (refer to the packing note or the invoice, if available).



2 Ensure that the equipment is in good condition.

3 If stored more than 1 month, connect the Vivo 45 LS to the power supply to recharge the internal battery.

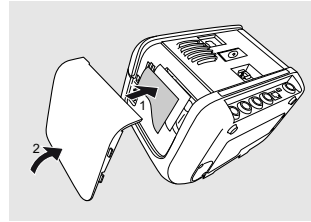4 Check that the grey and white air filters are installed.



### 4.2 Placing the Vivo 45 LS



### WARNING!



Read 2.3 *Environmental Conditions*, page 15 carefully to make sure all conditions are met and considered.

- 1 Place the Vivo 45 LS on a solid, flat, and clean surface.



The Vivo 45 LS should be placed lower than the patient in order to prevent the device from falling on the patient, as well as preventing condensed water from reaching the patient.

Overnight, the Vivo 45 LS should be placed close enough to the patient's bedside to allow movements during the sleep without pulling the Vivo 45 LS off its surface.

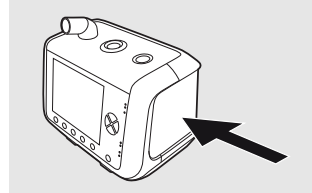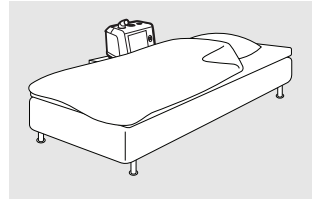- 2 Make sure that nothing can block the patient air inlet.



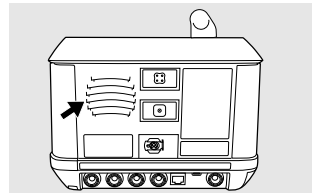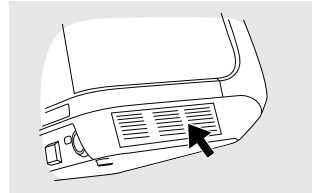- 3 Make sure that nothing can block the cooling air inlet or outlet.




- 4 Make sure that the controls are accessible for the operator.



### CAUTION!

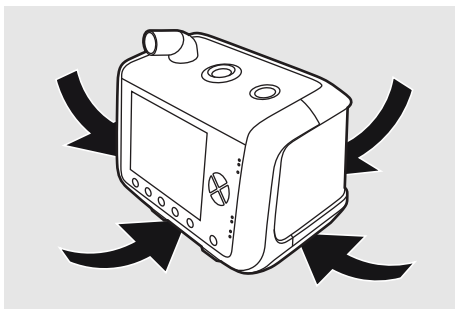- Do not place the Vivo 45 LS on a soft surface that will prevent the air flow underneath the device.  
Never cover the device.



- Always position the Vivo 45 LS so the power supply lays on a surface without strain to the power cord. The power supply shall be easy to disconnect, if required to isolate the Vivo 45 LS from the mains.


## 4.3 Connecting the Vivo 45 LS to Mains




### WARNING!

Read 2.2 *Electrical Safety*, page 15 carefully to make sure all conditions are considered and met.

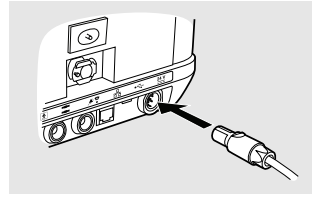
1 Plug the power supply into the power inlet of the Vivo 45 LS.



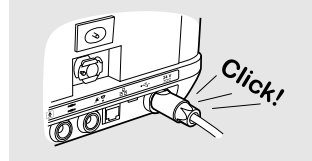2 Make sure a clicking sound is heard to ensure the power supply is completely inserted.



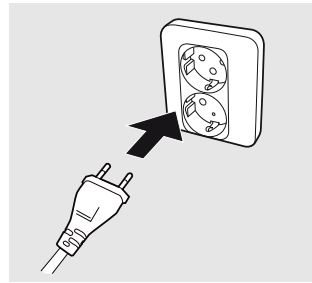3 Connect the power supply's power cord to the mains supply.



To isolate the Vivo 45 LS from the mains supply, disconnect the power supply.

## 4.4 Connecting the Patient Circuit



### **WARNING!**


Read 2.4 *Usage of Patient Circuit*, page 16 carefully to make sure all conditions are considered and met.

The Vivo 45 LS can be used with the following circuits:


- Single limb circuit with external active exhalation valve
- Single limb circuit with external leakage port
- Single limb circuit connected to patient interface with integrated leakage port
- Circuit with mouthpiece interface

## Connect the Patient Circuit

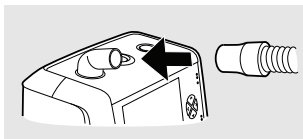
- 1 Inspect the circuit for damages.
- 2 Connect the patient circuit to the patient air outlet on the ventilator.



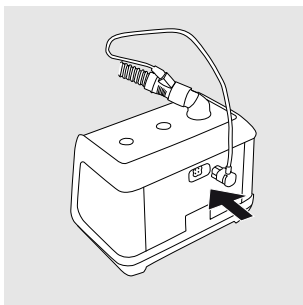- 3 If having a heated patient circuit, connect the heated wire electrical plug to the socket on the ventilator.



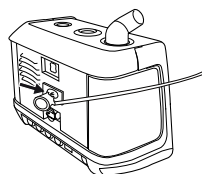- 4 If having an active exhalation valve patient circuit, connect the pilot pressure tube at the back of the ventilator.



- 5 Connect the other end of the patient circuit to the patient interface.

If having a leakage circuit and a patient interface without integrated leakage port: Make sure to use a leakage port between the circuit and the patient interface.

- 6 If using a leakage circuit or an active exhalation circuit, make the correct setting for the type of circuit. See .

For MPV circuits, the type of circuit is set when activating the MPV mode.

## 4.5 Inspecting the Vivo 45 LS before Use

### Inspection of Device

- Check that there is no visible damage.
- Check that the surface is clean.

### Inspection of Cables

- Check that all cables are recommended by Breas.
- Check that the cables are undamaged.
- Check that the cables are properly connected.

## Inspection of Placement

- The Vivo 45 LS shall be placed on solid flat surface below the patient level (see 4.2 *Placing the Vivo 45 LS*, page 32).
- Make sure that nothing can block the air inlet at the side.

## Inspection at Ventilator Startup

This procedure checks the ventilator's alarm handling and power source handling. If any check fails, take the ventilator out of use and contact your service provider.

- 1 Connect a patient circuit to the ventilator
- 2 Connect the power supply to the ventilator. If the power is supply is connected and the ventilator is turned off, press the Start/Stop button.

⇒The ventilator now turns on and enters stand by mode. If it is the first time the ventilator is turned on, you also have to select language.

- 3 If needed, perform a pre-use test.
- 4 Press and hold the Start/Stop button until the progress bar is filled to start the treatment.

⇒At the start of the treatment the ventilator performs an alarm test. Check that:

- The alarm LEDs flash
- The ventilator beeps

- 5 Disconnect the power supply for more than 5 seconds.

⇒The ventilator now switches to the internal battery (or click-in battery, if connected). Check that *Lost mains power* alarm is given.

- 6 Reconnect the power supply.

Check that the device switches to the mains supply, indicated by an information message and a beep.

- 7 Ensure that the treatment settings and alarm settings are set as prescribed before taking the ventilator to use.

## 4.6

### Adjusting the Vivo 45 LS Patient Settings

#### WARNING!



The configuration of the Vivo 45 LS therapy settings must always be prescribed by a licensed physician and carried out by an authorised health care professional.

For detailed information about the treatment parameters of the Vivo 45 LS, see 5.5 *Functions and Parameters in the Vivo 45 LS*, page 59.

Follow the instructions below when setting up the Vivo 45 LS:

- Adjust the settings to find the best possible breathing comfort for the patient.
- If you have changed the ventilation mode, press **Select Next** and review the settings before pressing **Confirm**.
- Always document the patient settings.
- The ventilator always starts in the mode and with the settings that were active when it was switched off.

## 4.7

### Performing the Pre-use Test

The pre-use test is used for detecting the type and characteristics of the patient circuit that is connected to the Vivo 45 LS. The resistance and compliance of the patient circuit are

measured and calculated. This will be used to compensate for pressure drop in the patient circuit and the compliance of the patient circuit.

The patient shall not be connected during the pre-use test.

#### NOTE

If the pre-use test has not been performed, the Vivo 45 LS will operate with default patient circuit compensation.



### Starting the Pre-use Test Manually

- 1 On the **Others** menu, Select **Pre-use Test** and then **Start Pre-use Test**.
- 2 Follow the instructions on the display.

### Activating the Pre-use Test Prompt

- 1 On the **Others** menu, Select **Pre-use Test**.
- 2 Set **Pre-use Test** to **On**.

### Pre-use Test Sequence

When performing a pre-use test, the instructions on the display will guide you through this sequence:

| Step | Action   |
|------|--|
| 1    | Start of pre-use test.   |
| 2    | Connect the patient circuit.   |
| 3    | Make sure that nothing is blocking the patient end of the circuit.   |
| 4    | Wait while the Vivo 45 LS is checking the patient circuit resistance. If the resistance is not within the limits, the test will end without performing the following steps. The result will be displayed for review. |
| 5    | Block the end of the patient circuit with an air tight object.   |
| 6    | Wait while the Vivo 45 LS is checking the patient circuit compliance and leakage.  |
| 7    | Test finished. Review the test result.   |

#### 4.7.1 Actions At Pre-Use Test Failure

At the end of the pre-use test the individual results for leakage, resistance and compliance are shown.

#### Failure Due To Incorrect Leakage

Indication: **Leakage: Fail**

- 1 Check the ventilator set-up (circuit, filters, humidifier, etc.) for leakage.
- 2 Ensure that all connectors are tightly fitting.
- 3 Run the pre-use test again.
- 4 Replace the circuit if the test is failed repeatedly.



## Failure Due To Incorrect Resistance or Compliance

Indication: **Resistance: Fail** or **Compliance: Fail**

**1** Check the ventilator set-up (circuit, filters, humidifier, etc.) for blockage or pinched tubing.

**2** Run the pre-use test again.

If the pre-use test is continually failed due to resistance or compliance, it is permitted to use the ventilator but be aware that the pressure (resistance) or volume (compliance) delivered to the patient may not meet with the specified accuracy.

The ventilator will apply the default values to compensate for circuit resistance and compliance. These values will deviate from the values for the circuit in use.

Ensure that the delivered ventilation is closely monitored.

## 5 How to Use the Vivo 45 LS



### WARNING!

Read 2 *Safety Information*, page 12 before using the ventilator.

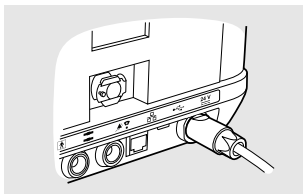
When the ventilator is handed over to the patient, the physician in charge or hospital staff must instruct the patient in how the unit works.

### 5.1 Switch the Vivo 45 LS On and Off

#### Switch On and Enter Standby Mode


1 If having access to mains power, connect the mains power supply.

⇒The Vivo 45 LS now starts up and enters standby mode.



2 If running the Vivo 45 LS on the internal battery or the click-in battery, press the Start/Stop button.

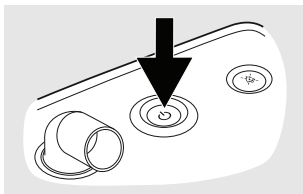
⇒The Vivo 45 LS now starts up and enters standby mode.



3 Select Yes/No when asked to “Perform Pre-use Test”. (The Vivo 45 LS may be configured to not ask for pre-use test at start up. See 5.2.9 *The Others Section* , page 49.

#### Start the Treatment

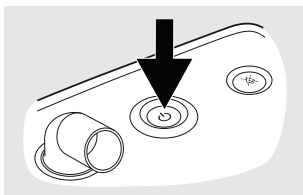1 Press and hold the Start/Stop button on the Vivo 45 LS.




2 Release the Start/Stop button when the progress bar is filled.

**Tip:** You can also press the start button shortly and then press the “Yes” button when asked to start the treatment.


## Stop Treatment and Switching Off

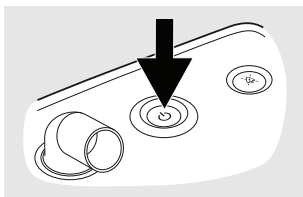1 To stop treatment and enter standby mode, first press and hold the Start/Stop button on the front panel.



2 Release the Start/Stop button when the progress bar is filled.

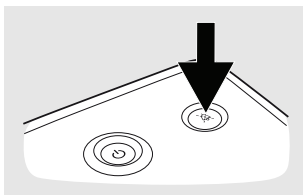
3 Press “Ok” or the Audio Pause button to stop the treatment. To cancel the stopping process and continue the treatment, press “Cancel”.




## 5.2 Using the Menu

### 5.2.1 Navigating with the Buttons

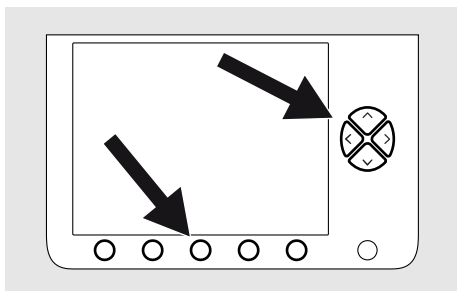Use the five navigation buttons and the up and down buttons on the panel in order to navigate the menu.



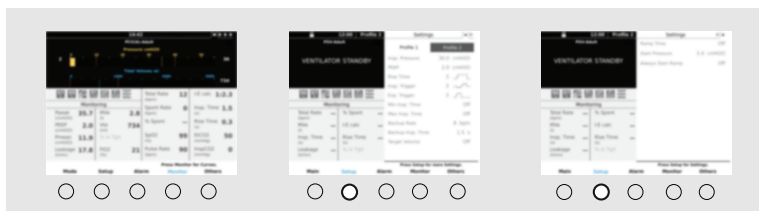The navigation buttons are used to view the different sections defined above each navigation button. The same navigation button can also be used to view additional information in some sections, or it can be designated a temporary function while an event window is active.




- The currently displayed section is indicated by a title in blue text.
- For sections with multiple pages, press the navigation button repeatedly to circle between the pages. Multiple pages are indicated by dotted tabs in the upper right corner.



## NOTE

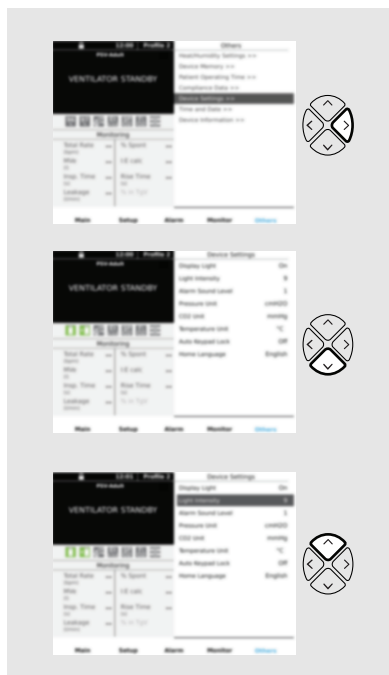Use the up or down button to enter the menu list.

Use the up and down buttons to navigate up or down in a menu list, or to select different parameters. Use the left and right buttons to alter parameters, or enter and exit sub sections.



## 5.2.2 Symbols Used in the Menu



**Internal battery.**

For detailed information about charge level indications, see 5.6.4 *Battery Icons* , page 70.



#### Click-in battery level

For detailed information about charge level indications, see 5.6.4 *Battery Icons* , page 70.



#### Home Mode activated



#### Connected to nurse call system



#### Remote alarm unit connected



#### Humidifier and/or heated circuit connected

The number in the drop symbol indicates the humidity setting. If the humidifier is connected but not activated, the drop is crossed over. The temperature above the heating waves symbol indicates the temperature setting for the heated circuit. If the heated circuit is connected but not activated, the temperature is replaced by “Off”.



#### Single limb circuit with leakage port selected (Leakage)



#### Single limb circuit with active exhalation valve selected (Exh. valve)



#### MPV breath mode selected.

Circuit with mouthpiece interface shall be used.



#### SpO2 sensor connected



#### FiO2 connected



#### EtCO2 connected



#### PtcCO2 connected



#### Effort belt connected

Connection failures are indicated by red coloured belts.



#### Multiple pages



#### High priority alarm event in history list

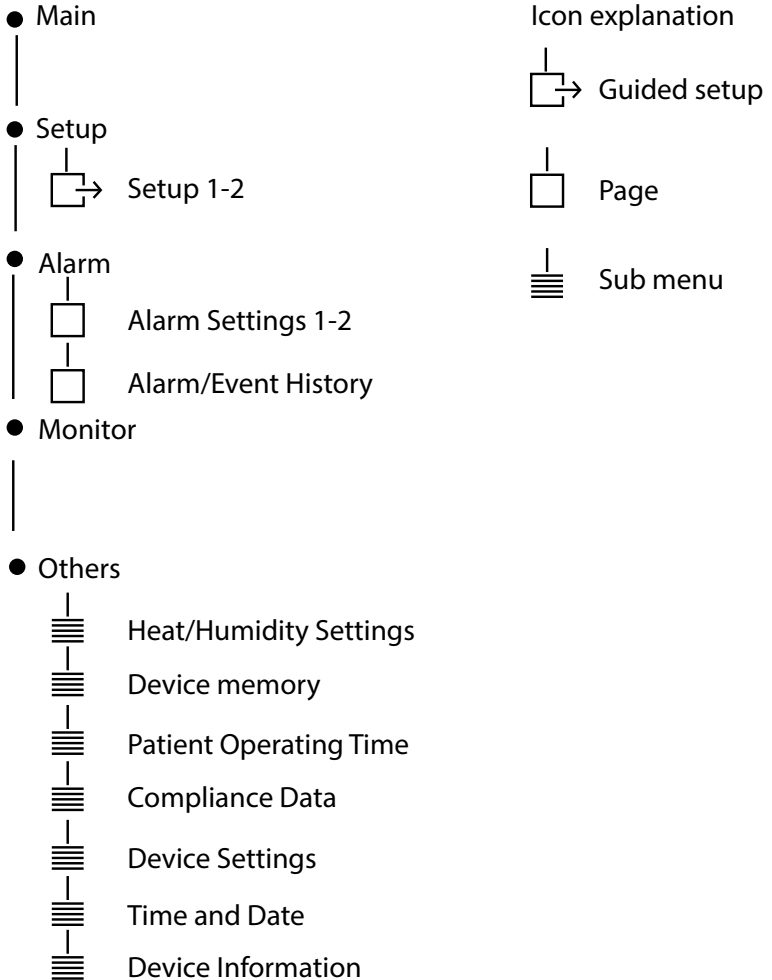


### Medium priority alarm event in history list

Size and colour of symbols may vary depending on displayed view.


## 5.2.3 Menu Overview

In *Home Mode*, the Vivo 45 LS menu has the following section layout:



## 5.2.4 Display Overview

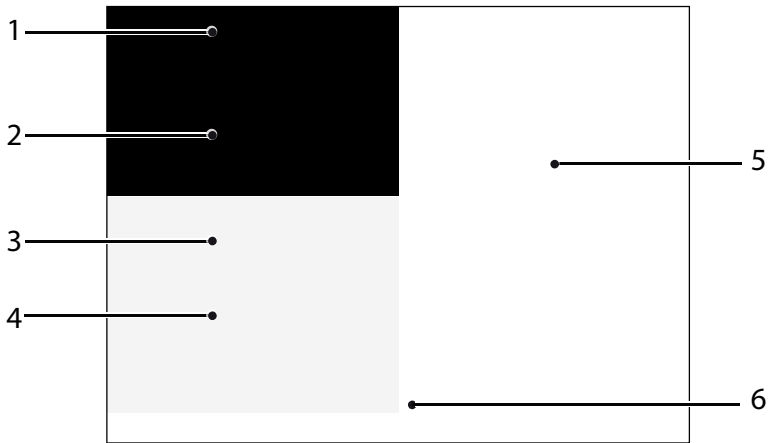This section describes the areas on the display.



### 1. Title and Icons

Ventilation mode, breath mode, patient mode and device mode are displayed. When more than 1 profile is used, the active profile is displayed.

### 2. Pressure and Volume Indicators

The bar graphs are used to display current pressure, PEEP, Ppeak, pressure and volume alarm limits, and inspired/expired tidal volume.

The red lines represent the low and high pressure alarms, the low and high tidal volume alarms.

### 3. Icons/Alarm Message

Information icons are presented here to give a quick overview of the Vivo 45 LS basic status (see 5.2.2 *Symbols Used in the Menu*, page 42).

### 4. Monitoring Field

This field is designed to give a brief overview of the key monitor values. Navigate to the Monitor section for the all available values.

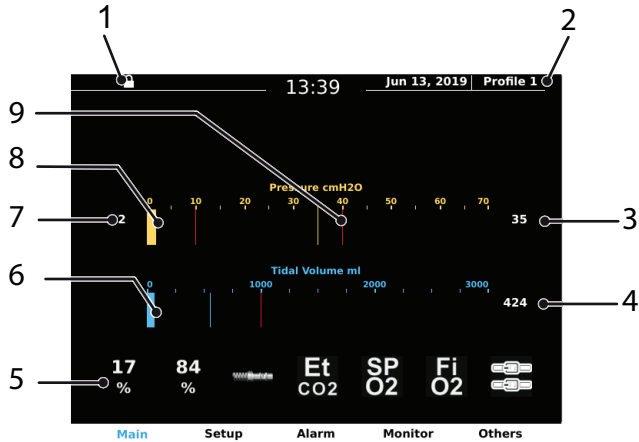
### 5. Screen Title and Context Area

This field displays settings and information for the section. Screen title, page number (if more than one sub page exists in the section) and time are also displayed.

### 6. Navigation Field

This field is mainly used for displaying the section layout of the menu, and determines the function of each navigation button. Depending on current operation, the navigation buttons can be assigned temporary functions such as “Yes” or “Next”, depending on which operation is active.

## 5.2.5 The Main Display, Home Mode



1. Indication of home mode lock status. See 5.2.2 *Symbols Used in the Menu*, page 42 for information about each icon.
2. When more than 1 profile is used, the active profile is displayed.
3. Ppeak (max pressure during last breath)
4. Tidal volume during last breath
5. Indication of battery status and connected sensors. See 5.2.2 *Symbols Used in the Menu*, page 42 for information about each icon.
6. Bar graph representing the tidal volume
7. PEEP value
8. Bar graph representing the current pressure
9. Red lines in the bar graphs represent alarm limits



5.2.6 The Setup Section, Home Mode



In the Setup section, treatment parameters can be viewed. See 5.5 *Functions and Parameters in the Vivo 45 LS*, page 59.

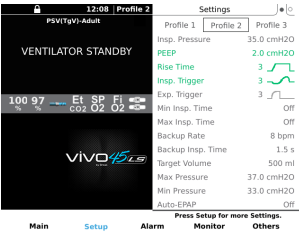
If Home Adjust is activated, some parameters might be adjustable within ranges defined by the prescribing physician. Adjustable parameters are indicated by green text.

If profiles with preset treatment settings have been created, the profiles are selected from the Setup section. The currently active profile is indicated by a frame.

Selecting a Profile

This procedure describes how to select between preset treatment profiles in home mode, if more than one profile has been created by the clinician.

- 1 Press the **Setup** button to display the Setup section.
- 2 Press the Down or Up arrow button to enter the Profile selection.
- 3 Press the Right or Left arrow to select a profile.  
The treatment settings for the profile are displayed below the profile selection.




- 4 Confirm the change of Profile by pressing the **Yes** button.

5.2.7 The Alarm Section



The Alarms section contains 3 pages: The two first are for viewing alarm settings and the third is for viewing the alarm and message history. Shift between the pages by pressing the Alarm button.

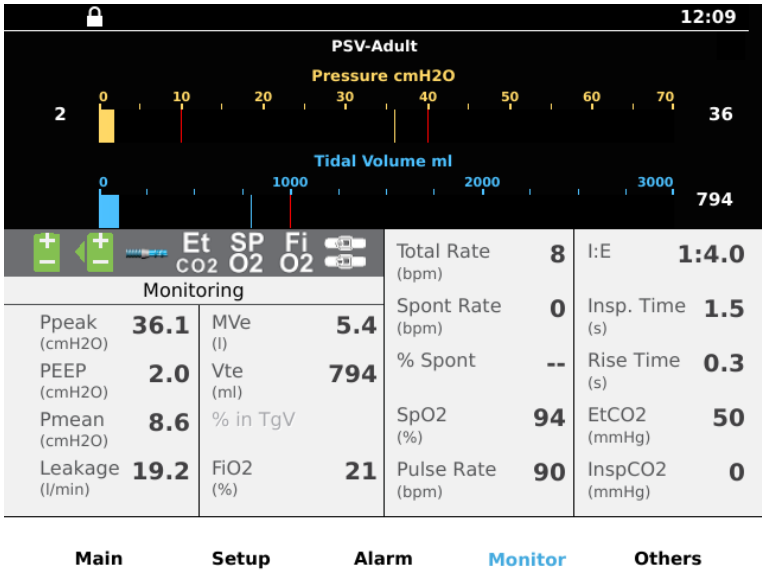5.2.8 The Monitor Section



The monitoring page contains bar graphs for pressure and volume, and a field with measurement readings for the monitored values that are applicable for the treatment mode. See 5.3 *Monitored Values in the Vivo 45 LS*, page 50 for a description on the monitored values.

5.2.9 The Others Section

5.2.9.1 Heat/Humidity Settings

- Humidifier Setting - Sets the level of humidity.
- Humidifier - Activates or deactivates the humidifier.
- Heated Circuit Temp - Sets the temperature of the heated circuit.
- Circuit Heating - Activates or deactivates the circuit heating.

5.2.9.2 Device Memory

The Device memory setting is used for transferring data to a memory card and for erasing the device data (requires Clinical mode).

5.2.9.3 Patient Operating Time

Shows the number of hours a patient has been using the Vivo 45 LS for breathing therapy.

5.2.9.4 Device Settings

General settings for the Vivo 45 LS:

- Display Light: On (will keep the display lit up regardless of use), Auto (will adjust the light intensity depending on the ambient light), Delayed (the display is dimmed after 30 seconds or more depending on the mode and battery setup. If any button is pressed or any alarm occurs, the display light will return to normal again).
- Light Intensity (setting range: 1-9, where 1 is the lowest and 9 is the highest light intensity setting. In cases where Display Light is set to “Auto”, the Light Intensity setting will not be available).  
If the light intensity is set too low, ambient light could cause difficulties in reading the alarm text.
- Auto Keypad Lock (On, Off).
- Home Language: Selects the language for Home Mode.

#### 5.2.9.5 Time and Date

- Time (set time: hours and minutes).
- Time Format (choose between 24 h or 12 h am/pm format).
- Date Format (choose between yyyy-mm-dd, dd/mm/yyyy, mm/dd/yyyy format).
- Date (set date: year, month and day).

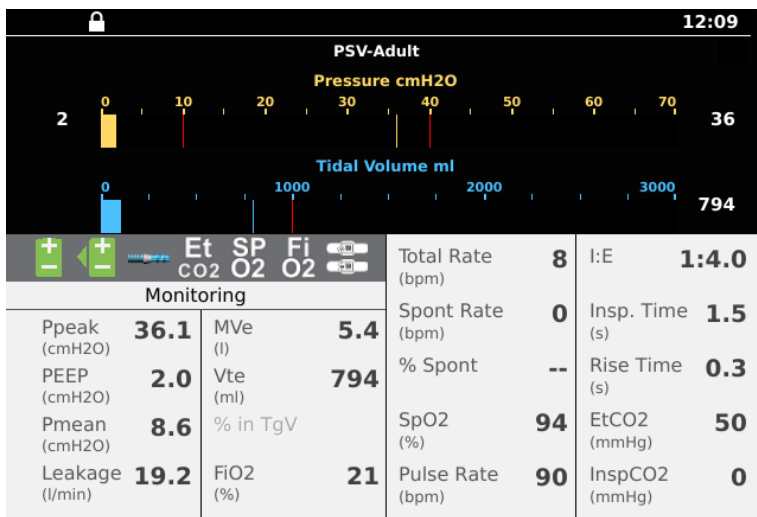
#### 5.2.9.6 Device Information

- Device Operating Time (hours)
- Firmware Version
- Language: Clinical Mode
- Language: Home Mode
- ~ AC (On/Off)
- External DC (V)
- Battery Health (%)
- Serial Number

### 5.3 Monitored Values in the Vivo 45 LS


Values monitored by the Vivo 45 LS stored in the data log and displayed on the screen:

1. The monitoring screen:



**Main**      **Setup**      **Alarm**      **Monitor**      **Others**

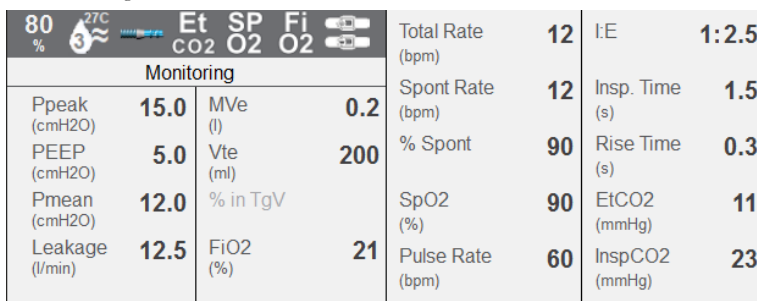2. The monitoring field:



### 5.3.1 Ppeak

Ppeak displays the highest pressure that is recorded during the inspiratory phase.

### 5.3.2 PEEP

PEEP displays the pressure at the end of the expiratory phase.

### 5.3.3 Pmean

Pmean displays the calculated mean value of pressure during a complete ventilatory cycle (inspiratory + expiratory phase).

### 5.3.4 Leakage

Leakage displays the total leakage (intentional and unintentional) as calculated at expiratory pressure level

### **5.3.5 MVe**

MVe displays the expired Minute Volume calculated as expired Tidal Volume multiplied with the Total Breath Rate.

MVe is displayed when a leakage circuit is used.

### **5.3.6 MVi**

MVi displays the inspired Minute Volume calculated as inspired Tidal Volume multiplied with the Total Breath Rate.

MVi is displayed when an MPV circuit or an active exhalation valve circuit is used.

### **5.3.7 Vte**

Vte displays the Expired Tidal Volume that the patient exhales during each breath.

Vte is displayed when using a leakage patient circuit . Vte is a calculated value.

### **5.3.8 Vti**

Vti displays the inspired Tidal Volume that the ventilator delivers during each breath.

Vti is displayed when an MPV circuit or an active exhalation valve circuit is used.

### **5.3.9 FiO<sub>2</sub>**

The FiO<sub>2</sub> displays the fraction of inspired oxygen as measured in the air channel of the ventilator. An FiO<sub>2</sub> sensor needs to be in place to measure and display this value (see the section on using the ventilator with the FiO<sub>2</sub> sensor.)

### **5.3.10 % in TgV**

% in TgV displays the percentage of breaths where the actual delivered Tidal Volume matches with the set Target Volume (not calculated until 100 breaths are registered).

### **5.3.11 Total Rate**

The Total Rate displays the actual total breath rate independent of whether the breaths are patient- or ventilator-triggered breaths.

### **5.3.12 Spont Rate**

The Spont Rate displays the actual spontaneous breath rate.

### **5.3.13 % Spont**

% Spont displays the percentage of spontaneous breaths calculated since the ventilator was last started (not calculated until 100 breaths are registered).

### 5.3.14 SpO<sub>2</sub> (Oxygen Saturation)

The SpO<sub>2</sub> displays the patient's oxygen saturation, if measured with the SpO<sub>2</sub> module accessory.

- When using the ventilator with the SpO<sub>2</sub> sensor accessory, the ventilator displays functional oxygen saturation measured by the sensor.
- The following information concerns the light emitted by the SpO<sub>2</sub> Sensor:
  1. Peak Wavelength (red): 660 nm
  2. Peak Wavelength (infrared): 905 nm
  3. Maximum Optical Output Power: 15 mW
  4. For more information regarding the oxygen probe's range of peak wavelengths, max optical power and usage, please refer to the respective probe manual.
- Environmental factors may influence the function or accuracy of the pulse oximeter, such as ambient light, physical movement, diagnostic testing, low perfusion, electromagnetic interference, dysfunctional haemoglobin, presence of certain dyes and inappropriate positioning of the pulse oximeter probe.
- A functional test cannot be used to assess the Accuracy of a Pulse Oximeter Probe or a Pulse Oximeter.

### 5.3.15 Pulse Rate

Pulse Rate displays the patient's pulse rate as measured with the SpO<sub>2</sub> module.

### 5.3.16 I:E

I:E displays the ratio between the length of the inspiration and the length of the expiration.

### 5.3.17 Insp. Time

Insp. Time displays the duration of the inspiratory cycle, measured from the start of inspiration to the start of expiration.

### 5.3.18 Rise Time

Rise Time displays the duration of the pressure or volume increase, measured from the start of inspiration until the set pressure or volume has been reached.

### 5.3.19 EtCO<sub>2</sub>

EtCO<sub>2</sub> displays the end-tidal carbon dioxide, measured on the last portion of the exhaled volume that is passing through the EtCO<sub>2</sub> sensor.

### 5.3.20 InspCO<sub>2</sub>

InspCO<sub>2</sub> displays the inspired carbon dioxide.

### 5.3.21 PtcCO<sub>2</sub>

PtcCO<sub>2</sub> displays transcutaneous CO<sub>2</sub> pressure from an external PtcCO<sub>2</sub> monitor, if connected by the PtcCO<sub>2</sub> cable accessory.

## 5.4 Modes in the Vivo 45 LS

In the Modes section of the Vivo 45 LS display, the operator selects the ventilation mode, breath mode, patient mode and device mode for the treatment.

### 5.4.1 Standby and Operating Mode

Standby mode is defined as the state of the Vivo 45 LS when it is powered on but without having treatment or pre-use test started.

Operating mode is defined as the state of the Vivo 45 LS when the blower is operating and producing an air flow.

Switch between operating and standby mode by starting/stopping the ventilator (see 5.1 *Switch the Vivo 45 LS On and Off*, page 40).

Some operations (such as setting time and date) are only available in standby mode.

### 5.4.2 Device Mode

The two device modes of the ventilator are used for controlling the user access to the ventilator settings.

Clinical mode allows full access to the ventilator treatment parameters for health care professionals.

Home mode is used to limit the access to the ventilator's settings for patients and lay persons.

If the Home Adjust function is activated, the patient is given the possibility to change a selection of patient parameters that is determined by the responsible care provider. The limited settings for Home Mode are defined in the Setup section



See 5.2.3 *Menu Overview*, page 44 for a chart with settings available in the ventilator's Home Mode.

### 5.4.3 Patient Mode

The Vivo 45 LS can be operated in either Adult mode or Paediatric mode.

In Paediatric mode, some ventilator parameters, for example Breath Rate, Inspiratory Time and Target Volume have special setting limits, in order to optimize the treatment for paediatric patients. The settings and alarm limits of the Vivo 45 LS are reset to default values when the patient mode is altered.

### 5.4.4 Ventilation and Breath Modes

The ventilation and breath modes are used for controlling the ventilation treatment with the ventilator. The ventilation mode selected can be either Pressure, Volume or CPAP.

Pressure mode is used in combination with the Support, Assist/Control, MPV or SIMV breath mode.

Volume mode is used in combination with the Assist/Control, MPV or SIMV breath mode.

The following combinations of ventilation and breath modes can be selected for the Vivo 45 LS:



- PSV (Pressure Support Ventilation). May be combined with:
  - TgV (Target Volume)
  - AE (Auto-EPAP)
- PCV (Pressure Controlled Ventilation). May be combined with:
  - A (Assist function — inspiration trigger)
  - TgV (Target Volume)
  - AE (Auto-EPAP)
  - MPV (Mouthpiece Ventilation)
  - SIMV (Synchronized Intermittent Mandatory Ventilation)
- VCV (Volume Controlled Ventilation). May be combined with:
  - A (Assist function — inspiration trigger)
  - MPV (Mouthpiece Ventilation)
  - SIMV (Synchronized Intermittent Mandatory Ventilation)
- CPAP (Continuous Positive Airway Pressure)

#### 5.4.4.1 Pressure Support Ventilation (PSV)

In the PSV mode, the patient's spontaneous breathing is assisted by the ventilator. The patient controls the start of inspiration through the inspiratory trigger and the start of exhalation by the expiratory trigger.

The set pressure is used as a target pressure, if the flow is decreased to the expiratory trigger level before the set pressure is reached, the expiration starts.

When an inspiration is started, either when the patient triggers a breath, or when the backup rate setting initiates an inspiration in case of a prolonged Apnea, the ventilator delivers a flow up to a certain preset pressure limit. In case of a patient initiated breath, the patient continues the breath for as long as they wish and cycle off when a percentage drop in peak inspiratory flow (expiratory trigger) has been reached.

Spontaneous breaths stop and an exhalation starts in three cases:

- The inspiration flow has dropped to the value set for expiratory trigger.
- The inspiration time is longer than the limit for maximal inspiration time or when inspiration time reaches 3 seconds.
- The limit for the high-pressure alarm is reached.

#### 5.4.4.2 Pressure Support Ventilation with Target Volume (PSV+ TgV)

The PSV+TgV mode acts as the PSV mode but with an added regulation of the pressure. Target volume is a feature that automatically adapts the pressure to make sure that the Vivo 45 LS delivers the desired set target volume to the patient. The delivered volume is compared to the set target volume on a breath by breath basis. The delivered pressure for the next breath will be increased or decreased depending on the difference between the delivered volume and the set target volume. Automatic pressure adjustments will be made in between two adjustable limits (min IPAP and max IPAP) in order to deliver the optimal support to the patient



See 5.5.16 *Target Volume*, page 66 for more information about Target Volume.

#### 5.4.4.3 Pressure Controlled Ventilation (PCV)

In the PCV mode the ventilation is controlled by the Vivo 45 LS. This is done at the preset pressure, breath rate, inspiratory time, and rise time settings set by the operator.

The inspiration stops and an expiration starts in two cases:

- The inspiration time expires.
- The limit for the high pressure alarm is reached.

#### 5.4.4.4 Pressure Controlled Ventilation with Target Volume (PCV+TgV)

The PCV+TgV mode acts as the PCV mode, but with an added regulation of the pressure. Target volume is a feature that automatically adapts the pressure to make sure that the Vivo 45 LS delivers the desired set target volume to the patient. For every breath, the delivered volume is compared to the set target volume. The delivered pressure for the next breath will be increased or decreased depending on the difference between the delivered volume and the set target volume. Automatic pressure adjustments will be made in between two adjustable limits (min pressure and max pressure) in order to deliver the optimal support to the patient.



See 5.5.16 *Target Volume*, page 66 for more information about Target Volume.

#### 5.4.4.5 PCV(A) – Assisted Pressure Controlled Ventilation

In the PCV(A) mode the ventilation is controlled by the ventilator, but the patient has the possibility to start a breath through the Inspiratory Trigger. This patient triggered breath will be delivered with the inspiratory time, rise time and pressure setting set by the operator.

#### 5.4.4.6 PCV(A+TgV) – Assisted Pressure Controlled Ventilation with Target Volume

The PCV(A+TgV) mode acts like the PCV(A) mode, but with an added regulation of the pressure. Target volume is a feature that automatically adapts the pressure to make sure that the ventilator delivers the desired set target volume to the patient. The delivered volume is compared to the set target volume on a breath by breath basis. The delivered pressure for the next breath will be increased or decreased depending on the difference between the delivered volume and the set target volume. Automatic pressure adjustments will be made in between two adjustable limits (min pressure and max pressure) in order to deliver the optimal support to the patient.

#### 5.4.4.7 PCV-SIMV – Pressure Controlled Ventilation with Synchronized Intermittent Mandatory Ventilation

In the PCV-SIMV mode, the ventilator synchronizes mandatory pressure-controlled breaths with the patient's breathing efforts.

In this mode, the ventilator delivers mandatory pressure-controlled breaths with a preset breath frequency defined as the SIMV Rate. The SIMV Rate setting determines the length of the SIMV cycle.

For each SIMV cycle, there is an initial mandatory period in which the patient may trigger one mandatory breath. This mandatory period is always 80% of the SIMV cycle time. If the patient does not trigger a breath during this period, the ventilator will deliver one mandatory breath automatically in the end of the period.

Between the mandatory breaths, the patient may trigger spontaneous breaths until the next SIMV cycle begins. The ventilator will respond to the patient's inspiration efforts with additional pressure support breaths. The inspiratory pressure of these support breaths is defined

by the SIMV Support Pressure, together with the settings for Rise Time and Expiratory Trigger. The default value for SIMV Support Pressure is the Inspiratory pressure in PCV.

After a mandatory breath, Vivo 45 LS will always wait at least one second before a new mandatory breath can be initiated in the following SIMV cycle.

#### **5.4.4.8 PCV-MPV– Pressure Controlled MouthPiece Ventilation**

##### **Mouthpiece Ventilation**

This mode is tailored specifically for those patients that use mouthpiece interface together with pressure controlled ventilation.

It is available in **Adult** patient mode only.

In MPV breath mode it is possible to set Breath Rate to zero so that breaths are only initiated when the patient triggers them using the mouthpiece. The PEEP setting is not available but always set to Off so that no air blows from the mouthpiece when no breath is being delivered.

As the patient is not always connected to the ventilator several of the surveillance functions are not working in the same way as in other modes:

- Disconnection alarm is not available in MPV mode.
- Low Pressure alarm is only active during breaths.
- Apnea alarm is possible to set to longer times and is an important mean of surveillance to make sure that the patient gets ventilation support regularly.

When switching to MPV mode all alarms, except from High/Low Pressure Alarms, are automatically switched off to avoid false alarm triggering (when changing between profiles the alarm settings does not change to Off, but stays as defined in the profiles)

- The alarm levels must be assessed and adjusted considering the patient's condition and treatment settings.
- MPV breath mode shall be used with the mouthpiece circuit including the mouthpiece interface only.
- The MPV circuit including the mouthpiece interface shall be used in MPV breath mode only and is not for use in any other breath mode.
- MPV breath mode shall not be used with ventilator dependent patients.

#### **5.4.4.9 VCV – Volume Controlled Ventilation**

In the VCV mode the ventilation is controlled by the Vivo 45 LS. This is done with the pre-set tidal volume, breath rate, inspiratory time, and rise time settings set by the operator.

The inspiration stops and an expiration starts in two cases:

- The inspiratory time expires.
- The limit for the high pressure alarm is reached.

#### **5.4.4.10 VCV(A) – Assisted Volume Controlled Ventilation**

In the VCV(A) mode the ventilation is controlled by the Vivo 45 LS, but the patient has the possibility to start a breath through the inspiratory trigger. This patient triggered breath will be delivered with the inspiratory time, rise time and tidal volume setting set by the operator.

#### **5.4.4.11 VCV-SIMV – Volume Controlled Ventilation with Synchronized Intermittent Mandatory Ventilation**

In the VCV-SIMV mode, the Vivo 45 LS synchronizes mandatory volume-controlled breaths with the patient's breathing efforts.

In this mode, the ventilator delivers mandatory volume-controlled breaths with a preset breath frequency defined as the SIMV Rate. The SIMV Rate setting determines the length of the SIMV cycle.

For each SIMV cycle, there is an initial mandatory period in which the patient may trigger one mandatory breath. This mandatory period is always 80% of the SIMV cycle time. If the patient does not trigger a breath during this period, the ventilator will deliver one mandatory breath automatically in the end of the period.

Between the mandatory breaths, the patient may trigger spontaneous breaths until the next SIMV cycle begins. The ventilator will respond to the patient's inspiration efforts with additional pressure support breaths. The inspiratory pressure of these support breaths is defined by the SIMV Support Pressure, together with the settings for Rise Time and Expiratory Trigger. The default value for SIMV Support Pressure is the Inspiratory pressure set in PCV or PSV.

After a mandatory breath, Vivo 45 LS will always wait at least one second before a new mandatory breath can be initiated in the following SIMV cycle.

#### 5.4.4.12 VCV-MPV– Volume Controlled MouthPiece Ventilation

This mode is tailored specifically for those patients that use mouthpiece interface together with pressure controlled ventilation.

It is available in **Adult** patient mode only.

In MPV breath mode it is possible to set Breath Rate to zero so that breaths are only initiated when the patient triggers them using the mouthpiece. The PEEP setting is not available but always set to Off so that no air blows from the mouthpiece when no breath is being delivered.

As the patient is not always connected to the ventilator several of the surveillance functions are not working in the same way as in other modes:

- Disconnection alarm is not available in MPV mode.
- Low Pressure alarm is only active during breaths.
- Apnoea alarm is possible to set to longer times and is an important mean of surveillance to make sure that the patient gets ventilation support regularly.

When switching to MPV mode all alarms, except from High/Low Pressure Alarms, are automatically switched off to avoid false alarm triggering (when changing between profiles the alarm settings does not change to Off, but stays as defined in the profiles).

#### **WARNING!**

- The alarm levels must be assessed and adjusted considering the patient's condition and treatment settings.
- MPV breath mode shall be used with the mouthpiece circuit including the mouthpiece interface only.
- The MPV circuit including the mouthpiece interface shall be used in MPV breath mode only and is not for use in any other breath mode.
- MPV breath mode shall not be used with ventilator dependent patients.



#### 5.4.4.13 CPAP – Continuous Positive Airway Pressure

In CPAP mode the Vivo 45 LS is applying a continuous positive pressure to the airways. The flow will automatically be adjusted to maintain the set CPAP level.

5.5 Functions and Parameters in the Vivo 45 LS

The parameters used for controlling the breathing by the Vivo 45 LS are described below. See also 8.2.3 *Parameter Specifications*, page 131 for detailed parameter data.

Depending on the setup for Home mode adjustment, parameters may not be available or have limited ranges when the Vivo 45 LS is in Home mode.

5.5.1 Parameters' Availability in Different Modes

This table offers an overview of the available parameters in each mode.

| Modes<br>Parameter  | PSV | PSV<br>(TgV) | PCV | PCV<br>(TgV) | PCV (A) | PCV<br>(A<br>+TgV) | VCV(A) | PCV<br>SIMV | VCV-<br>SIMV | PCV-<br>MPV | VCV-<br>MPV | CPAP |
|---|-----|--------------|-----|--------------|---------|--------------------|--------|-------------|--------------|-------------|-------------|------|
| Insp. Pressure  | x   |              | x   |              | x       |                    |        | x           |              | x           |             |      |
| PEEP  | x   | x            | x   | x            | x       |                    | x      | x           |              |             |             |      |
| Auto-EPAP;<br>MIN EPAP<br>MAX EPAP<br>Pressure limit<br>EPAP step<br>Relax Time | x   | x            | x   | x            | x       |                    |        |             |              |             |             |      |
| Breath Rate   |     |              | x   | x            | x       |                    | x      |             |              |             |             |      |
| Insp. Time  |     |              | x   | x            | x       |                    | x      | x           |              | x           |             |      |
| Backup Insp.<br>Time  | x   | x            |     |              |         |                    |        |             |              |             |             |      |
| Sigh  | x   | x            | x   | x            | x       |                    | x      |             |              |             |             |      |
| Rise Time   | x   | x            | x   | x            | x       |                    | x      | x           |              | x           |             |      |
| Insp. Trigger   | x   | x            |     |              | x       |                    | x      | x           |              | x           |             |      |
| Exp. Trigger  | x   | x            |     |              |         |                    |        | x           |              |             |             |      |
| Min Insp.<br>Time   | x   | x            |     |              |         |                    |        |             |              |             |             |      |
| Max Insp.<br>Time   | x   | x            |     |              |         |                    |        |             |              |             |             |      |
| Backup Rate   | x   | x            |     |              |         |                    | x      |             |              | x           |             |      |
| Target Volume:<br>Max Pressure<br>Min Pressure                                  |     | x            |     | x            | x       |                    |        |             |              |             |             |      |
| Tidal Volume  |     |              |     |              |         |                    | x      | x           |              |             |             |      |
| Flow pattern  |     |              |     |              |         |                    | x      |             |              |             |             |      |
| SIMV rate   |     |              |     |              |         |                    |        | x           |              |             |             |      |

| Modes<br>Parameter       | PSV | PSV<br>(TgV) | PCV | PCV<br>(TgV) | PCV (A) | PCV<br>(A<br>+TgV) | VCV(A) | PCV<br>SIMV | VCV-<br>SIMV | PCV-<br>MPV | VCV-<br>MPV | CPAP |
|--------------------------|-----|--------------|-----|--------------|---------|--------------------|--------|-------------|--------------|-------------|-------------|------|
| SIMV support<br>pressure |     |              |     |              |         |                    |        | x           | x            |             |             |      |
| CPAP                     |     |              |     |              |         |                    |        |             |              |             |             | x    |

### 5.5.2 Insp. Pressure

The Insp. Pressure setting is used to define the airway pressure during the inspiratory phase. Minimum/maximum working IPAP is limited/achieved by a software control of blower speed vs. measured pressure.

| Unit               | Min | Max | Default |
|--------------------|-----|-----|---------|
| cmH <sub>2</sub> O | 4   | 50  | 15      |

### 5.5.3 PEEP

The PEEP setting is used to define the airway pressure at the end of the expiratory phase.

| Unit               | Min                     | Max               | Default                              |
|--------------------|-------------------------|-------------------|--------------------------------------|
| cmH <sub>2</sub> O | 2<br>Off <sup>(1)</sup> | 20 <sup>(2)</sup> | 5 <sup>(A)</sup><br>2 <sup>(P)</sup> |

(1)= Off only for patient circuits with active exhalation valve

(2)= For pressure ventilation modes: The max setting is also limited by Insp. Pressure - 2 cmH<sub>2</sub>O and Min Pressure -2 cmH<sub>2</sub>O.

(A)= Adult mode, (P)= Paediatric mode

### 5.5.4 Breath Rate

The Breath Rate setting defines the minimum number of breaths the ventilator will deliver as long as no inspiratory trigger effort from the patient is detected. The cycles will be ventilator-initiated breaths.

The combination of the Breath Rate and Inspiratory Time setting is limited by the I:E ratio range 1:9.9 to 2:1.

| Unit                                  | Min                                  | Max                                    | Default                                |
|---------------------------------------|--------------------------------------|--|--|
| 1/breath                              | 4 <sup>(A)</sup><br>6 <sup>(P)</sup> | 40 <sup>(A)</sup><br>60 <sup>(P)</sup> | 12 <sup>(A)</sup><br>20 <sup>(P)</sup> |
| (A)= Adult mode, (P)= Paediatric mode |                                      |  |  |

### 5.5.5 Backup Rate

The Backup Rate setting defines the minimum number of breaths the ventilator will deliver as long as no inspiratory trigger effort from the patient is detected. The cycles will be ventilator-initiated breaths.

The combination of the Breath Rate and Inspiratory Time setting is limited by the I:E ratio range 1:9.9 to 2:1.

| Unit                                  | Min   | Max  | Default   |
|---------------------------------------|---|--|---|
| breath/<br>minute                     | 4 <sup>(A)</sup><br>6 <sup>(P)</sup><br>0 (MPV) | 40 <sup>(A)</sup><br>60 <sup>(P)</sup><br>40 (MPV) | 12 <sup>(A)</sup><br>20 <sup>(P)</sup><br>0 (MPV) |
| (A)= Adult mode, (P)= Paediatric mode |   |  |   |

### 5.5.6 SIMV Rate

The SIMV Rate setting is used in the SIMV ventilation modes, for defining the minimum frequency of mandatory, ventilator-controlled breaths. The mandatory breaths can be either triggered by an inspiratory effort from the patient, or ventilator-initiated. The SIMV Rate setting determines the SIMV cycle time.



The combination of the SIMV Rate and Inspiratory Time setting is limited by the I:E ratio 2:1.

| Unit                                  | Min                                  | Max                                    | Default                                |
|---------------------------------------|--------------------------------------|--|--|
| 1/breath                              | 4 <sup>(A)</sup><br>6 <sup>(P)</sup> | 40 <sup>(A)</sup><br>60 <sup>(P)</sup> | 12 <sup>(A)</sup><br>20 <sup>(P)</sup> |
| (A)= Adult mode, (P)= Paediatric mode |                                      |  |  |

### 5.5.7 Insp. Time (Inspiratory Time)

The Inspiratory Time setting defines the length of each inspiration from start of inspiration to cycling off to expiration.

The combination of the Inspiratory Time and Breath Rate or SIMV Rate settings is limited by the I:E ratio 2:1.

| Unit                                  | Min | Max                                  | Default                                |
|---------------------------------------|-----|--------------------------------------|--|
| s                                     | 0.3 | 5 <sup>(A)</sup><br>2 <sup>(P)</sup> | 1.5 <sup>(A)</sup><br>1 <sup>(P)</sup> |
| (A)= Adult mode, (P)= Paediatric mode |     |                                      |  |

### 5.5.8 Backup Insp. Time (Backup Inspiratory Time)

The Backup Inspiratory Time setting defines the length of each inspiration delivered during ventilator-triggered backup ventilation, initiated by the set Backup Rate.

The combination of the Backup Inspiratory Time and Backup Rate setting is limited by the I:E ratio 2:1.

| Unit                                  | Min | Max                                  | Default                                |
|---------------------------------------|-----|--------------------------------------|--|
| s                                     | 0.3 | 5 <sup>(A)</sup><br>2 <sup>(P)</sup> | 1.5 <sup>(A)</sup><br>1 <sup>(P)</sup> |
| (A)= Adult mode, (P)= Paediatric mode |     |                                      |  |

### 5.5.9 Sigh Parameters

With the Sigh function, the ventilator will periodically deliver extended breaths.

#### NOTE

In pressure modes (during the sigh breath), the high pressure alarm will automatically be set 10 cmH<sub>2</sub>O above set sigh pressure .

When the sigh function is on, the following parameters are available:

| Sigh Rate  |     |     |         |
|--|-----|-----|---------|
| The Sigh rate sets the frequency of which breaths with an increased pressure or volume are delivered to the patient.<br>If the High Pressure alarm or the High Tidal Volume alarm is given, the Sigh function will be disabled as long as the alarm condition remains. |     |     |         |
| Unit   | Min | Max | Default |
| 1/ breath  | 10  | 250 | 50      |

| Sigh %   |     |     |         |
|--|-----|-----|---------|
| Sigh % sets the increased % of the set pressure is delivered to the patient. |     |     |         |
| Unit   | Min | Max | Default |
| %  | 125 | 200 | 125     |

| Sigh Inspiratory Time  |   |                                      |  |
|--|---|--------------------------------------|--|
| Sigh inspiratory time sets the inspiratory time during sigh breaths. |   |                                      |  |
| Unit   | Min   | Max                                  | Default                                |
| s  | Current <i>Inspiratory Time</i><br>or<br><i>Backup Inspiratory Time</i> | 5 <sup>(A)</sup><br>2 <sup>(P)</sup> | 1.5 <sup>(A)</sup><br>1 <sup>(P)</sup> |
| (A)= Adult mode, (P)= Paediatric mode                                |   |                                      |  |

### 5.5.10 Rise Time

The Rise Time setting controls the speed of the pressure/volume increase from start of inspiration to the set pressure or volume.

A low setting will give a faster increase and therefore a longer plateau at the set value. A high setting will give a slower increase and therefore a shorter plateau.

| Rise Time in Pressure mode |                       |                       |                       |
|----------------------------|-----------------------|-----------------------|-----------------------|
| Unit                       | Min                   | Max                   | Default               |
| Step                       | 1 (PSV and PCV modes) | 9 (PSV and PCV modes) | 3 (PSV and PCV modes) |

| Rise Time in Volume mode   |                               |                                      |         |
|--|-------------------------------|--------------------------------------|---------|
| Unit   | Min                           | Max                                  | Default |
| In the Volume Control mode (VCV), the rise time for the support is set as a percentage of the inspiratory time.  |                               |                                      |         |
| %  | 50% of set <i>Insp Time</i> . | 90% of set <i>Insp Time</i> .<br>Off | Off     |
| In the VCV-SIMV mode, the rise time set as percentage if the inspiratory time is for the mandatory breaths. The rise time for the support breath triggered by the patient is set in steps. |                               |                                      |         |
| Step   | 1                             | 9                                    | 3       |

### 5.5.11 Insp. Trigger (Inspiratory Trigger)

The inspiratory trigger defines the patient's effort required to initiate a ventilator assisted breath. When the patient starts a breath, an increasing flow is created in the patient circuit. If the patient's effort reaches the set inspiratory trigger level an inspiration is initiated.

If the patient cannot trigger a breath, the ventilator will deliver breaths according to the set Backup Rate, Breath Rate or SIMV Rate.

| Unit   | Min              | Max                     | Default                              |
|--|------------------|-------------------------|--------------------------------------|
| Step   | 1 <sup>(1)</sup> | 9<br>Off <sup>(2)</sup> | 3 <sup>(A)</sup><br>2 <sup>(P)</sup> |
| (A)= Adult mode, (P)= Paediatric mode<br>(1)= Low value is easy to trigger, high value is harder to trigger.<br>(2)= Off turns off the assisted breath function (Control mode only). |                  |                         |                                      |

### 5.5.12 Sup. Pressure (SIMV mode)

The Support Pressure setting is used in the SIMV ventilation modes, for defining the inspiratory pressure for the support breaths triggered by the patient.

| Unit               | Min            | Max | Default |
|--------------------|----------------|-----|---------|
| cmH <sub>2</sub> O | Current PEEP+2 | 50  | 15      |

### 5.5.13 Exp. Trigger (Expiratory Trigger)

The Expiratory Trigger setting defines the moment when the ventilator will cycle from the inspiratory to the expiratory phase.

| Unit  | Min              | Max              | Default |
|---|------------------|------------------|---------|
| Step  | 1 <sup>(1)</sup> | 9 <sup>(1)</sup> | 3       |
| (1)= Low value is easy to trigger, high value is harder to trigger. |                  |                  |         |

### 5.5.14 Max Insp. Time (Maximum Inspiratory Time)

The Maximum Inspiratory Time setting defines a maximum length for each inspiration. If the Maximum Inspiratory Time is set to Off, the length of the inspiration and/or minimum inspiratory time is dependent on the set Expiratory Trigger.

| Unit                                  | Min | Max   | Default |
|---------------------------------------|-----|---|---------|
| s                                     | 0.3 | 5 <sup>(A)</sup><br>2 <sup>(P)</sup><br>Off | Off     |
| (A)= Adult mode, (P)= Paediatric mode |     |   |         |

### 5.5.15 Min Insp. Time (Minimum Inspiratory Time)

The Minimum Inspiratory Time setting defines a minimum length for each inspiration. If the Minimum Inspiratory Time is set to Off, the length of the inspiration is dependent on the set Expiratory Trigger.

| Unit                                  | Min        | Max                                  | Default |
|---------------------------------------|------------|--------------------------------------|---------|
| s                                     | 0.3<br>Off | 3 <sup>(A)</sup><br>2 <sup>(P)</sup> | Off     |
| (A)= Adult mode, (P)= Paediatric mode |            |                                      |         |

**NOTE**

If Target Volume is used with a patient circuit with an active exhalation valve, leakage may be misinterpreted by the ventilator as an increase of tidal volume. This will lead to a decrease of the Inspiratory Pressure (the Inspiratory Pressure will not be lower than the set Min Pressure). This may result in hypoventilation as the true delivered tidal volume will decrease both as a result of the leakage and the decrease in Inspiratory Pressure. This does not occur if a patient circuit with leakage port is used.

The Target Volume setting defines the tidal volume that the ventilator will aim for while ventilating the patient in a pressure mode. To aim for the preset volume, the ventilator will adapt the Inspiratory Pressure between two adjustable pressure limits: Min Pressure and Max Pressure.

When Target Volume is active, the mode field on the ventilator display will indicate "(IgV)".

| Unit | Min  | Max                                       | Default |
|------|--|---|---------|
| ml   | Off<br>300 <sup>(A)</sup><br>50 <sup>(P)</sup> | 2000 <sup>(A)</sup><br>500 <sup>(P)</sup> | Off     |

(A)= Adult mode, (P)= Paediatric mode

**Target Volume Parameters**

When target volume is on, the following parameters are enabled:

| Max Pressure   |                             |                             |         |
|--|-----------------------------|-----------------------------|---------|
| Max Pressure defines the upper pressure limit up to where the ventilator can increase the pressure to reach the set Target Volume. If Target Volume is not reached at Max Pressure, the ventilator will continue to ventilate at this Max Pressure setting.                  |                             |                             |         |
| Unit   | Min                         | Max                         | Default |
| cmH <sub>2</sub> O   | Current <i>Min Pressure</i> | 50                          | 15      |
| Min Pressure   |                             |                             |         |
| Min Pressure defines the lower pressure limit down to where the ventilator can decrease the pressure to maintain the set Target Volume. If the actual volume is above Target Volume at Min Pressure, the ventilator will continue to ventilate at this Min Pressure setting. |                             |                             |         |
| Unit   | Min                         | Max                         | Default |
| cmH <sub>2</sub> O   | 4                           | Current <i>Max Pressure</i> | 15      |

**Tidal Volume**

The Tidal Volume setting defines the volume that will be delivered by the Vivo 45 LS for each breath when using volume control modes.

In VCV-SIMV mode, this setting is applicable to the mandatory breaths that are controlled by the ventilator.

| Unit | Min                                     | Max                                       | Default                                  |
|------|---|---|--|
| ml   | 300 <sup>(A)</sup><br>50 <sup>(P)</sup> | 2000 <sup>(A)</sup><br>500 <sup>(P)</sup> | 500 <sup>(A)</sup><br>150 <sup>(P)</sup> |

(A)= Adult mode, (P)= Paediatric mode

### 5.5.18 Flow Pattern

Flow pattern sets the characteristics of the air flow in VCV modes.

| Unit | Min   | Max | Default |
|------|---|-----|---------|
| —    | Square (constant flow)<br>Decelerating (Flow is decreased linearly, may prevent air hunger) |     | Square  |

### 5.5.19 CPAP

The CPAP setting defines the pressure that will be applied to the airways in CPAP mode.

| Unit                                  | Min | Max | Default                               |
|---------------------------------------|-----|-----|---------------------------------------|
| cmH <sub>2</sub> O                    | 4   | 20  | 10 <sup>(A)</sup><br>8 <sup>(P)</sup> |
| (A)= Adult mode, (P)= Paediatric mode |     |     |                                       |

### 5.5.20 Humidifier

Humidifier allows the user to start or stop the heated humidification.

The click-in water chamber needs to be connected before the setting can be turned On.

| Unit | Min | Max | Default |
|------|-----|-----|---------|
| —    | Off | On  | Off     |

### 5.5.21 Humidifier Setting

The Humidifier Setting defines the level of humidity of the air delivered to the patient.

| Unit | Min | Max | Default |
|------|-----|-----|---------|
| Step | 1   | 5   | 3       |

### 5.5.22 Heated Circuit Temp

Heated Circuit Temp setting will define the temperature of the heated circuit.

| Unit  | Min   | Max   | Default |
|-------|-------|-------|---------|
| °C/°F | 16/61 | 30/86 | 27/81   |

### 5.5.23 Circuit Heating

Circuit Heating allows the user to start or stop the heating of the circuit.

The Heated Circuit needs to be connected before the setting can be turned On.

| Unit | Min | Max | Default |
|------|-----|-----|---------|
| —    | Off | On  | Off     |

## 5.6 Using Batteries

Since all batteries, in general, degenerate over time, the recommendations below will ensure that the battery capacity of the ventilator is maximized during its lifetime.

The internal and click-in batteries in the ventilator are of the Lithium-ion type, which is a high performance battery. It has long expected lifetime, low weight in relation to its capacity and low self discharge.

### **5.6.1 Power Source Priority**

1. AC power (Mains)
2. External DC
3. Click-in battery
4. Internal battery

When a power source fails, the ventilator will switch to either the external DC (if installed), the click-in battery (if attached) or the internal battery and show a message in the display window.

### **5.6.2 Power Switchover Tests**

This section describes methods for testing the ventilator's switch over from one power source to another.

## Internal Battery Switchover Test

### Test Conditions

- AC power connected
  - No click-in battery installed
  - Treatment is running
- 1 Disconnect the AC power cord.
  - 2 Check that the:
    - Internal battery power source LED is illuminated
    - Medium priority “Lost Mains Power” alarm is triggered
    - Information message “Switched to Internal Battery” is posted

## Click-in Battery Switchover Test

### Test Conditions

- AC power connected
  - Click-in battery installed
  - Treatment is running
- 1 Disconnect the AC power cord.
  - 2 Check that the:
    - Click-in battery power source LED is illuminated
    - Medium priority “Lost Mains Power” alarm is triggered
    - Information message “Switched to Click-In Battery” is posted

## External DC Switchover Test

### Test Conditions

- AC power connected
  - External DC connected
  - Treatment is running
- 1 Disconnect the AC power cord.
  - 2 Check that the:
    - External DC power source LED is illuminated
    - Medium priority “Lost Mains Power” alarm is triggered
    - Information message “Switched to External DC” is posted

### 5.6.3

## Charging the Batteries

### CAUTION!

Do not charge the Vivo 45 LS while placed in the carry bag or other types of closed or non-ventilated spaces.



The internal and click-in batteries are automatically charged when connecting the ventilator to the mains supply. To ensure that the batteries are fully charged, a maintaining charging cycle will be performed. The batteries are not charged when connecting the ventilator to an external DC supply. While charging, the battery level will be animated. The batteries are only charged if the internal temperatures are between 0 to 45°C (32 to 113°F). High power

consuming settings in combination with high ambient temperatures may make the battery temperature rise above 45°C (113°F).



## NOTE

Charging of batteries is only started when the state of charge is below 95%

## Behaviour of the ventilator while Internal Battery or Click-in Battery is Charging

The ventilator can operate as normal while a battery is charging. The battery icon will be animated (filling in from bottom to top) while charging.

## Charging Times

| Battery          | Charger    | Time |
|------------------|------------|------|
| Internal battery | Vivo 45 LS | 2 h  |
| Click-in battery | Vivo 45 LS | 4 h  |

\* Times are based on charging empty batteries.

## 5.6.4

## Battery Icons

When running on battery, the battery status is indicated by the following symbols:

| Symbols | Battery Status   |
|---------|--|
|         | High state of charge (over 50 %)   |
|         | Medium State of charge (20 % – 50 %)   |
|         | Low state of charge (below 20 %)   |
|         | Malfunctioning battery   |
|         | Information box displaying state of charge (alters with battery symbol above, when in standby mode)                      |
|         | Information box displaying estimated remaining treatment time (alters with battery symbol above, when in treatment mode) |

## 5.6.5

## Internal Battery

The internal battery is intended as a backup power source if the primary power source fails. It can also be used as a temporary power source. For example during transportation between one stationary power source to another.



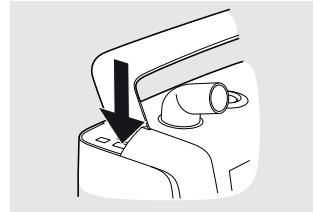
## 5.6.6 Click-in Battery

The click-in battery is intended as a power source during transportation, or if the primary AC power source fails.


The click-in battery can be replaced during treatment, provided that the internal battery is charged.

### Connect the Click-in Battery


1 Release the side cover by pressing the button under the handle.




2 Open and remove the side panel.

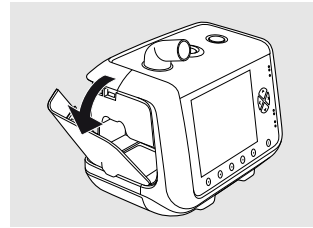

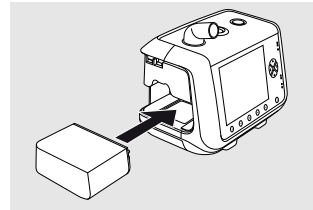3 Insert the click-in battery. Make sure the latch at the bottom of the click-in compartment is engaged.



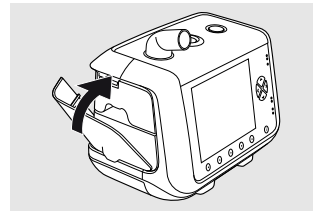4 Close the side panel. Make sure there is a clicking sound to secure the side panel.



When removing the battery, press down the latch at the bottom of the battery compartment and tilt the ventilator sideways. Make sure to close the side panel after removing the click-in battery.

### 5.6.7 Battery Operating Time (Internal and Click-in)

The operation time is dependent on the battery condition, its capacity, the ambient air temperature and the ventilator pressure setting. These data are based on new and fully charged batteries.

| Condition                       | Value                     |
|---------------------------------|---------------------------|
| <b>Environmental Conditions</b> |                           |
| Ambient temperature             | 20°C (68°F)               |
| <b>Ventilator settings</b>      |                           |
| Mode                            | PCV                       |
| Pressure*                       | 20 cmH <sub>2</sub> O     |
| PEEP                            | 4 cmH <sub>2</sub> O      |
| Breath Rate*                    | 20 bpm                    |
| Insp. Time*                     | 1.0 s                     |
| I:E                             | 1:2                       |
| Insp. Trigger                   | Off                       |
| Rise Time                       | 1                         |
| Target Volume                   | Off                       |
| Display Light*                  | Off                       |
| Light Intensity*                | -                         |
| <b>Other</b>                    |                           |
| Tidal Volume                    | 800 ml                    |
| Resistance                      | 5 hPa (l/s) <sup>-1</sup> |
| Compliance                      | 50 ml (hPa) <sup>-1</sup> |

\*: These ventilator settings affect the operating time significantly.

| Battery          | Operating Time |
|------------------|----------------|
| Internal Battery | 2.5 h          |
| Click-in Battery | 6.5 h          |

### 5.6.8 Storing the Internal Battery and the Click-in Battery

Storage longer than 1 month should be initiated with half-charged batteries in order to maintain maximum capacity. Optimal storage temperature is 5 to 30°C (41 to 86°F).

#### 5.6.9 External DC **WARNING!**



Do not connect the Vivo 45 LS to a wheelchair unless the operating manual for the wheelchair permits this as this can affect the Vivo 45 LS performance and consequently result in patient death.

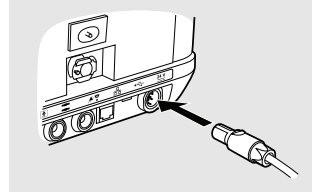
The ventilator can be operated from:

- The Breas XPAC using the XPAC battery cable.
- A 12 V external DC source using the 12-24 VDC car adapter cable.
- A 24 V external DC source using the external battery cable.
- Both AC power supply and external DC using the Y-cable.

With both AC power supply and external DC source connected, the ventilator will automatically switch over to the external DC source if the AC power supply fails. The external DC voltage level is shown under “Others”, “Device Information” in the menu

**1** Connect the external DC cable to the ventilator.

Make sure that it is fitted correctly.




**2** Connect the other end of the cable to the DC source.


## 5.7 Using Accessories
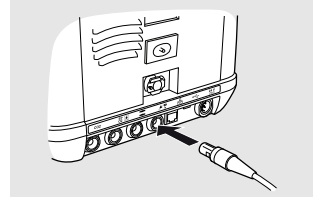
### 5.7.1 Connecting and Disconnecting the Cables


- 1 Insert the cable in the appropriate port.



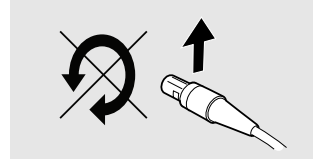- 2 Make sure to insert the connector with the marking pointing upwards.



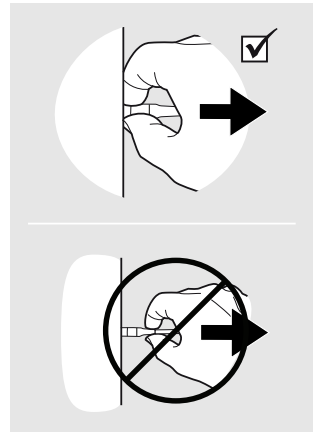- 3 Pull the connector sleeve, not the cable itself or cable restrainer to release the connector.




### 5.7.2 Using the ventilator with a Nurse Call System

The ventilator can be connected to a nurse call system using the nurse call cable. When connected, the ventilator alarms will also be forwarded to the nurse call system.

5.7.2.1 Connect the ventilator to a Nurse Call System

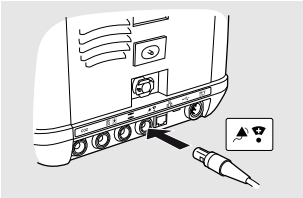1 Connect the nurse call cable at the back of the ventilator.



2 Test the connection by triggering an alarm on the ventilator and verify that the nurse call system activates.

5.7.3 Using the ventilator with the FiO<sub>2</sub> Sensor

The FiO<sub>2</sub> sensor can be used to monitor and store FiO<sub>2</sub> measurements. The FiO<sub>2</sub> sensor measures the fraction of inspired oxygen in the air channel of the ventilator. The FiO<sub>2</sub> measurements will be stored in the data memory which can be downloaded to a PC and viewed in Breas software.

| Usage                   | Time   |
|-------------------------|--|
| Operating temperature   | 10 to 40°C (50 to 104°F)   |
| Operating pressure      | 700 to 1250 mbar   |
| Expected operating life | <3 years (in ambient air) or 500,000 Vol.% h, whichever comes first. |
| Shelf life              | < 6 months (recommended)   |


**CAUTION!**

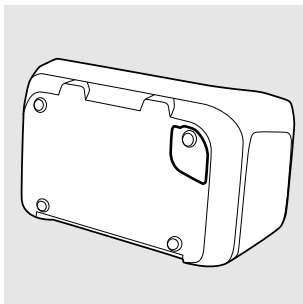Note that the operating conditions for the FiO<sub>2</sub> sensor is different from the ventilator system conditions. If the sensor is used outside its operating conditions the FiO<sub>2</sub> measurements might deviate.




### 5.7.3.1 Installing the FiO<sub>2</sub> Sensor

- 1 Place the ventilator so the bottom is accessible.
- 2 Remove the hatch for FiO<sub>2</sub> sensor. Use a torx TX10 screwdriver.



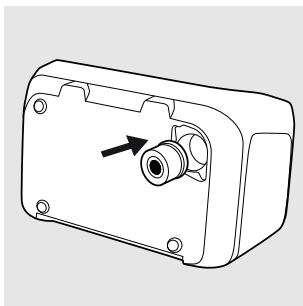- 3 Insert the FiO<sub>2</sub> sensor with the electrical contact side in.



- 4 Reinstall the hatch.
- 5 Calibrate the FiO<sub>2</sub> sensor in the advanced settings of the main menu.



When installed, the ventilator automatically detects the sensor, also after powering off/on and after power failure.

### 5.7.3.2 Calibrating the FiO<sub>2</sub> sensor

The FiO<sub>2</sub> sensor should be calibrated when first used and then at least once a month.



FiO<sub>2</sub> calibration can be performed from the “FiO<sub>2</sub>/CO<sub>2</sub> Calibration” page under the Advanced settings section of the main menu.

### 5.7.4 Using the ventilator with the Remote Alarm



Information about safety, warnings, product description, installation, usage, cleaning, maintenance and technical specifications can be found in the user instruction for Remote Alarm.

The Remote Alarm enables care providers and clinical personnel to monitor the ventilator alarms remotely. The Remote Alarm forwards alarms from the ventilator. When an alarm sounds, the care provider or clinical personnel must attend to the patient quickly.

When installing a remote alarm system, check that it operates as intended before starting the treatment.

## 5.7.5 Using the ventilator with the EtCO<sub>2</sub> Sensor

The EtCO<sub>2</sub> sensor can be connected to the patient breathing circuit and to a Vivo 45 LS in order to monitor and store CO<sub>2</sub> measurements. The CO<sub>2</sub> measurements will be stored in the ventilator data memory which can be downloaded to a PC and viewed in the ventilator PC software.



More information about safety, warnings, product description, installation, usage, cleaning, maintenance and technical specifications can be found in the user instruction for the EtCO<sub>2</sub> sensor.

### 5.7.5.1 Safety Information

#### **WARNING!**

Read the instructions thoroughly so that you completely understand how the EtCO<sub>2</sub> sensor is operated before taking it into use, to ensure correct usage and maximum performance.

Breas Medical reserves the right to make changes to this product without any prior notification.

Do not use a damaged CO<sub>2</sub> sensor or adapter.

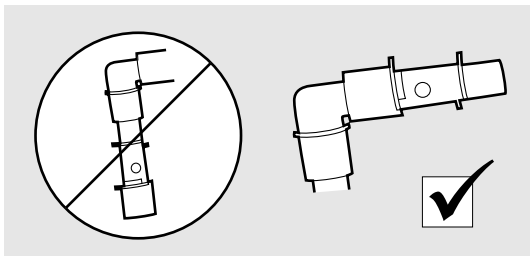
The CO<sub>2</sub> sensor is intended to be used by authorized and trained medical personnel only.

The CO<sub>2</sub> sensor is intended only as an adjunct in patient assessment. It must be used in conjunction with other assessments of clinical signs and symptoms.

Masks' dead-space, patient's volumes and unintentional leakage may influence the CO<sub>2</sub> measurements.

Used airway adapters shall be disposed of in accordance with local regulations for medical waste.

Measurements can be affected by mobile and RF communications equipment. It should be assured that the CO<sub>2</sub> sensor is used in the electromagnetic environment specified in 8.3 *Emission and Immunity Declaration*, page 138.



Do not place the airway adapter between the endotracheal tube and an elbow as this may allow patient secretions to block the adapter windows and result in incorrect operation.

Incorrect CO<sub>2</sub> zeroing will result in false gas readings.

Replace the airway adapter if rain-out/condensation occurs inside the airway adapter.

Only use airway adapters distributed by Breas Medical.

Do not apply tension to the CO<sub>2</sub> sensor cable.



To keep secretions and moisture from pooling on the windows, always position the CO<sub>2</sub> sensor in a vertical position with the green LED pointing upwards.



#### **WARNING!**

Disposable airway adapters shall not be reused. Reuse of the single use adapter can cause cross infection.



#### **CAUTION!**

If an intentional leakage port is used, make sure that the CO<sub>2</sub> sensor is placed between the patient interface and the leakage port.

If a patient interface with integrated leakage is used, the monitored CO<sub>2</sub> values may be influenced.



The CO<sub>2</sub> sensor should be placed as close to the patient interface as possible. However, a HME (if used) shall be placed between the patient interface and the CO<sub>2</sub> sensor. This will protect the airway adapter from secretions and effects of water vapour and eliminates the need of changing the airway adapter.



#### **NOTE**

The CO<sub>2</sub> monitoring automatically compensates for changes in ambient barometric pressure. The CO<sub>2</sub> monitor shall fulfil the ISO 80601-2-55 standard (Medical electrical equipment – Particular requirements for the basic safety and essential performance of respiratory gas monitors).

### **5.7.5.2 How to Connect the EtCO<sub>2</sub> Sensor**



#### **WARNING!**

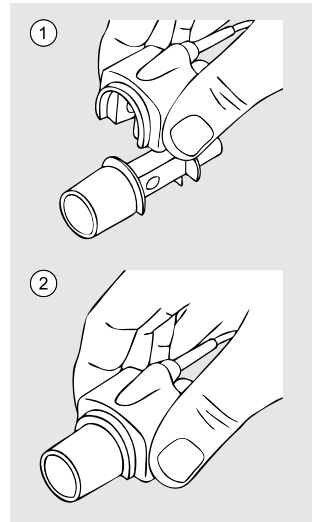
The CO<sub>2</sub> sensor is not intended to be in contact with the patient body.



1 Connect the CO<sub>2</sub> sensor cable to the CO<sub>2</sub> connection port on the ventilator (according to the instruction 5.7.1 *Connecting and Disconnecting the Cables*, page 74).


A green LED indicates that the CO<sub>2</sub> sensor is ready to use.

2 Snap the CO<sub>2</sub> sensor probe on top of the airway adapter. It will click into place when properly sealed.




3 Perform a CO<sub>2</sub> zeroing procedure.

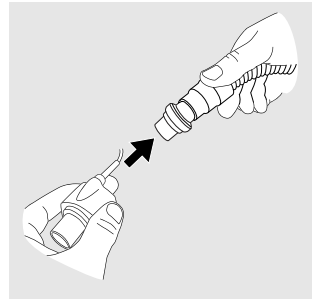4 Connect the airway adapter to the patient circuit.



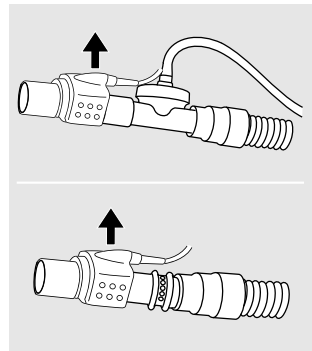5 Make sure to position the CO<sub>2</sub> sensor with the LED pointing upwards.



When installed, the ventilator automatically detects the sensor, also after powering off/on or after power failure.

## CO<sub>2</sub> Zeroing

CO<sub>2</sub> zeroing is recommended when changing the airway adapter. Besides from that, zeroing only needs to be performed when an offset in monitored CO<sub>2</sub> values is observed, or when a CO<sub>2</sub> sensor accuracy unspecified message is displayed.

| LED Status           | Description         |
|----------------------|---------------------|
| Steady green light   | System OK           |
| Flashing green light | Zeroing in progress |
| Steady red light     | Sensor error        |
| Flashing red light   | Check adapter       |

## Maintenance

No periodical maintenance is required for the CO<sub>2</sub> sensor.

To verify the CO<sub>2</sub> sensor readings, a gas span check shall be performed every year, preferably when the ventilator is sent for service.



See the ventilator service manual for how to perform the gas span check.



### WARNING!

Do not under any circumstances attempt to service or repair the CO<sub>2</sub> sensor yourself. If you do so, the manufacturer will no longer be responsible for the performance and safety of the CO<sub>2</sub> sensor.



## Cleaning

### WARNING!

- Always be careful when cleaning to ensure that you do not damage any equipment.
- Fluid must not be allowed to enter into the CO<sub>2</sub> sensor.
- Always clean the T-piece with plug when to be used by a new patient. All parts that come into contact with the respiration gas must be cleaned.
- Remove the airway adapter before cleaning.
- Do not sterilise the CO<sub>2</sub> sensor.
- Do not autoclave the CO<sub>2</sub> sensor.



Clean the outside of the CO<sub>2</sub> sensor using a lint-free cloth moistened, but not wet, with ethanol or isopropyl alcohol (< 70%).

## Disposal

The CO<sub>2</sub> sensor must be disposed of in accordance with the local environmental regulations regarding the disposal of used equipment and waste.

## 5.7.6 Using the ventilator with the PtcCO<sub>2</sub> Cable

An external monitor for Transcutaneous CO<sub>2</sub> Pressure (PtcCO<sub>2</sub>) may be connected to the ventilator by an accessory PtcCO<sub>2</sub> cable. For information about available PtcCO<sub>2</sub> cables, see 9 *Accessories and Parts*, page 144.



### NOTE

Both the PtcCO<sub>2</sub> cable and the EtCO<sub>2</sub> sensor connects at the ventilator's yellow CO<sub>2</sub> port. Only one CO<sub>2</sub> measuring device can be connected at a time.

When connected, the ventilator will:

- Display the monitored values and include them in trend views.
- Store monitored values in the internal memory. The PtcCO<sub>2</sub> values will also be included in the data that can be downloaded and analysed with Breas PC software.
- Repeat CO<sub>2</sub> alarms from the external PtcCO<sub>2</sub> monitor.
- Automatically detect the sensor, also after powering off/on or after power failure

## 5.7.7 Using the Ventilator with the SpO<sub>2</sub> module



Information about safety, warnings, product description, installation, usage, cleaning, maintenance and technical specifications can be found in the user instruction for SpO<sub>2</sub> module.

The SpO<sub>2</sub> module enables connection to an SpO<sub>2</sub> sensor for measuring of functional oxygen saturation of arterial haemoglobin (SpO<sub>2</sub>) and pulse rate. The SpO<sub>2</sub> module can be connected to the Vivo 45 LS in order to monitor and store SpO<sub>2</sub> measurements.

The SpO<sub>2</sub> measurements will be stored in the data memory which can be downloaded to a PC and viewed in the Breas PC software.

When installed, the Vivo 45 LS automatically detects the sensor, also after powering off/on and after power failure.

## 5.7.8 Using the ventilator with the Remote Start/Stop



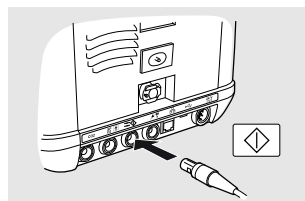
### NOTE

The effort belt communication box and the Remote Start/Stop use the same port on the ventilator. Only one of the accessories can be connected at a time.

### 5.7.8.1 Connecting the Remote Start/Stop

1 Connect the Remote Start/Stop cable to the ventilator.

Information about safety, warnings, product description, installation, usage, cleaning, maintenance and technical specifications can be found in the user instruction for Remote Start/Stop.




## 5.7.9 Using the ventilator with the Protective Cover

The protective cover is intended for additional protection of the ventilator during transportation, and in hospital, institutional and home care environments. It can be used while the

ventilator is operating, for example mounted on a wheelchair, in a personal vehicle, or carried by hand.

The protective cover protects the ventilator from environmental impact such as shock, water spill, sunlight, dust and dirt, under normal handling.

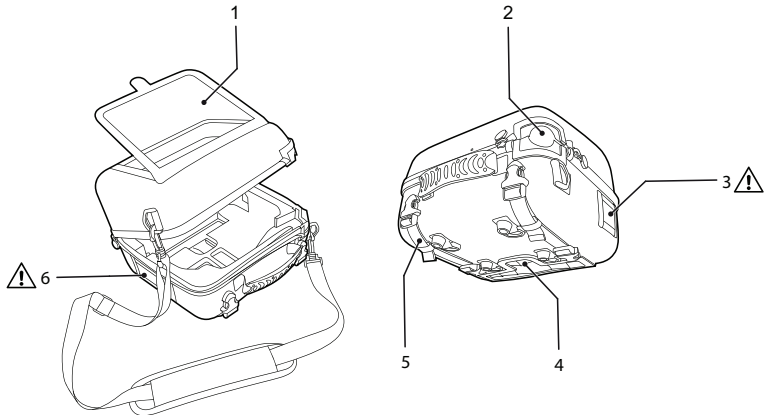The protective cover has the following functions:



1. Transparent window, for accessing front panel and buttons
2. Port for patient circuit
3. Cooling air inlet
4. Port for cables and O<sub>2</sub> inlet
5. Mounting straps
6. Patient air inlet

#### **CAUTION!**

Do not cover the air inlets or outlets.



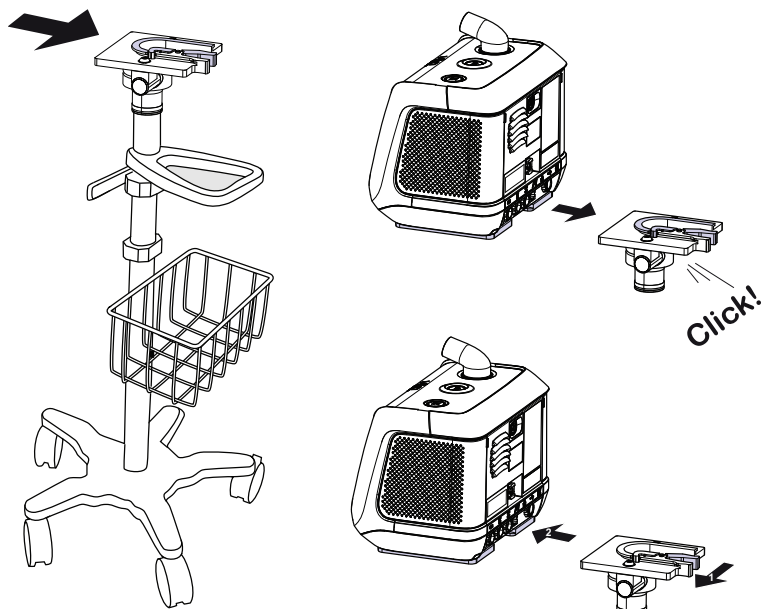
### **5.7.10 Using the Vivo 45 LS with the Trolley**

#### **Intended Use**

The intended use of the trolley system is to allow patient mobility while receiving ventilator treatment. The trolley shall only be used in indoor, hospital environment. The trolley system consists of a trolley base and a mounting bracket.

This section describes how to use the Vivo 45 LS and a trolley with mounting bracket.

Mount and dismount the Vivo 45 LS as shown in the picture:



The bottom plate is mounted to the trolley using two screws.

Be careful when handling the trolley with the ventilator mounted, in order to avoid any risk of the trolley falling. The trolley can be tipped 10° and return to vertical position, when loaded in accordance with the weight specifications below.

#### **WARNING!**



The maximum total weight of the trolley and added accessories is 37 kg (82 lbs). (Trolley base weight = 12 kg (26 lbs), maximum externally added load = 25 kg (55 lbs).)

- The maximum load of the trolley basket is 0.9 kg (2 lbs).
- The maximum load of the IV-pole is 3 kg (6.5 lbs).
- The maximum load of the trolley rail is 9 kg (20 lbs).
- The maximum load of the E-cylinder holder is 7.9 kg (17.5 lbs).

No maintenance is required.

#### **5.7.11**

#### **Using the Click-in Humidifier**

#### **WARNING!**



Read the chapter 2.6 *Humidification*, page 19 before using the Vivo 45 LS with the humidifier.

#### **CAUTION!**



The click-in humidifier and the circuit heating operates on the AC power source only. If the AC power source fails and the internal or the external battery activates, the click-in humidifier and the circuit heating will be turned off automatically.

The click-in humidifier is intended to humidify the patient air. It is intended for non-invasive use only. The click-in water chamber is for single patient use only. Reusing a water chamber for a new patient might cause a risk of cross-contamination. The Vivo 45 LS shall not be moved with a filled water chamber installed.

About the Click-in Humidifier

The information in the table below is applicable to the recommended breathing system configuration, which is the click-in humidifier and the heated circuit.

| Property   | Value   |
|--|---|
| Humidifier classification                        | ISO 80601-2-74:2021, Class 2  |
| Rated Flow                                       | Max 50 l/min  |
| Operating Conditions                             | +5°C to +40°C. Humidity: max 90% RH, non-condensing                                     |
| Max humidification output                        | > 10 mg/l   |
| Duration of Operation between Humidifier refills | Default setting (3): 16 hours and 40 minutes<br>Max setting (5): 8 hours and 40 minutes |
| Static temperature stability *                   | ±2°C  |
| Measurement uncertainty                          | ±0.5°C  |

\* The static temperature stability have been measured at the patient port, when using the attachable humidifier. The measurement conforms to ISO 80601–2–74:2021 and discloses the value for the worst case breathing gas pathway configuration.

5.7.11.1 Adding Water to the Water Chamber

CAUTION!



Use only distilled or sterilised water or boiled, chilled tap water in the humidifier water chamber. This is to reduce mineral deposits and maximize the life of the water chamber.



Do not fill the water chamber with hot water.



Do not overfill the water chamber. Fill only the water chamber to the maximum level indicated on the water chamber.



Always ensure the lid with seal is properly mounted after filling and reassembling the water chamber. Also check that the water chamber is correctly docked in place and locked to the ventilator.



Avoid to remove the seal from the lid at normal, daily usage.

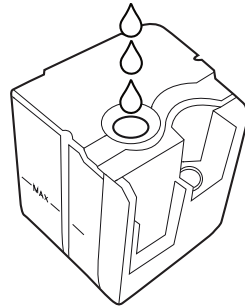


Make sure all parts are dry before the ventilator is connected to the mains and put into operation.

- 1 Detach the water chamber, see 5.7.11.4 *Detaching the Water Chamber*, page 87.
- 2 Inspect the water chamber for damages, dirt or deposits. Clean if required, see 5.7.11.6 *Cleaning the Water Chamber*, page 88. If the water chamber is damaged, replace it before use.
- 3 Fill water to the chamber, by filling through one of the airway connections.

Make sure not to fill above the Max indication. A water chamber filled to the maximum level contains approximately 350 ml

You can also remove the lid and fill water through the top of the chamber.



### 5.7.11.2 Installing the Water Chamber



#### **CAUTION!**

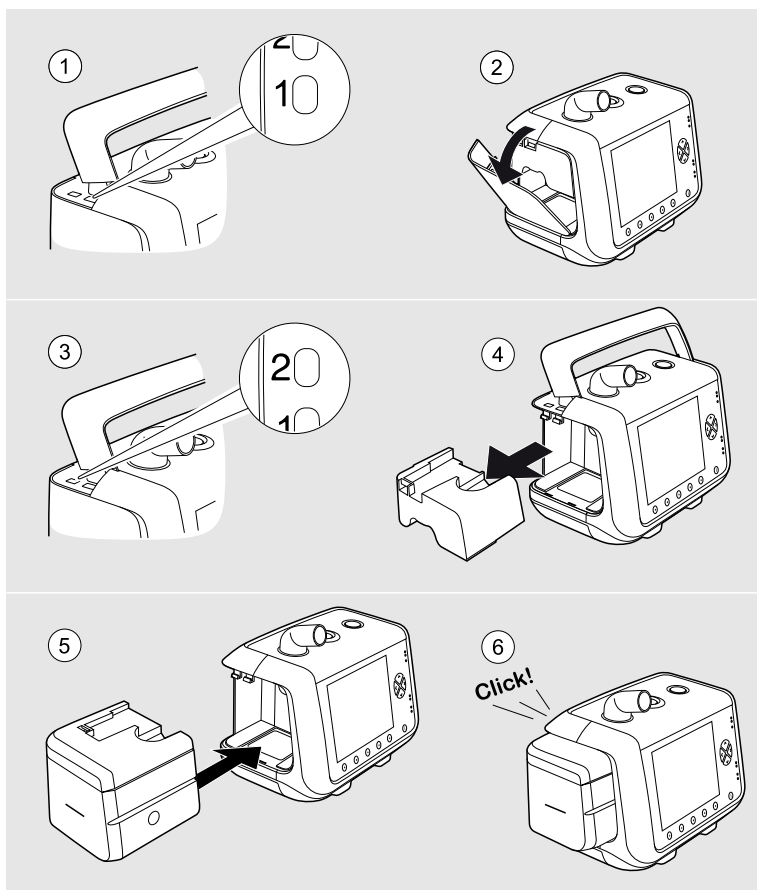
Do not switch on the humidifier without a filled water chamber in order to avoid burn or damage to the humidifier's electronics.



#### **NOTE**

If the ventilator is equipped with click-in battery, remove it before installing the water chamber.

Follow the instructions in the illustration below to install the water chamber to the ventilator



### CAUTION!

Always make sure the water chamber is in correct position before use.  
Store the airway bypass unit in a clean and dust free environment.

### 5.7.11.3 Activating the Humidification

The water chamber must be installed in order to access the humidifier setting on the ventilator menu, both in clinical and home mode. If the water chamber is disconnected and reconnected after usage, the ventilator will remember the humidity setting used.

The click-in humidifier only operates during treatment. When the ventilator is in standby mode, the humidification is paused.

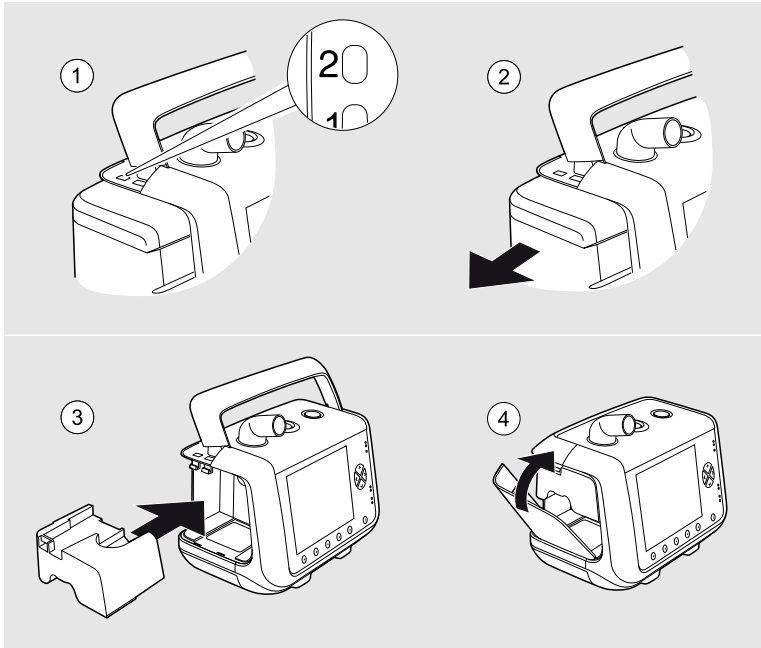
#### Prerequisites

- The water chamber shall be filled with water and attached.
- The ventilator shall be connected to the mains power supply



1. On the **Others** page, select **Heat/Humidity Settings**.
2. Select **Humidifier Setting** and set the level of humidification. 1 is the lowest level 5 is the highest level.
3. Select **Humidifier** and set it to **On**.
4. The humidifier is now activated and will start to operate when the treatment starts.

#### 5.7.11.4 Detaching the Water Chamber



#### CAUTION!

Always insert the air bypass unit after disconnecting the water chamber.



#### WARNING!

Always stop treatment before detaching or attaching the water chamber. Make sure the Vivo 45 LS with the attached water chamber is placed lower than the patient and on a flat surface. This is to prevent personal injury from accidental spillage or from excess water or condensation flowing down the patient tube and into the patient's interface.



Never add or pour out water from the water chamber when it is attached to the ventilator. If there is water outside of the water chamber after filling, dry it using a lint-free cloth before reconnecting it to the ventilator.

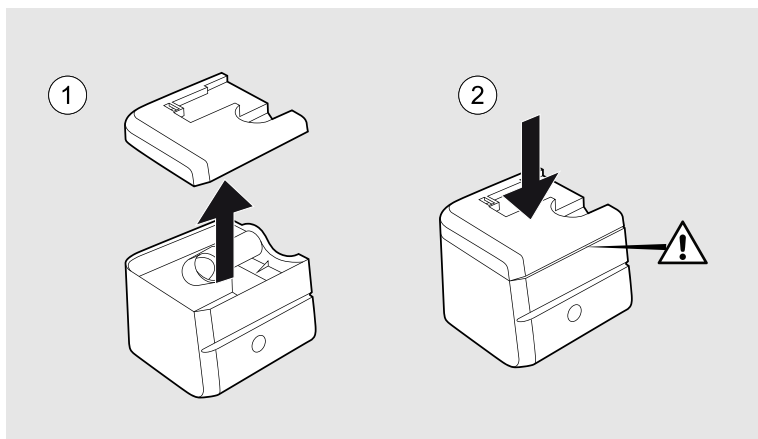


## WARNING!

To avoid burn injury, be careful not to touch the heater plate or the heated water in the water chamber when the humidifier is switched on or has not yet cooled down. Wait 10 minutes for the heater plate and water to cool.

### 5.7.11.5 Opening the Water Chamber

The water chamber lid shall be opened when manually emptying or cleaning the water chamber



## CAUTION!

Always make sure the lid of the water chamber is totally sealed.

### 5.7.11.6 Cleaning the Water Chamber

The cleaning and disinfection intervals should be established by the care provider, based on the care provider's infection control procedures.

1 Open the water chamber as described in 5.7.11.5 *Opening the Water Chamber*, page 88.

2 Clean the parts of the water chamber either by hand using a mild detergent or in a dishwasher without dishwashing detergent. Max. temperature: 60°C (140°F).

3 If there are mineral deposits inside the water chamber, dissolve them using warm water and citric acid for 30 minutes.

For disinfecting the water chamber, use any of the agents listed below. Follow the provider's instructions. The water chamber will withstand at least 20 disinfections without degradation.

| Disinfection Agent       | Duration   |
|--------------------------|------------|
| Gigasept® FF 5% solution | 15 minutes |
| Steranios 2% solution    | 10 minutes |

### 5.7.12 Using the Patient Circuit with Heated Wire

The ventilator may be used with the accessory *Patient Circuit, Heated Wire with Cable Connector*.

When the heated circuit is used, the time for the patient air temperature to reach the set temperature from a starting temperature of (23±2)°C may be up to 3 minutes.

## Prerequisites

The wire heating only operates during treatment. When the ventilator is in standby mode, the wire heating is paused.



Read the User Instruction for the Patient Circuit, Heated Wire with Cable Connector before using the patient circuit.

### 5.7.12.1 Connecting the Patient Circuit

Connect the circuit as described in 4.4 *Connecting the Patient Circuit*, page 35.

When the circuit is connected, continue with activating the circuit heating.

### 5.7.12.2 Activating the Circuit Heating

The ventilator shall be connected to the mains power supply

- 1 On the **Others** page, select **Heat/Humidity Settings**.
- 2 Select **Heated Circuit Temp** and set the temperature according to the respiratory therapist's prescription.
- 3 Select **Circuit Heating** and set it to **On**.

The circuit heating is now activated and will start to operate when the treatment starts.

## 6

## Alarms

### WARNING!



The adjustable alarm settings should be re-evaluated whenever a change in settings is made on the ventilator.

### CAUTION!



Never leave a patient unattended during an alarm condition.



Setting alarm limits to extreme values could put the patient at risk.  
Permitted distributed alarm systems are Vivo 45 LS remote alarm with cable and Vivo 45 LS nurse call cables provided by Breas Medical only.

### NOTE



The alarm settings are maintained during an extended power failure.

This chapter describes the alarm functions used for the ventilator.

## 6.1

## Alarm Function

The alarm function of the ventilator consists of the alarm LEDs on the front panel, an audible alarm, and messages on the display (see the front panel section for an overview of the position of the LEDs).

### 6.1.1

### Alarm Indication

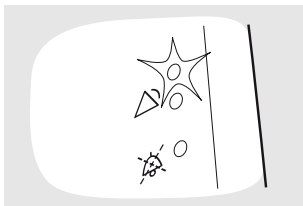
As soon as an alarm condition is detected, the ventilator main unit and the remote alarm unit (if connected) will alarm without delay.

When an alarm condition arises, the alarm is indicated in three ways:

#### Colour LED on the panel

Indicates the priority of the active alarm condition.

- High priority: red colour, flashing twice per second.
- Medium priority: yellow colour, flashing every 2 seconds.



#### Alarm text in display

Displays the name of the active alarm condition.

If several alarm conditions have been reached, the alarm descriptions are rolling in the display, displaying the alarms with highest priority first. A ">>" symbol is indicating that more than 1 alarm is set.

For detailed information about an alarm, press the info button at the lower left corner of the front.

Rebreathing

## Audible signals

- **High priority:** 3 signals followed by 2 more. The signal sequence is repeated with a 0.5 second pause and thereafter a 3 second pause.



- **Function failure:** Same signal as the high priority alarm or a constant signal, depending on the kind of function failure.

- **Medium priority:** 3 signals, with a lower frequency than the high priority alarm. The signal sequence repeats after a 6 second pause.



- **Information:** 1 signal with a low frequency. The signal is repeated after a 5 second pause and stopped after 5 sequences.



The power failure alarm sounds in the case of power failure.

If the external DC falls below the warning limit and it is the last power source, the Low External DC warning is displayed.

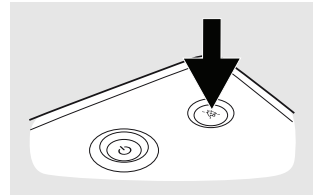
If a battery that is the last power source falls below the warning limit, the Low Last Power Source alarm is set.

### 6.1.2

#### Audible Signal Pause

The audible signal of an active alarm can be paused for 60 seconds by pressing the Audio Pause button. The audible signal can be reactivated by pressing the Audio Pause button again.

If a new alarm condition occurs during the audio pause period, the audible signal will be reactivated.



### 6.1.3

#### Audible Signal Presilence

The audible signal can be turned off for the coming 2 minutes.

##### CAUTION!

During the presilence period, any new alarms will only be indicated by the visual signals, the audible signal will not be activated.

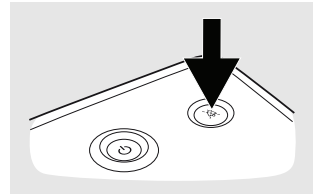
##### NOTE

Power Failure and Function Failure alarms are not affected by presilence and, if triggered, will sound during the presilence period.



1 Press and hold the Audio Pause button for about 3 seconds.

⇒ A confirmation request is displayed.



2 Press OK to confirm.

#### 6.1.4 Alarm Reset

An alarm will automatically be reset once the cause of the alarm has been corrected.

In the alarm descriptions, read the *Possible cause* information and perform corrective actions, if applicable.



#### **WARNING!**

If an alarm condition cannot be corrected, discontinue use and refer the ventilator for service.

### 6.2 Operator's Position

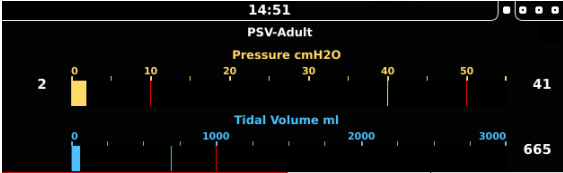
To receive the audible part of an alarm, the operator's position should be within audible range from the ventilator, depending on the set audible alarm level.

To receive the visual part of an alarm and its priority, the operator's position should be within a distance of 4 metres (13 feet) from the ventilator, and within an angle of 30° to the normal of the ventilator display.


### 6.3 Physiological Alarms

The ventilator only enables the alarms that are relevant for the used treatment. If changing modes or treatment settings, review the alarm settings.

6.3.1 High Pressure Alarm

| Property           | Description  |
|--------------------|--|
| Alarm text         | High Pressure  |
| Priority           | High   |
| Alarm condition    | A High Pressure alarm will be given when the patient pressure reaches the set High Pressure alarm limit for three consecutive breaths.<br>It will also be given if pressure exceeds 75 cmH <sub>2</sub> O. |
| Possible cause     | <ul style="list-style-type: none"><li>• Mismatch between pressure setting and alarm setting.</li><li>• Coughing during inspiration.</li><li>• Changes in airway resistance and or compliance.</li></ul>    |
| Reset criteria     | A full breath is performed with maximum pressure below the alarm limit.  |
| Ventilator action  | The Vivo 45 LS will continue treatment according to the current settings. The actual breath is however terminated if the High Pressure alarm limit is reached.   |
| Setting range      | <ul style="list-style-type: none"><li>• 5 cmH<sub>2</sub>O to 70 cmH<sub>2</sub>O</li></ul> <p>Note that the High pressure alarm cannot be set lower than the value set for the Low pressure alarm.</p>    |
| Setting resolution | Below 10 cmH <sub>2</sub> O: 0.5 cmH <sub>2</sub> O<br>Above 10 cmH <sub>2</sub> O: 1.0 cmH <sub>2</sub> O   |
| Setting Display    | The High Pressure alarm setting is displayed by a red line in the pressure bar graph. <div></div>                        |

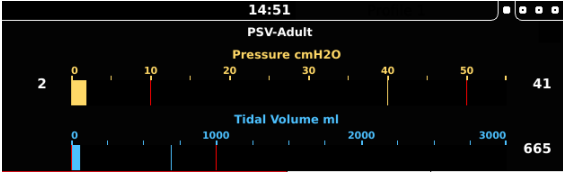6.3.2 Low Pressure Alarm

| Property           | Description  |
|--------------------|--|
| Alarm text         | Low Pressure   |
| Priority           | High   |
| Alarm condition    | A Low Pressure alarm will be given when the Vivo 45 LS pressure fails to reach the low pressure alarm limit for 15 seconds.  |
| Possible cause     | <ul style="list-style-type: none"><li>Disconnection of patient circuit.</li><li>Mismatch between pressure setting and alarm setting.</li><li>Leakage from the mask or other components of the patient circuit.</li></ul> |
| Reset criteria     | The pressure rises above the alarm limit.  |
| Ventilator action  | The Vivo 45 LS will continue treatment according to the current settings.  |
| Setting range      | <ul style="list-style-type: none"><li>1 cmH<sub>2</sub>O to 50 cmH<sub>2</sub>O</li></ul> <p>Note that the Low pressure alarm cannot be set higher than the value set for the High pressure alarm.</p>                   |
| Setting resolution | Below 10 cmH <sub>2</sub> O: 0.5 cmH <sub>2</sub> O<br>Above 10 cmH <sub>2</sub> O: 1.0 cmH <sub>2</sub> O   |
| Setting display    | The Low Pressure alarm setting is displayed by a red line in the pressure bar graph. <div></div>  |



### 6.3.3

#### High PEEP Alarm

| Property          | Description   |
|-------------------|---|
| Alarm text        | <b>High PEEP</b>  |
| Priority          | Medium  |
| Alarm condition   | A High PEEP alarm will be given when the measured PEEP is 30% above the set value for more than 15 seconds  |
| Possible cause    | <ul style="list-style-type: none"> <li>Blocked leakage port.</li> <li>Too short expiratory time.</li> <li>Changes in airway resistance and or compliance.</li> <li>Malfunction of the exhalation valve.</li> <li>Blocked exhalation valve.</li> </ul> |
| Reset criteria    | PEEP has gone below the alarm limit (lower than 30% above the set value).   |
| Ventilator action | The Vivo 45 LS will continue treatment according to the current settings.   |
| Setting range     | <ul style="list-style-type: none"> <li>On</li> <li>Off</li> </ul>   |

### 6.3.4

#### Low PEEP Alarm

| Property          | Description  |
|-------------------|--|
| Alarm text        | <b>Low PEEP</b>  |
| Priority          | Medium   |
| Alarm condition   | A Low PEEP alarm will be given when the measured PEEP is 30% below the set value for more than 60 seconds          |
| Possible cause    | <ul style="list-style-type: none"> <li>Excessive leakage.</li> <li>Malfunction of the exhalation valve.</li> </ul> |
| Reset criteria    | PEEP has gone above the alarm limit (higher than 30% below the set value).   |
| Ventilator action | The Vivo 45 LS will continue treatment according to the current settings.  |
| Setting range     | <ul style="list-style-type: none"> <li>On</li> <li>Off</li> </ul>  |

### 6.3.5 High Vt<sub>i</sub> (High Inspired Tidal Volume Alarm)

| Property           | Description  |
|--------------------|--|
| Alarm text         | <b>High Vti</b>  |
| Priority           | Medium   |
| Alarm Condition    | A High Inspired Tidal Volume alarm will be given when the monitored Inspired Tidal Volume exceeds the set limit for the High Inspired Tidal Volume alarm for 15 seconds.   |
| Possible cause     | <ul style="list-style-type: none"><li>• Mismatch between Inspired Tidal Volume and alarm setting.</li><li>• Pressure settings causing the Inspired Tidal Volume to exceed the set alarm level.</li><li>• Leakage from the mask or other components of the patient circuit.</li><li>• Mismatch between selected and used patient circuit.</li></ul> |
| Reset criteria     | When inspired tidal volume is below set alarm limit  |
| Setting range      | <ul style="list-style-type: none"><li>• Adult mode: 150 ml to 2500 ml</li><li>• Paediatric mode: 30 ml to 600 ml</li><li>• Off</li></ul>   |
| Setting resolution | 10 below 600 ml, 100 above 600 ml  |
| Ventilator action  | The ventilator will continue treatment with the same settings.   |

### 6.3.6 Low Vt<sub>i</sub> Alarm (Low Inspired Tidal Volume)

| Property           | Description   |
|--------------------|---|
| Alarm text         | <b>Low Vti</b>  |
| Priority           | High  |
| Alarm Condition    | A Low Inspired Tidal Volume alarm will be given when the monitored Inspired Tidal Volume fails to reach the set limit for the Low Inspired Tidal Volume alarm for 15 seconds. |
| Possible cause     | <ul style="list-style-type: none"><li>• Mismatch between Inspired Tidal Volume and Alarm setting.</li><li>• Changes in airway resistance and or compliance.</li></ul>         |
| Setting range      | <ul style="list-style-type: none"><li>• Adult mode: 100 ml to 2000 ml</li><li>• Paediatric mode: 20 ml to 500 ml</li><li>• Off</li></ul>                                      |
| Setting resolution | 10 below 600 ml, 100 above 600 ml   |
| Reset criteria     | A full breath above set alarm limit   |
| Ventilator action  | The ventilator will continue treatment with the same settings.  |

## 6.3.7

**High MV<sub>i</sub> Alarm (High Inspired Minute Volume Alarm)**

| Property           | Description  |
|--------------------|--|
| Alarm text         | <b>High MV<sub>i</sub></b>   |
| Priority           | Medium   |
| Alarm condition    | A High Inspired Minute Volume alarm will be given when the monitored inspired minute volume exceeds the set limit for 15 seconds.  |
| Possible cause     | <ul style="list-style-type: none"> <li>• Mismatch between Breath Rate, Inspired Tidal Volume settings and the alarm setting.</li> <li>• Increased Breath Rate.</li> <li>• Leakage around the mask or within one of the components of the circuit.</li> </ul> |
| Reset criteria     | When inspired minute volume is below set alarm limits  |
| Setting range      | <ul style="list-style-type: none"> <li>• Adult mode: 1.0 to 40 l/min</li> <li>• Paediatric mode: 1.0 to 20 l/min</li> <li>• Off</li> </ul>   |
| Setting resolution | 0.5 l/min  |
| Ventilator action  | The ventilator will continue treatment with the same settings.   |

## 6.3.8

**Low MV<sub>i</sub> (Low Inspired Minute Volume Alarm)**

| Property           | Description   |
|--------------------|---|
| Alarm text         | <b>Low MV<sub>i</sub></b>   |
| Priority           | High  |
| Alarm condition    | A Low Inspired Minute Volume alarm will be given when the monitored minute volume does not reach the alarm limit for 15 seconds.  |
| Possible cause     | <ul style="list-style-type: none"> <li>• Mismatch between Breath Rate and Inspired Tidal Volume settings and the alarm setting.</li> <li>• Changes in airway resistance and or compliance.</li> <li>• Decreased Breath Rate.</li> </ul> |
| Setting range      | <ul style="list-style-type: none"> <li>• 1.0 l/min to 30 l/min (Adult mode)</li> <li>• 0.1 l/min to 10 l/min (Paediatric mode)</li> <li>• Off</li> </ul>  |
| Setting resolution | 0.1 l up to 1.0 l, 0.5 l above 1.0 l.   |
| Ventilator action  | The ventilator will continue treatment with the same settings.  |

### 6.3.9 High Vt<sub>e</sub> Alarm (High Expired Tidal Volume)

| Property           | Description  |
|--------------------|--|
| Alarm text         | <b>High Vte</b>  |
| Priority           | Medium   |
| Alarm condition    | A High Expired Tidal Volume alarm will be given when the monitored Expired Tidal Volume exceeds the alarm limit for 15 seconds.  |
| Possible cause     | <ul style="list-style-type: none"><li>• Mismatch between Expired Tidal Volume and alarm setting.</li><li>• Mismatch between selected and used patient circuit.</li><li>• Pressure settings causing the Expired Tidal Volume to exceed the set alarm level.</li></ul> |
| Setting range      | <ul style="list-style-type: none"><li>• 150 ml to 2500 ml (Adult mode)</li><li>• 30 ml to 600 ml (Paediatric mode)</li><li>• Off</li></ul>   |
| Setting resolution | 10 below 600 ml, 100 above 600 ml  |
| Ventilator action  | The ventilator will continue treatment with the same settings.   |

### 6.3.10 Low Vt<sub>e</sub> Alarm (Low Expired Tidal Volume)

| Property           | Description  |
|--------------------|--|
| Alarm text         | <b>Low Vte</b>   |
| Priority           | High   |
| Alarm Condition    | A Low Expired Tidal Volume alarm will be given when the monitored Expired Tidal Volume fails to reach the set limit for the Low Expired Tidal Volume alarm for 15 seconds.   |
| Possible cause     | <ul style="list-style-type: none"><li>• Mismatch between Expired Tidal Volume and Alarm setting.</li><li>• Changes in airway resistance and or compliance.</li><li>• Leakage around the mask or within one of the components of the circuit.</li></ul> |
| Reset criteria     | Full breath above set alarm limit  |
| Setting range      | <ul style="list-style-type: none"><li>• 100 ml to 2000 ml (Adult mode)</li><li>• 20 ml to 500 ml (Paediatric mode)</li><li>• Off</li></ul>   |
| Setting resolution | 10 below 600 ml, 100 above 600 ml  |
| Ventilator action  | The ventilator will continue treatment with the same settings.   |

### 6.3.11 High MV<sub>e</sub> (High Expired Minute Volume Alarm)

| Item               | Description   |
|--------------------|---|
| Alarm text         | <b>High MVe</b>   |
| Priority           | Medium  |
| Alarm condition    | A High Expired Minute Volume alarm will be given when the monitored expired minute volume exceeds the alarm limit for 15 seconds.                             |
| Possible cause     | <ul style="list-style-type: none"><li>• Mismatch between Breath Rate, Tidal Volume settings and the alarm setting.</li><li>• Increased Breath Rate.</li></ul> |
| Setting range      | <ul style="list-style-type: none"><li>• Adult mode: 1.0 to 40 l/min</li><li>• Paediatric mode: 1.0 to 20 l/min</li><li>• Off</li></ul>                        |
| Setting resolution | 0.5 l/min   |
| Ventilator action  | The ventilator will continue treatment with the same settings.  |

### 6.3.12 Low MV<sub>e</sub> Alarm (Low Expired Minute Volume)

| Property           | Description  |
|--------------------|--|
| Alarm text         | <b>Low MVe</b>   |
| Priority           | High   |
| Alarm condition    | A Low Expired Minute Volume alarm will be given when the monitored minute volume is below the alarm limit for more than 15 seconds.  |
| Possible cause     | <ul style="list-style-type: none"><li>• Mismatch between Breath Rate and Tidal Volume settings and the alarm setting.</li><li>• Changes in airway resistance and or compliance.</li><li>• Decreased Breath Rate.</li><li>• Leakage around the mask or within one of the components of the circuit.</li></ul> |
| Setting range      | <ul style="list-style-type: none"><li>• 1.0 l/min to 30 l/min (Adult mode)</li><li>• 0.1 l/min to 10 l/min (Paediatric mode)</li><li>• Off</li></ul>   |
| Setting resolution | 0.1 l up to 1.0 l, 0.5 l above 1.0 l.  |
| Ventilator action  | The ventilator will continue treatment with the same settings.   |

### 6.3.13 High Breath Rate Alarm

| Property           | Description  |
|--------------------|--|
| Alarm text         | <b>High Breath Rate</b>  |
| Priority           | Medium   |
| Alarm condition    | A High Breath Rate alarm will be given when the alarm limit has been exceeded for 15 seconds.  |
| Possible cause     | <ul style="list-style-type: none"><li>• Mismatch between the Breath Rate setting and the alarm setting.</li><li>• Increased Breath Rate.</li><li>• Too sensitive setting of the inspiratory trigger setting.</li></ul> |
| Reset criteria     | The breath rate goes below the alarm limit.  |
| Ventilator action  | The Vivo 45 LS will continue treatment according to the current settings.  |
| Setting range      | <ul style="list-style-type: none"><li>• 10 bpm to 70 bpm (Adult mode)</li><li>• 10 bpm to 99 bpm (Paediatric mode)</li><li>• Off</li></ul>   |
| Setting resolution | 1 bpm.   |

### 6.3.14 Low Breath Rate Alarm

| Property           | Description   |
|--------------------|---|
| Alarm text         | <b>Low Breath Rate</b>  |
| Priority           | High  |
| Alarm condition    | A Low Breath Rate alarm will be given when the delivered total breath rate is below the alarm limit for 15 seconds.   |
| Possible cause     | <ul style="list-style-type: none"><li>• Mismatch between the Breath Rate setting and the alarm setting.</li><li>• The patient cannot trigger breaths because the inspiratory trigger setting is too high.</li><li>• Decrease in the patient's spontaneous breathing.</li><li>• Circuit disconnection.</li></ul> |
| Reset criteria     | The breath rate goes above the alarm limit.   |
| Ventilator action  | The Vivo 45 LS will continue treatment according to the current settings.   |
| Setting range      | <ul style="list-style-type: none"><li>• 4 bpm to 30 bpm (Adult mode, non-MPV modes)</li><li>• 1 bpm to 30 bpm (Adult mode MPV modes)</li><li>• 6 bpm to 50 bpm (Paediatric mode)</li><li>• Off</li></ul>  |
| Setting resolution | 1 bpm.  |

6.3.15 Apnea Alarm

| Property           | Description   |
|--------------------|---|
| Alarm text         | <b>Apnea</b>  |
| Priority           | High  |
| Alarm condition    | An Apnea alarm will be given when no patient-triggered breath is detected for the set period of time.<br>The Apnea alarm is only available if the Inspiratory trigger is activated.                             |
| Possible cause     | <ul style="list-style-type: none"><li>• Patient stopped breathing.</li><li>• Patient decreases spontaneous breathing.</li><li>• Circuit disconnection.</li><li>• Inspiratory Trigger is set too high.</li></ul> |
| Reset criteria     | Inspiratory effort detected by the Vivo 45 LS.  |
| Ventilator action  | The Vivo 45 LS will continue treatment according to the current settings.   |
| Setting range      | <ul style="list-style-type: none"><li>• 5 to 60 s. (Non MPV mode)</li><li>• 15 to 900 s. (MPV mode)</li><li>• Off</li></ul>   |
| Setting resolution | 5 s below 15 s.<br>15 s above 15 s.<br>MPV mode:<br>15 s below 60 s.<br>60 s above 60 s.  |

6.3.16 Disconnection Alarm  
**CAUTION!**



No single alarm can reliably detect all disconnections due to the number of possible combinations of therapy settings, circuit configurations and patient interfaces. To verify that patient disconnection can be detected, including if the patient interface becomes accidentally detached from the patient, it is advised to test the functionality of the Disconnection Alarm upfront with the complete set-up as used during treatment, including items such as filters, circuit, connectors, interface (mask, cannula etc.)

| Property          | Description  |
|-------------------|--|
| Alarm text        | <b>Disconnection</b>   |
| Priority          | High   |
| Alarm condition   | A Disconnection alarm will be given when the measured flow exceeds the expected leakage flow at the set Pressure for 15 seconds.   |
| Possible cause    | <ul style="list-style-type: none"> <li>• Too high leakage in the patient circuit.</li> <li>• The patient has removed the mask.</li> <li>• Circuit disconnection.</li> <li>• Pilot pressure tube disconnection</li> </ul> |
| Reset criteria    | The leakage is back within limits.   |
| Ventilator action | The Vivo 45 LS will continue treatment according to the current settings   |
| Setting range     | <ul style="list-style-type: none"> <li>• On</li> <li>• Off</li> </ul>  |



### 6.3.17 Rebreathing Alarm

| Property          | Description   |
|-------------------|---|
| Alarm text        | <b>Rebreathing</b> (with leakage circuit)<br><b>Rebreathing</b> (with active exhalation valve circuit)  |
| Priority          | High (with leakage circuit)<br>Medium (with active exhalation valve circuit)  |
| Alarm condition   | <b>Leakage Circuit</b><br>A Rebreathing alarm will be given if the intentional leakage is too low for more than 15 seconds.<br><b>Exhalation valve circuit</b><br>A Rebreathing alarm will be given if the exhalation valve is obstructed for more than 10 consecutive breaths.<br><b>MPV circuit</b><br>A Rebreathing alarm will be given if air returns into the ventilator for more than 10 consecutive breaths. |
| Possible cause    | <ul style="list-style-type: none"> <li>Obstructed or occluded patient circuit.</li> <li>Incorrect patient circuit.</li> <li>Patient exhales through mouthpiece.</li> <li>Obstructed or removed CO<sub>2</sub> port from leakage circuit.</li> <li>For active exhalation valve circuits: Disconnected pilot pressure line.</li> </ul>  |
| Reset criteria    | The leakage is back within limits.<br>The bias flow is restored   |
| Ventilator action | The Vivo 45 LS will continue treatment according to the current settings.   |
| Setting range     | <ul style="list-style-type: none"> <li>On</li> <li>Off</li> </ul>   |

### 6.3.18 Obstruction Alarm

| Property          | Description   |
|-------------------|---|
| Alarm text        | <b>Obstruction</b>  |
| Priority          | High  |
| Alarm condition   | An Obstruction alarm will be given if the inspiratory breathing tube becomes blocked and remains blocked for 2 consecutive breaths.   |
| Ventilator action | With each breath cycle, upon detection of an obstruction the ventilator will reduce the airway pressure to the set PEEP. Treatment will resume with the start of the next breath cycle. |
| Reset Criteria    | When the monitored compliance and resistance become normal after a breath.  |
| Setting Range     | <ul style="list-style-type: none"> <li>High</li> <li>Low</li> <li>Off</li> </ul>  |

6.3.19 High FiO<sub>2</sub> Alarm

| Property           | Description  |
|--------------------|--|
| Alarm text         | <b>High FiO2</b>   |
| Priority           | Medium   |
| Alarm condition    | A High FiO <sub>2</sub> alarm will be given when the measured FiO <sub>2</sub> exceeds the alarm limit for 30 seconds. |
| Possible cause     | <ul style="list-style-type: none"><li>Increased oxygen inflow.</li><li>Decreased minute ventilation.</li></ul>         |
| Reset criteria     | FiO <sub>2</sub> goes below the alarm limit  |
| Setting range      | <ul style="list-style-type: none"><li>21% to 100%</li><li>Off</li></ul>  |
| Setting resolution | 1%   |
| Ventilator action  | The ventilator will continue treatment with the same settings.   |

6.3.20 Low FiO<sub>2</sub> Alarm

| Property           | Description  |
|--------------------|--|
| Alarm text         | <b>Low FiO2</b>  |
| Priority           | High   |
| Alarm condition    | A Low FiO <sub>2</sub> alarm will be given when the measured FiO <sub>2</sub> is below the alarm limit for 30 seconds.   |
| Possible cause     | <ul style="list-style-type: none"><li>Decreased oxygen inlet.</li><li>Disconnection of oxygen inlet.</li><li>Increased minute ventilation.</li><li>High leakage.</li></ul> |
| Setting range      | <ul style="list-style-type: none"><li>21% to 100%</li><li>Off</li></ul>  |
| Setting resolution | 1%   |
| Ventilator action  | The ventilator will continue treatment with the same settings.   |

### 6.3.21 High SpO<sub>2</sub> Alarm

| Property           | Description  |
|--------------------|--|
| Alarm text         | <b>High SpO2</b>   |
| Priority           | Medium   |
| Alarm condition    | A High SpO <sub>2</sub> alarm will be given when the measured SpO <sub>2</sub> exceeds the alarm limit for 30 seconds. |
| Possible cause     | Too high flow of bleed-in oxygen.  |
| Reset criteria     | The SpO <sub>2</sub> value goes back below the alarm limit.  |
| Ventilator action  | The Vivo 45 LS will continue treatment according to the current settings.  |
| Setting range      | <ul style="list-style-type: none"><li>90 % to 100 %</li><li>Off</li></ul>  |
| Setting resolution | 1 %  |

This alarm requires a connected SpO<sub>2</sub> sensor.

### 6.3.22 Low SpO<sub>2</sub> Alarm

| Property           | Description   |
|--------------------|---|
| Alarm text         | <b>Low SpO2</b>   |
| Priority           | High  |
| Definition         | A Low SpO <sub>2</sub> alarm will be given when the measured SpO <sub>2</sub> is below the alarm limit for 30 seconds.  |
| Possible cause     | <ul style="list-style-type: none"><li>Too low flow of bleed-in oxygen.</li><li>Oxygen inlet is disconnected.</li><li>Delivered tidal volumes are too small.</li></ul> |
| Setting range      | 85% to 100%   |
| Setting resolution | 1%  |
| Ventilator action  | The ventilator will continue treatment with the same settings.  |

This alarm requires a connected SpO<sub>2</sub> sensor.

### 6.3.23 High EtCO<sub>2</sub> Alarm

| Property           | Description  |
|--------------------|--|
| Alarm text         | <b>High EtCO<sub>2</sub></b>   |
| Priority           | High   |
| Alarm condition    | A High EtCO <sub>2</sub> alarm will be given when the measured EtCO <sub>2</sub> exceeds the alarm limit for 30 seconds.   |
| Possible cause     | <ul style="list-style-type: none"><li>• Alarm limit is set too low.</li><li>• Breath Rate too low.</li><li>• Delivered Tidal Volume too low.</li><li>• Excessive dead space between patient and exhalation valve/ leakage port.</li><li>• Exhalation valve / leakage port is occluded.</li></ul> |
| Setting range      | 1 to 99 mmHg<br>Off  |
| Setting resolution | 1 mmHg   |
| Ventilator action  | The ventilator will continue treatment with the same settings.   |

This alarm requires a connected EtCO<sub>2</sub> sensor

### 6.3.24 Low EtCO<sub>2</sub> Alarm

| Property           | Description  |
|--------------------|--|
| Alarm text         | <b>Low EtCO<sub>2</sub></b>  |
| Priority           | Medium   |
| Alarm condition    | A Low EtCO <sub>2</sub> alarm will be given when the measured EtCO <sub>2</sub> is below the alarm limit for 30 seconds.   |
| Possible cause     | <ul style="list-style-type: none"><li>• Alarm limit is set too high.</li><li>• Ventilator disconnection.</li><li>• Excessive leakage in the Patient circuit/Interface.</li><li>• Partial obstruction of the airways.</li><li>• Breath Rate too high.</li><li>• Delivered Tidal Volume too high.</li><li>• Self triggering of the ventilator.</li></ul> |
| Setting range      | 1 to 99 mmHg<br>Off  |
| Setting resolution | 1 mmHg   |
| Ventilator action  | The ventilator will continue treatment with the same settings.   |

This alarm requires a connected EtCO<sub>2</sub> sensor

### 6.3.25 High InspCO<sub>2</sub> Alarm (High Inspired CO<sub>2</sub>)

| Property           | Description  |
|--------------------|--|
| Alarm text         | <b>High InspCO<sub>2</sub></b>   |
| Priority           | High   |
| Alarm condition    | A High Inspired CO <sub>2</sub> alarm will be given when the measured inspired CO <sub>2</sub> exceeds the alarm limit for 30 seconds.   |
| Possible cause     | <ul style="list-style-type: none"><li>• Alarm limit is set too low.</li><li>• Excessive dead space between patient and exhalation valve/leakage port.</li><li>• Leakage port/valve occluded.</li></ul> |
| Setting range      | 1 to 99 mmHg<br>Off  |
| Setting resolution | 1 mmHg   |
| Ventilator action  | The ventilator will continue treatment with the same settings.   |

This alarm requires a connected EtCO<sub>2</sub> sensor

### 6.3.26 High Pulse Rate Alarm

| Property           | Description   |
|--------------------|---|
| Alarm text         | <b>High Pulse Rate</b>  |
| Priority           | Medium  |
| Alarm condition    | A High Pulse Rate alarm will be given when the measured pulse rate exceeds the alarm limit for 15 seconds.  |
| Possible cause     | <ul style="list-style-type: none"><li>• Insufficient ventilatory support.</li><li>• Too low flow of bleed-in oxygen.</li><li>• The PEEP value is set too high.</li><li>• Bad positioning of the finger probe.</li></ul> |
| Reset criteria     | The pulse rate goes back below the alarm limit.   |
| Ventilator action  | The Vivo 45 LS will continue treatment according to the current settings.   |
| Setting range      | 30 to 230 bpm (beats per minute)<br>Off   |
| Setting resolution | 5 bpm (beats per minute)  |

This alarm requires a connected SpO<sub>2</sub> sensor.

### 6.3.27 Low Pulse Rate Alarm

| Property           | Description   |
|--------------------|---|
| Alarm text         | <b>Low Pulse Rate</b>   |
| Priority           | High  |
| Alarm condition    | A low pulse rate alarm will be given when the measured pulse rate goes below the alarm limit for 15 seconds.  |
| Possible cause     | <ul style="list-style-type: none"><li>• Bad positioning of the finger probe.</li><li>• Too low flow of bleed-in oxygen.</li><li>• Insufficient ventilatory support.</li></ul> |
| Reset criteria     | The pulse rate goes back above the alarm limit.   |
| Ventilator action  | The Vivo 45 LS will continue treatment according to the current settings.   |
| Setting range      | 30 to 230 bpm (beats per minute)<br>Off   |
| Setting resolution | 5 bpm (beats per minute)  |

This alarm requires a connected SpO<sub>2</sub> sensor.

### 6.3.28 PtcCO<sub>2</sub> Alarm

| Property          | Description  |
|-------------------|--|
| Alarm text        | <b>PtcCO<sub>2</sub> Outside Limits</b>  |
| Priority          | High   |
| Alarm Condition   | A PtcCO <sub>2</sub> alarm will be given when PtcCO <sub>2</sub> is outside alarm limits.<br>Check the PtcCO <sub>2</sub> monitor.   |
| Possible cause    | <ul style="list-style-type: none"><li>• External PtcCO<sub>2</sub> monitor is outside its alarm limits.</li><li>• Breath Rate needs adjustment</li><li>• Delivered Tidal Volume needs adjustment.</li><li>• Excessive dead space between patient and exhalation valve/leakage port.</li><li>• Leak port/valve occluded.</li><li>• Ventilator disconnection.</li><li>• Excessive leakage in the Patient circuit/Interface.</li><li>• Partial obstruction of the airways.</li><li>• Self triggering of the ventilator.</li></ul> |
| Ventilator action | The ventilator will continue treatment with the same settings.   |

This alarm requires a connected PtcCO<sub>2</sub> sensor.

## 6.4 Technical Alarms

### 6.4.1 Power Fail Alarm

| Property          | Description   |
|-------------------|---|
| Alarm text        | The alarm is given audibly with a tone and the display is blinking with the alarm message<br><b>Power Fail</b>  |
| Priority          | High  |
| Alarm condition   | The Power Fail alarm is given if the last power source fails to provide enough power for running the ventilator.  |
| Possible cause    | The last available power source cannot deliver power to the ventilator.<br>Battery discharged or battery failure.   |
| Reset criteria    | External power supply connected to ventilator.  |
| Ventilator action | The Vivo 45 LS stops the treatment, and gives the Power Fail alarm for at least 2 minutes.<br>If power is restored within the alarm time, the ventilator will automatically resume treatment with current settings.<br>When powered up again, the power failure will be logged. |

### 6.4.2 High Patient Air Temp. (High Patient Air Temperature)

| Property          | Description  |
|-------------------|--|
| Alarm text        | <b>High Patient Air Temp</b>   |
| Priority          | High   |
| Alarm condition   | A High Patient Air Temperature alarm will be given when the patient air temperature exceeds 43°C (109.4°F).                                    |
| Possible cause    | <ul style="list-style-type: none"><li>Blocked air inlets.</li><li>Blocked cooling air outlets.</li><li>Too high ambient temperature.</li></ul> |
| Ventilator action | The Vivo 45 LS will continue treatment.<br>If a heated circuit or the click-in humidifier is used, these will be turned off.                   |
| Reset criteria    | The temperature goes below the limit again.  |

### 6.4.3 Low Patient Air Temp. (Low Patient Air Temperature Alarm)

| Property          | Description   |
|-------------------|---|
| Alarm text        | <b>Low Patient Air Temp</b>   |
| Alarm condition   | A Low Patient Air Temperature alarm will be given when the patient air temperature is below the preset limit -30°C (-22°F). |
| Priority          | High  |
| Possible cause    | Too low ambient temperature   |
| Ventilator action | The ventilator will continue treatment with the same settings.  |

#### 6.4.4 Low Last Power Source Alarm

| Property          | Description   |
|-------------------|---|
| Alarm text        | <b>Low Last Power Source</b>  |
| Priority          | Medium  |
| Alarm condition   | This alarm will be given when the last battery source (internal battery) has 15 minutes of operating time left with current settings. |
| Ventilator action | The Vivo 45 LS will continue treatment according to the current settings.   |

#### 6.4.5 Crit. Low Last Power Source Alarm

| Property          | Description   |
|-------------------|---|
| Alarm text        | <b>Crit. Low Last Power Source</b>  |
| Alarm condition   | A Crit. Low Last Power Source alarm will be given when the last battery source (internal battery or click-in battery) has 5 minutes of operating time left with current settings. |
| Priority          | High  |
| Ventilator action | The ventilator will continue treatment with the same settings.  |
| Reset             | Connection of “higher” power source.  |

#### 6.4.6 Lost Mains Alarm

| Property          | Description  |
|-------------------|--|
| Alarm text        | <b>Lost Mains Power</b>  |
| Alarm condition   | A Mains Power Lost alarm will be given when the ventilator switched from AC power (Mains) to another power source due to AC Power (Mains) is lost. |
| Priority          | Medium   |
| Ventilator action | The ventilator will continue treatment with the same settings. An information message will be shown on the screen.                                 |
| Reset             | Confirmation by user or AC power (Mains) reconnected.  |



#### 6.4.7 Exhalation Valve Control Error Alarm

| Property        | Description   |
|-----------------|---|
| Alarm text      | <b>Exhalation Valve Control Error</b>   |
| Alarm condition | An Exhalation Valve Control Error alarm will be given when the ventilator fails to control the internal /external exhalation valve.   |
| Priority        | High  |
| Possible cause  | <ul style="list-style-type: none"><li>Exhalation valve occluded</li><li>Exhalation valve control tube disconnected</li><li>Internal function failure of the exhalation valve controls</li></ul> |
| Reset           | The pilot pressure gets a normal value.   |

#### 6.4.8 SpO<sub>2</sub> Disconnected (SpO<sub>2</sub> Sensor Failure/Disconnection Alarm)

| Property          | Description   |
|-------------------|---|
| Alarm text        | <b>SPO2 Disconnected</b>  |
| Alarm condition   | An SpO <sub>2</sub> Sensor Failure/Disconnection alarm will be given when an error signal or no signal from the SpO <sub>2</sub> sensor has been detected for 2 seconds. Check the SpO <sub>2</sub> sensor. |
| Priority          | High  |
| Possible cause    | The SpO <sub>2</sub> electronics cable has been disconnected and subsequently no communication (possibly due to disconnection) for 2 seconds.<br>Failure in the SpO <sub>2</sub> sensor.                    |
| Ventilator action | The ventilator will continue treatment with the same settings.  |
| Reset             | Confirmation by user or reconnected/changed.  |

#### 6.4.9 SpO<sub>2</sub> Signal Lost Alarm

| Property          | Description  |
|-------------------|--|
| Alarm text        | <b>SPO2 Signal Lost</b>  |
| Alarm condition   | SpO <sub>2</sub> signal lost.  |
| Priority          | High   |
| Possible cause    | Signal lost reported by SpO <sub>2</sub> electronics (due to patient removing the probe from finger, or sensor detached from SpO <sub>2</sub> electronics. |
| Ventilator action | The ventilator will continue treatment with the same settings.   |
| Reset             | User presses OK or electronics cable is disconnected by the user, or the sensor is reconnected to the finger.  |

#### 6.4.10 Poor SpO<sub>2</sub> Signal

| Property          | Description   |
|-------------------|---|
| Alarm text        | <b>Poor SPO2 Signal</b>   |
| Alarm condition   | A Poor SpO <sub>2</sub> signal alarm will be given when the SpO <sub>2</sub> signal is not correct. Check the SpO <sub>2</sub> sensor.                |
| Priority          | High  |
| Possible cause    | Artifact or low perfusion reported by SpO <sub>2</sub> electronics  |
| Ventilator action | The ventilator will continue treatment with the same settings.  |
| Reset             | OK message from SpO <sub>2</sub> electronics or SpO <sub>2</sub> electronics disconnected by user or SpO <sub>2</sub> Signal Lost alarm is triggered. |

#### 6.4.11 CO<sub>2</sub> Disconnected (CO<sub>2</sub> Sensor Failure/Disconnection Alarm)

| Property          | Description  |
|-------------------|--|
| Alarm text        | <b>CO2 Sensor Disconnected</b>   |
| Alarm condition   | A CO <sub>2</sub> Sensor Failure/Disconnection alarm will be given when communication between the ventilator and the CO <sub>2</sub> sensor has been lost for 2 seconds. Check the CO <sub>2</sub> sensor. |
| Priority          | High   |
| Possible cause    | <ul style="list-style-type: none"><li>• CO<sub>2</sub> Sensor disconnected.</li><li>• Failure in the CO<sub>2</sub> sensor.</li></ul>  |
| Ventilator action | The ventilator will continue treatment with the same settings.   |
| Reset             | Confirmation by user or reconnected/changed.   |

#### 6.4.12 CO<sub>2</sub> Accuracy Error Alarm

| Property          | Description  |
|-------------------|--|
| Alarm text        | <b>CO2 Accuracy Error</b>  |
| Alarm condition   | A CO <sub>2</sub> Accuracy Error alarm will be given when an accuracy error in the CO <sub>2</sub> measurement has occurred. |
| Priority          | High   |
| Ventilator action | The ventilator will continue treatment with the same settings.   |
| Reset             | OK message from sensor or sensor disconnected by user.   |

#### 6.4.13 Check CO<sub>2</sub> Adapter Alarm

| Property          | Description  |
|-------------------|--|
| Alarm text        | <b>Check CO2 Adapter</b>   |
| Alarm condition   | A Check CO <sub>2</sub> Adapter alarm will be given when the airway adapter is not attached correctly to the CO <sub>2</sub> sensor. Check/replace the airway adapter. |
| Priority          | High   |
| Ventilator action | The ventilator will continue treatment with the same settings.   |
| Reset             | OK message from sensor or sensor disconnected by user.   |

#### 6.4.14 CO<sub>2</sub> Sensor Error Alarm

| Property          | Description   |
|-------------------|---|
| Alarm text        | <b>CO2 Sensor Error</b>   |
| Alarm condition   | A CO <sub>2</sub> Sensor Error alarm will be given when an error in the CO <sub>2</sub> sensor has occurred.<br>Replace the CO <sub>2</sub> sensor. CO <sub>2</sub> monitoring cannot be performed in this condition. |
| Priority          | High  |
| Ventilator action | The ventilator will continue treatment with the same settings.  |
| Reset             | OK message from sensor or sensor disconnected by user.  |

#### 6.4.15 FiO<sub>2</sub> Disconnected (FiO<sub>2</sub> Sensor Failure/Disconnection Alarm)

| Property          | Description   |
|-------------------|---|
| Alarm text        | <b>FiO2 Disconnected</b>  |
| Alarm condition   | An FiO <sub>2</sub> Sensor Failure/Disconnection alarm will be given when no signal from the FiO <sub>2</sub> sensor has been detected for 2 seconds.<br>Check the FiO <sub>2</sub> sensor. |
| Priority          | High  |
| Possible cause    | <ul style="list-style-type: none"><li>• FiO<sub>2</sub> Sensor disconnected.</li><li>• Communication with the FiO<sub>2</sub> sensor failed.</li></ul>                                      |
| Ventilator action | The ventilator will continue treatment with the same settings.  |
| Reset             | Confirmation by user or reconnected/changed.  |

#### 6.4.16 Ambient Pressure Compensation Lost Alarm

| Property          | Description   |
|-------------------|---|
| Alarm text        | <b>Pressure Comp Lost</b>   |
| Priority          | Medium  |
| Alarm condition   | An Ambient Pressure Compensation Lost alarm will be given when the automatic ambient pressure compensation functionality is out of order.   |
| Ventilator action | The Vivo 45 LS will continue treatment according to the current settings.<br>Normal atmospheric pressure at sea level will be used as approximation for the temporary ambient pressure compensation. If used at other altitude, delivered and measured pressures may deviate. |
| Reset             | Reset of ventilator.  |

#### 6.4.17 Temperature Comp. Lost (Ambient Temperature Compensation Lost Alarm)

| Property          | Description   |
|-------------------|---|
| Alarm text        | <b>Temperature Comp. Lost</b>   |
| Alarm condition   | An Ambient Temperature Compensation Lost alarm will be given when the automatic ambient temperature compensation is out of order.<br>There is no communication with the air temperature sensor or the value is out of range (less than -30°C (-22°F) or more than 70°C (158°F). |
| Priority          | Medium  |
| Ventilator action | The ventilator will continue treatment with the same settings.<br>The accuracy of the volume measurement may be impaired.   |
| Reset             | Ambient temperature inside valid range.   |

#### 6.4.18 Humidity Comp. Lost (Humidity Compensation Lost Alarm)

| Property          | Description   |
|-------------------|---|
| Alarm text        | <b>Humidity Comp. Lost</b>  |
| Alarm condition   | An Humidity Compensation Lost alarm will be given when the automatic humidity compensation is out of order.<br>50% relative humidity is used for temporary compensation. If the ventilator is used at other humidities, delivered and measured pressure and flow may deviate. |
| Priority          | Medium  |
| Ventilator action | The ventilator will continue treatment with the same settings.<br>The accuracy of the volume measurement may be impaired.   |
| Reset             | Air humidity sensor values (RH and temperature) inside valid range.   |

#### 6.4.19 LED Failure Alarm

| Property          | Description  |
|-------------------|--|
| Alarm text        | <b>LED Failure</b>   |
| Alarm condition   | A LED Failure alarm will be given when one or more LED indicators on the front panel are broken. |
| Priority          | Medium   |
| Ventilator action | The ventilator will continue treatment with the same settings.                                   |
| Reset             | Power-on reset of ventilator (or repair).  |

#### 6.4.20 Low Alarm Battery Alarm

| Property          | Description  |
|-------------------|--|
| Alarm text        | <b>Low Alarm Battery</b>   |
| Alarm condition   | An alarm for <i>Low Alarm Battery</i> will be given if the alarm battery is not charged enough to have power for a <i>Power Fail</i> alarm for at least 2 minutes. |
| Priority          | Medium   |
| Ventilator action | The ventilator will continue treatment with the same settings and start charging the alarm batteries.  |
| Reset             | When alarm energy storage level is sufficient to give an alarm for at least 2 minutes.   |

#### 6.4.21 Alarm Battery Error Alarm

| Property          | Description   |
|-------------------|---|
| Alarm text        | <b>Alarm Battery Error</b>  |
| Alarm condition   | Unable to communicate with super capacitor and read super capacitor status. |
| Priority          | Medium  |
| Ventilator action | The ventilator will continue treatment with the same settings.              |
| Reset             | When triggering condition is removed.                                       |

## 6.4.22 Internal/Click-In Battery Hot Alarm

| Property          | Description  |
|-------------------|--|
| Alarm text        | <b>Internal Battery — Internal Battery Hot<br/>Click-In Battery — Click-In Battery Hot</b>   |
| Alarm condition   | An alarm for Internal/Click-In Battery Overheat in Discharge will be given when the internal or click-in battery reaches 55°C (131°F).   |
| Priority          | High   |
| Ventilator action | The ventilator will continue treatment with the same settings. Battery discharging will be disabled (by the battery electronics) once the temperature gets to 60°C (140°F). (If the battery is last power source, the ventilator will stop running). |



### NOTE

The battery electronics by manufacture stops discharge at 60°C (140°F).

## 6.4.23 Heated Circuit Temp. Alarm

| Property          | Description  |
|-------------------|--|
| Alarm text        | <b>Heated Circuit Temp.</b>  |
| Alarm condition   | A Heated Circuit temp alarm will be given when the measured temperature of the heated wire is outside the tolerance. |
| Priority          | Medium   |
| Ventilator action | The ventilator will continue treatment with the same settings.   |
| Priority          | Medium   |
| Reset             | Heated wire measured temp tolerance is inside limits.  |

## 6.4.24 High Humidifier Temp. Alarm

| Property          | Description   |
|-------------------|---|
| Alarm text        | <b>High Humidifier Temp.</b>  |
| Alarm condition   | A Humidifier High Temperature alarm will be given if the humidifier heater plate temperature exceeds 76°C (169°F) for more than 2 seconds.  |
| Priority          | Medium  |
| Ventilator action | The ventilator will turn off the click-in humidifier and then continue treatment with the same settings. A message with option to turn on the humidifier again will be displayed. |
| Reset             | The alarm is dismissed when the humidifier temperature drops below 76°C (169°F), set humidifier temperature).   |

6.4.25 Humidifier Fault Alarm

| Property          | Description   |
|-------------------|---|
| Alarm text        | <b>Humidifier Fault</b>   |
| Alarm condition   | <ul style="list-style-type: none"><li>• All humidifier enabling conditions have been satisfied for 10 minutes, and</li><li>• No humidifier setting changes have been made for 10 minutes, and</li><li>• Heater plate temperature &lt; 50°C (122 °F)</li><li>• Humidifier set temperature &gt; Ambient temperature ,and</li><li>• The heater plate temperature is more than 5°C (41 °F) below the set temperature, or the heater plate temperature &lt; -20°C (68 °F) or greater than 400°C (752 °F)</li></ul> |
| Priority          | Medium  |
| Ventilator action | The ventilator will turn off the humidifier and continue treatment with the same settings.<br>The humidifier must be restarted manually when the cause of the alarm is resolved.  |

6.4.26 Heated Circuit Fault Alarm

| Property          | Description   |
|-------------------|---|
| Alarm text        | <b>Heated Circuit Fault</b>   |
| Alarm condition   | A Heated Circuit Fault alarm will be given if a fault in the heated circuit electronics or temperature sensor is detected.  |
| Priority          | Medium  |
| Ventilator action | The ventilator will turn off the heated circuit and continue treatment with the same settings.<br>The heated circuit must be restarted manually when the cause of the alarm is resolved.            |
| Reset             | The alarm is dismissed when the heated circuit setting is changed to OFF, or the treatment is stopped.<br>The power to the heated circuit is re-enabled when all enabling conditions are satisfied. |

#### 6.4.27 Internal Function Failure

| Property          | Description  |
|-------------------|--|
| Alarm text        | <b>Int. Function Failure</b>   |
| Priority          | High   |
| Alarm condition   | Failure of internal function that prevents treatment or normal operation of the ventilator.<br>The error code that follows the alarm text indicates the kind of function failure.<br>All Internal Function Failure alarm error codes are defined and explained in the ventilator Service Manual. |
| Reset criteria    | Restart the ventilator.  |
| Ventilator action | The ventilator will stop the treatment and shut down.  |
| Action to take    | Restart the Vivo 45 LS.<br>If the alarm persists or reoccurs: Take a note of the error code and contact your supplier of the Vivo 45 LS .  |

#### 6.4.28 Air Temp. Sensor Fail Alarm

| Property          | Description  |
|-------------------|--|
| Alarm text        | <b>Air Temp Sensor Fail</b>  |
| Alarm condition   | The alarm is given in case of swivel boot temperature sensor communication failure or sensor reporting temperatures out of range (below -30°C (-22°F) or above 60°C (140°F). |
| Priority          | Medium   |
| Ventilator action | The ventilator will continue treatment with the same settings.   |

#### 6.4.29 Internal Error Alarm

| Property          | Description  |
|-------------------|--|
| Alarm text        | <b>Internal Error</b>  |
| Priority          | High   |
| Alarm Condition   | An internal Error alarm will be given when the ventilator has an internal error, followed by an error code for the specific failure.<br>All Internal Error alarm error codes are defined and explained in the ventilator Service Manual. |
| Ventilator action | The ventilator will continue treatment with the same settings.   |
| Reset action      | Power off and restart the ventilator.  |



#### 6.4.30 Database Integrity Fail Alarm

| Property          | Description  |
|-------------------|--|
| Alarm text        | <b>Database Integrity Failed</b>                               |
| Priority          | High   |
| Alarm Condition   | This alarm is given when the database integrity check fails.   |
| Ventilator action | The ventilator will continue treatment with the same settings. |
| Reset action      | Rebuild the database and restart the ventilator.               |

#### 6.4.31 Cooling Fan Error Alarm

| Property          | Description  |
|-------------------|--|
| Alarm text        | <b>Cooling Fan Error</b>   |
| Alarm Condition   | The Cooling Fan Error alarm shall be given when the cooling fan runs too slow. |
| Priority          | High   |
| Ventilator action | The ventilator will continue treatment with the same settings.                 |
| Reset             | When cooling fan speed is above 275 rpm.                                       |

#### 6.4.32 Clock Failure Alarm

| Property          | Description   |
|-------------------|---|
| Alarm text        | <b>Clock Failure</b>  |
| Priority          | High  |
| Alarm condition   | The alarm shall be given when the real time clock value is invalid. |
| Ventilator action | The ventilator will continue treatment with the same settings.      |
| Reset action      | Restart the ventilator.   |

#### 6.4.33 Internal Temp High Alarm

| Property          | Description  |
|-------------------|--|
| Alarm text        | <b>Internal Temp High</b>  |
| Priority          | High   |
| Alarm condition   | The Internal High Temp alarm shall be given when the ventilator internal temperature is high.<br>The internal temp high alarm is triggered when PTU/Sensor board temperature is higher than 65°C (149°F), or main board temperature is higher than 65°C (149°F), or motor temperature is higher than 85°C (185°F). |
| Ventilator action | The ventilator will continue treatment with the same settings.   |
| Reset             | When the triggering conditions are resolved.   |

6.4.34 Humidifier/Bypass Loose Alarm

| Property          | Description   |
|-------------------|---|
| Alarm text        | <b>Humidifier/Bypass Loose</b>  |
| Priority          | Medium  |
| Alarm condition   | The Humidifier/Bypass Loose alarm shall be given when the air bypass/humidifier latch is stuck in the down position for 5 secs. |
| Ventilator action | The ventilator will continue treatment with the same settings.  |
| Reset action      | Reinsert the air bypass unit/humidifier and make sure the latch closes.   |

6.5 Alarm Test

6.5.1 Alarm Signal Test

When starting treatment, an automatic alarm signal test is performed. Check that the test is performed successfully, this is indicated by:

- A short beep indicating functional audio signaling.
- The alarm LED first lights yellow, then red, indicating functional visual signaling.
- The audio pause LED lights yellow.
- In about a second, both LEDs are turned off.

If the test fails, do not use the Vivo 45 LS. Contact your supplier of the Vivo 45 LS for a technical check.

6.5.2 Mandatory Alarm Tests

This alarm test should be performed every 24th month or if the ventilator’s function needs to be checked for any other reason.

The alarm test should be included in the regular inspections during maintenance.

To perform the alarm test, follow the instructions below:

## Alarm Test Preparation

- 1 Connect the ventilator patient circuit to a test lung.
- 2 Connect the ventilator to Mains power supply.
- 3 Start the ventilator.
- 4 Adjust the settings as follows:

| Setting           | Value                              |
|-------------------|------------------------------------|
| Ventilation Mode  | Pressure Support Ventilation (PSV) |
| Patient Mode      | Adult                              |
| Insp. Pressure    | 15 cmH <sub>2</sub> O              |
| PEEP              | 5 cmH <sub>2</sub> O               |
| Rise Time         | 9                                  |
| Insp. Trigger     | 9                                  |
| Exp. Trigger      | 3                                  |
| Min Insp. Time    | Off                                |
| Max Insp. Time    | Off                                |
| Backup Rate       | 12 bpm                             |
| Backup Insp. Time | 2.0 s                              |
| Target Volume     | Off                                |

- 5 All alarm settings shall be set to Off if possible.
- 6 Start the treatment.

### 6.5.2.1 High and Low Pressure Alarm Tests

- 1 Set the high pressure alarm to 10 cmH<sub>2</sub>O.  
⇒ The high pressure alarm shall be given.
- 2 Set the high pressure alarm to 55 cmH<sub>2</sub>O.
- 3 Set the low pressure alarm to 20 cmH<sub>2</sub>O.  
⇒ The low pressure alarm shall be given.
- 4 Set the low pressure alarm to 1.0 cmH<sub>2</sub>O.

### 6.5.2.2 Expiratory Tidal Volume Alarm (V<sub>t<sub>e</sub></sub>) Tests

This alarm test applies if having a patient circuit with intentional leakage.

- 1 Set up the ventilator as described in *Alarm Test Preparation*, page 121.
- 2 Set the high V<sub>t<sub>e</sub></sub> alarm to 150 ml.  
⇒ The high V<sub>t<sub>e</sub></sub> alarm shall be given.
- 3 Set the high V<sub>t<sub>e</sub></sub> alarm to Off.
- 4 Set the low V<sub>t<sub>e</sub></sub> alarm to 400 ml.  
The low V<sub>t<sub>e</sub></sub> alarm shall be given.

### 6.5.2.3 Inspiratory Tidal Volume Alarm ( $V_{Ti}$ ) Tests

This alarm test applies if having a patient circuit with exhalation valve or a patient circuit with mouthpiece.

- 1 Set up the ventilator as described in *Alarm Test Preparation*, page 121.
- 2 Set the high  $V_{Ti}$  alarm to 150 ml.  
⇒ The high  $V_{Ti}$  alarm shall be given.
- 3 Set the high  $V_{Ti}$  alarm to Off.
- 4 Set the low  $V_{Ti}$  alarm to 400 ml.  
The low  $V_{Ti}$  alarm shall be given.

### 6.5.2.4 EtCO<sub>2</sub> Related Alarm Test

This alarm test applies if the EtCO<sub>2</sub> accessory is used.

- 1 Connect the EtCO<sub>2</sub> sensor with an attached airway adapter to the Vivo 45 I.S.
- 2 Disconnect the airway adapter from the CO<sub>2</sub> sensor.  
⇒ The check CO<sub>2</sub> adapter alarm shall be given.
- 3 Connect the airway adapter to the CO<sub>2</sub> sensor again.

### 6.5.2.5 SpO<sub>2</sub> Related Alarm Tests

These tests apply if the SpO<sub>2</sub> accessory is used,

- 1 Connect SpO<sub>2</sub> sensor to device and to your finger.
- 2 Set the low SpO<sub>2</sub> alarm to 85%.
- 3 Set the high SpO<sub>2</sub> alarm to be 90%.
- 4 Start treatment and wait 30 s.  
⇒ High SpO<sub>2</sub> alarm should be given.
- 5 Stop treatment.
- 6 Set the high SpO<sub>2</sub> alarm to off.
- 7 Set the low SpO<sub>2</sub> alarm to be 100%.
- 8 Start treatment and wait 30 s.  
⇒ Low SpO<sub>2</sub> alarm should be given.
- 9 Stop Treatment.
- 10 Set the low SpO<sub>2</sub> alarm to 85%.
- 11 Set the low pulse rate alarm to off.
- 12 Set the high pulse rate alarm to 30 bpm.
- 13 Start treatment and wait 30 s.  
⇒ High pulse rate alarm should be given.
- 14 Stop treatment.
- 15 Set the high pulse rate alarm to off.
- 16 Set the low pulse rate alarm to be 230 bpm.
- 17 Start treatment and wait 30 s.  
⇒ Low pulse rate alarm should be given.
- 18 Stop Treatment.
- 19 Set the low pulse rate alarm to off.

**6.5.3 Optional Alarm Tests**

In this chapter, methods for additional alarm tests are described. These tests are optional and not needed to ensure safe use of the ventilator.

**6.5.3.1 High PEEP Alarm**

- 1 Connect the ventilator patient circuit to a test lung and a CPAP device.
- 2 Set the CPAP device treatment pressure to 10 cmH<sub>2</sub>O.
- 3 Adjust the ventilator settings as follows:

| Setting          | Value                 |
|------------------|-----------------------|
| Ventilation Mode | Pressure              |
| Breath Mode      | Assist/Control        |
| Insp. Pressure   | 15 cmH <sub>2</sub> O |
| PEEP             | 5 cmH <sub>2</sub> O  |
| Breath Rate      | 12 bpm                |
| Insp. Time       | 1.5 s                 |
| Rise Time        | 5                     |
| Insp. Trigger    | Off                   |
| Target Volume    | Off                   |

- 4 Start treatment on both the ventilator and the CPAP device.
- 5 Wait approximately 15 seconds before the High PEEP alarm shall be given.
- 6 Stop treatment. Test completed.

**6.5.3.2 Low Pressure and Disconnection Alarms**

- 1 Start treatment and disconnect the patient circuit.
- 2 Wait 15 seconds.
- 3 The Low Pressure Alarm and/or the Disconnection Alarm will be given.
- 4 Stop treatment. Test completed.

**6.5.3.3 Disconnection Alarm Test**

- 1 Set the disconnection alarm to On.
- 2 Disconnect the patient circuit.  
⇒ The disconnection alarm shall be given.
- 3 Set the disconnection alarm to Off.

**6.5.3.4 Obstruction Alarm**

- 1 Start treatment; block the patient circuit completely to simulate an obstruction.
- 2 Wait approximately 10 seconds.
- 3 The Obstruction Alarm will be given.
- 4 Stop treatment. Test completed.

## 7 Cleaning and Maintenance



### WARNING!

The Vivo 45 LS should be subjected to maintenance, service and control and any applicable upgrades, in accordance with Breas service instructions.



The Vivo 45 LS shall only be repaired or modified in accordance with Breas service manuals, technical bulletins, and any special service instructions, by service technicians that have been authorised after Breas Vivo 45 LS service training.



Do not under any circumstances attempt to service or repair the ventilator yourself. If you do so, the manufacturer will no longer be responsible for the performance and safety of the ventilator.

Deviation from these service instructions may lead to risk of personal injury!

The patient-connected parts and the filter must be cleaned and replaced regularly to ensure correct function of the ventilator. All replaced parts must be disposed of in accordance with local environmental regulations regarding the disposal of used equipment and waste.

### 7.1 Cleaning the Vivo 45 LS



#### WARNING!

To avoid electrical shock, disconnect the power supply to the ventilator before cleaning. Do not immerse the ventilator into any fluids.



#### CAUTION!

Always be careful when cleaning to ensure that you do not damage any equipment.



Fluid must not be allowed to enter the ventilator.



Never apply any liquids directly on the ventilator by spraying, splashing or pouring. Use a moistened lint-free cloth when cleaning.



Do not use an excessive amount of liquid when cleaning the ventilator.



Do not autoclave the ventilator.

#### 7.1.1 Main Unit

- 1 Switch off the Vivo 45 LS and disconnect the power supply.
- 2 Remove the patient circuit.
- 3 Disconnect all electric cables.
- 4 Clean the outside of the Vivo 45 LS using a lint-free cloth with a mild soap solution, and/or ethanol 70% for surface disinfection.
- 5 If the click-in humidifier is used, clean it as described in 5.7.11.6 *Cleaning the Water Chamber*, page 88.
- 6 Reconnect the patient circuit. Make sure all parts are dry before the ventilator is put into operation.

7.1.2 Air Pathway Disinfection

The table below lists the parts that might get contaminated by exhaled gases or bodily fluids during normal use or single fault condition.

| Condition               | Parts   |
|-------------------------|---|
| With bacteria filter    | <ul style="list-style-type: none"><li>• Patient circuit</li><li>• EtCO<sub>2</sub> airway adapter (if used)</li><li>• Bacteria filter</li></ul>   |
| Without bacteria filter | <ul style="list-style-type: none"><li>• Patient circuit</li><li>• EtCO<sub>2</sub> airway adapter (if used)</li><li>• FiO<sub>2</sub> sensor (if used)</li><li>• Patient air outlet/Pneumatic unit</li><li>• Air bypass unit/water chamber</li><li>• Blower/Inlet silencer</li><li>• Air inlet with filters</li></ul> |

In case of contamination, the internal air pathways of the Vivo 45 LS may be disinfected up to 5 times by a maximum 60 minute long validated ozone gas process.

Low resistance bacteria filter, if used, should be replaced every 24 hours.

7.1.3 Patient Circuit



The patient circuit should be cleaned and replaced in accordance with the manufacturer's instructions and care provider's instructions, where applicable. For safety information, read 2.4 *Usage of Patient Circuit*, page 16.

Always clean the patient circuit or use a new patient circuit when the ventilator is used by a new patient.

Check the patient circuit regularly for damage. In case of damage, replace the circuit

**CAUTION!**



Appropriate personnel should determine the duration of use for the patient circuit based on accepted infection control procedures.

## 7.2 Cleaning and Replacing the Filters

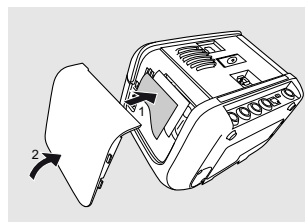
### Patient air filters

#### NOTE

- Coarse filter (grey): This is a washable filter, wash the filter at least once a week and replace once a year. See 7.2.1 *Washing a coarse filter*, page 127 for washing instructions.
- Fine filter (white): This is a disposable filter that not shall be washed or reused. Replace the fine filter at least every month, or more frequently when used in environments rich of pollen, pet hair or other particles that might clog the filter.

The patient air filters are located in the filter cassette at the side of the ventilator.

- 1 Turn off the ventilator and place it on a dust free surface.
- 2 First Place the filters in the air inlet compartment, with the coarse filter outside the fine filter.



- 3 Close the side panel carefully for not displacing the filters while closing. For detailed information about closing the side panel, see 3.3.1 *Detaching and Reattaching the Side Panels*, page 27.


### Cooling Air Filter
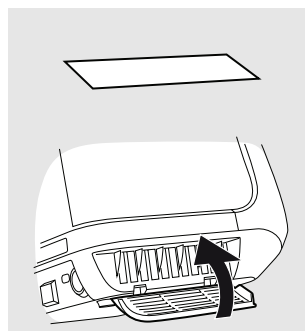
#### NOTE

The filter shall be washed at least once a week and replaced every second year. See 7.2.1 *Washing a coarse filter*, page 127 for washing instructions.

The cooling air inlet filter is located at the bottom left side of the ventilator.

- 1 Open the cooling air filter compartment by pulling at the top of the lid.



- 2 Remove the filter and wash or replace it.
- 3 Put back the filter and close the lid.



### 7.2.1 Washing a coarse filter

- 1 Wash the filter using warm water and a mild soap.
- 2 Rinse thoroughly.
- 3 Dry the filter by squeezing it out in a towel. Do not wring the filter.
- 4 Make sure the filter is completely dry before inserting.

### 7.3 Change of Patients

If the ventilator is used in a clinic by several patients, a low resistance bacterial filter may be used between the air outlet and the patient tube to prevent patient cross-contamination.

- 1 Follow the instructions in 7.1.1 *Main Unit*, page 124, steps 1 to 5.
- 2 Replace the patient filters according to 7.2 *Cleaning and Replacing the Filters*, page 126.
- 3 If a low resistance bacterial filter is used, it shall be replaced. To avoid cross-contamination when no bacterial filter has been used, a validated ozone-disinfection process may be used, see the section on disinfecting the main unit internally.
- 4 Use a new patient circuit when the ventilator is used by a new patient.

### 7.4 Regular Maintenance

Regular maintenance inspections and checks shall be carried out at least every 24 months, according to the ventilator Service Manual.

#### **WARNING!**

Do not use the device and contact your responsible care provider for an inspection of the device in the event of:

- Unexpected patient symptoms during treatment.
- Unexplainable or sudden pressure, performance or sound changes during operation.
- Suspected damage to the device, including the occurrence of Internal Functional Failure alarms.
- Suspected damage to the click-in battery, including evidence of battery cell leakage.



### 7.5 Service and Repair

The service and repair of the ventilator must only be carried out by authorised service personnel in accordance with Breas service instructions. Service inspections must always be carried out following any repairs to the device.

#### **WARNING!**

Do not use the device and contact your responsible care provider for an inspection of the device in the event of:

- Unexpected patient symptoms during treatment.
- Unexplainable or sudden pressure, performance or sound changes during operation.
- Suspected damage to the device, including the occurrence of Internal Functional Failure alarms.
- Suspected damage to the click-in battery, including evidence of battery cell leakage.





Authorised service workshops can order the ventilator Service Manual that contains all technical documentation required for the maintenance and service of the ventilator.

## 7.6 Storage

Store the ventilator in a dark room, where the temperature range is within  $-20$  to  $+60^{\circ}\text{C}$  ( $-4$  to  $+140^{\circ}\text{F}$ ).

### CAUTION!

The ventilator must not be stored in a warm place, such as direct sunlight or close to a radiator. The time required for the device to cool from the maximum storage temperature of  $+60^{\circ}\text{C}$  ( $+140^{\circ}\text{F}$ ) until it is ready for use in ambient temperature of  $+20^{\circ}\text{C}$  ( $+68^{\circ}\text{F}$ ) is 30 minutes.

If stored in a cold environment, let the ventilator adapt to room temperature before using the device. The time required for the device to warm from the minimum storage temperature of  $-20^{\circ}\text{C}$  ( $-4^{\circ}\text{F}$ ) until it is ready for use in ambient temperature of  $+20^{\circ}\text{C}$  ( $+68^{\circ}\text{F}$ ) is 30 minutes.



## 7.7 Disposal

The ventilator, any accessories and all replaced parts must be disposed of and recycled in accordance with the local environmental regulations regarding the disposal of used equipment and waste. Contact your service provider for information regarding the disposal procedure.

### NOTE



Batteries used with the ventilator shall be recycled in accordance with the local environmental regulations.




## 8 Technical Specifications

### 8.1 System Description

#### Active Exhalation Valve Configuration


This diagram provides an overview of the ventilator system when used with an active exhalation valve patient circuit.

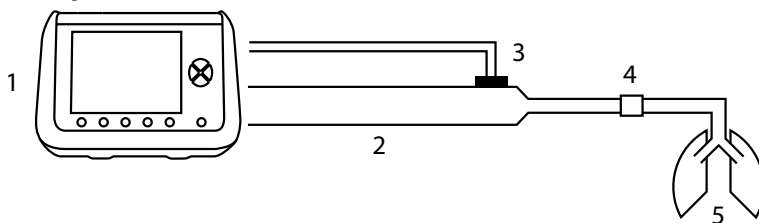

1. Vivo 45 LS
2. Tube
3. Active Exhalation valve
4. Patient interface connection
5. Patient


#### Leakage Port Configuration

This diagram provides an overview of the ventilator system when used with a leakage port patient circuit.

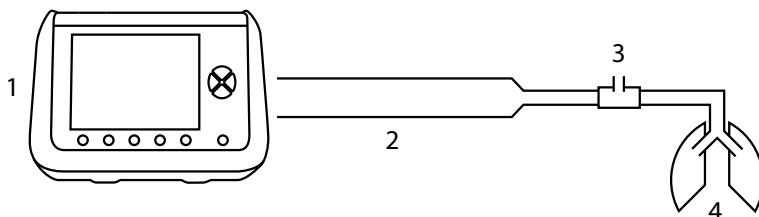

1. Vivo 45 LS
2. Tube
3. Leakage port / Patient interface connection
4. Patient

### 8.1.1 Pneumatic Diagram for the ventilator

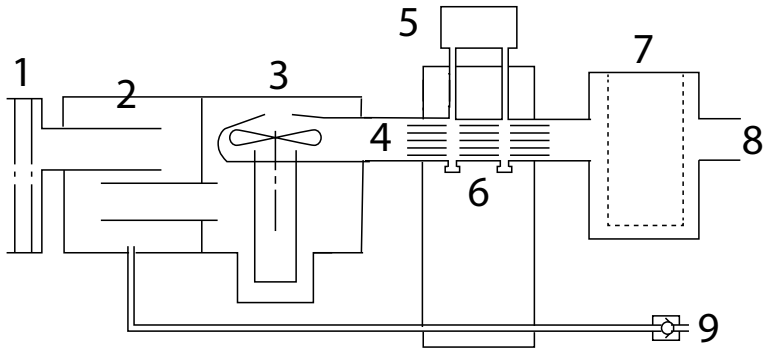

## 8.2 Data

### 8.2.1 Worst Case Accuracy

#### Pressure Control Modes

The worst case Vivo 45 LS configuration is the 15 mm patient circuit with HCH humidifier, bacterial filter and EtCO<sub>2</sub> sensor.

#### Volume control Modes

The worst case Vivo 45 LS configuration is the 15mm circuit with or without HCH humidifier, bacterial filter, FiO<sub>2</sub> sensor and EtCO<sub>2</sub> sensor.

### 8.2.2 Modes Specifications

This section lists the settings that can be made for the ventilator's modes.

#### Ventilation modes

- PSV (Pressure Support Ventilation). May be combined with:
  - TgV (Target Volume)
  - AE (Auto-EPAP)
- PCV (Pressure Controlled Ventilation). May be combined with:
  - A (Assist function — inspiration trigger)
  - TgV (Target Volume)
  - AE (Auto-EPAP)
  - MPV (Mouthpiece Ventilation)
  - SIMV (Synchronized Intermittent Mandatory Ventilation)
- VCV (Volume Controlled Ventilation). May be combined with:
  - A (Assist function — inspiration trigger)
  - MPV (Mouthpiece Ventilation)
  - SIMV (Synchronized Intermittent Mandatory Ventilation)
- CPAP (Continuous Positive Airway Pressure)

### **Device modes**

- Clinical
- Home

### **Patient modes**

- Adult
- Paediatric

## **8.2.3 Parameter Specifications**

This section lists the characteristics for the ventilator's parameters.

All stated tolerances includes measurement uncertainty. The accuracies have been tested with all allowed configurations. Stated tolerances only disclose the maximum tolerance. If a parameter's tolerance is described with both absolute and relative measures, the greater one applies.

| Setting                   | Unit               | Min   | Max  | Default   | Resolution           | Tolerance                            |
|---------------------------|--------------------|---|--|---|----------------------|--------------------------------------|
| Insp. Pressure            | cmH <sub>2</sub> O | 4   | 50   | 15  | 0.5 < 10<br>1.0 ≥ 10 | ±0.5 cmH <sub>2</sub> O or ±5%       |
| CPAP                      | cmH <sub>2</sub> O | 4   | 20   | 10 <sup>(A)</sup><br>8 <sup>(P)</sup>             | 0.5 < 10<br>1.0 ≥ 10 | ±0.5 cmH <sub>2</sub> O or ±5%       |
| PEEP                      | cmH <sub>2</sub> O | 2<br>Off <sup>(2)</sup>                         | 20 <sup>(3)</sup>                                  | 5 <sup>(A)</sup><br>2 <sup>(P)</sup>              | 0.5 < 10<br>1.0 ≥ 10 | ±0.5 cmH <sub>2</sub> O or ±5%       |
| Breath Rate               | bpm <sup>(4)</sup> | 4 <sup>(A)</sup><br>6 <sup>(P)</sup>            | 40 <sup>(A)</sup><br>60 <sup>(P)</sup>             | 12 <sup>(A)</sup><br>20 <sup>(P)</sup>            | 1                    | ±2%                                  |
| Backup Rate               | bpm <sup>(4)</sup> | 4 <sup>(A)</sup><br>6 <sup>(P)</sup><br>0 (MPV) | 40 <sup>(A)</sup><br>60 <sup>(P)</sup><br>40 (MPV) | 12 <sup>(A)</sup><br>20 <sup>(P)</sup><br>0 (MPV) | 1                    | ±2%                                  |
| Backup Insp. Time         | s                  | 0.3   | 5 <sup>(A)</sup><br>2 <sup>(P)</sup>               | 1.5 <sup>(A)</sup><br>1 <sup>(P)</sup>            | 0.1                  | ± (20 ms + 5% of setting) or ±0.1 s, |
| Inspiratory Time          | s                  | 0.3   | 5 <sup>(A)</sup><br>2 <sup>(P)</sup>               | 1.5 <sup>(A)</sup><br>1 <sup>(P)</sup>            | 0.1                  | ± (20 ms + 5% of setting) or ±0.1 s, |
| Min Insp. Time            | s                  | 0.3<br>Off                                      | 3 <sup>(A)</sup><br>2 <sup>(P)</sup>               | Off   | 0.1                  | ± (20 ms + 5% of setting) or ±0.1 s, |
| Max Insp. Time            | s                  | 0.3   | 5 <sup>(A)</sup><br>2 <sup>(P)</sup><br>Off        | Off   | 0.1                  | ± (20 ms + 5% of setting) or ±0.1 s, |
| Inspiratory Trigger       | Step               | 1   | 9<br>Off <sup>(5)</sup>                            | 3 <sup>(A)</sup><br>2 <sup>(P)</sup>              | 1                    | -                                    |
| Expiratory Trigger        | Step               | 1   | 9 <sup>(5)</sup>                                   | 3   | 1                    | -                                    |
| Rise Time (Pressure mode) | Step               | 1   | 9  | 3   | 1                    | -                                    |
| Rise Time (Volume mode)   | s                  | 50% of set <i>Insp Time.</i>                    | 90% of set <i>Insp Time.</i><br>Off                | Off   | 10%                  | ± (20 ms + 5% of setting) or ±0.1 s, |

(A)= Adult mode, (P)= Paediatric mode, (1)= 0.5 cmH<sub>2</sub>O, 1.0 cmH<sub>2</sub>O, (2)= Only with active exhalation valve circuit, (3)= For pressure modes, also limited by Insp. Pressure- 2 cmH<sub>2</sub>O, (4)= breath per minute, (5)= Off is only available in Control/Assist mode, 6) Values >10 might initially be set with a 0.5 step when turning on Auto-EPAP. When changing the value, whole numbers will be used.

| Setting               | Unit               | Min   | Max                           | Default            | Resolution               | Tolerance                            |
|-----------------------|--------------------|---|-------------------------------|--------------------|--------------------------|--------------------------------------|
| Max Pressure          | cmH <sub>2</sub> O | Current <i>Min Pressure</i>                                       | 50                            | 15                 | 0.5 < 10<br>1.0 ≥ 10     | ±0.5 cmH <sub>2</sub> O or ±5%       |
| Min Pressure          | cmH <sub>2</sub> O | 4   | Current <i>Max Pressure</i>   | 15                 | 0.5 < 10<br>1.0 ≥ 10     | ±0.5 cmH <sub>2</sub> O or ±5%       |
| Target Volume         | ml                 | Off<br>300 (A)<br>50 (P)  | 2000 (A)<br>500 (P)           | Off                | 10 < 500<br>50 ≥ 500     | ±12 ml or ±10%,                      |
| Tidal Volume          | ml                 | 300 (A)<br>50 (P)   | 2000 (A)<br>500 (P)           | 500 (A)<br>150 (P) | 10 < 500<br>50 ≥ 500     | ±12 ml or ±10%,                      |
| Sigh                  | -                  | Off   | On                            | Off                | -                        | -                                    |
| Sigh Rate             | 1/Breath           | 10  | 250                           | 50                 | 10                       | -                                    |
| Sigh Inspiratory Time | s                  | Current <i>Inspiratory Time</i> or <i>Backup Inspiratory Time</i> | 5 (A)<br>2 (P)                | 1.5 (A)<br>1 (P)   | 0.1                      | ± (20 ms + 5% of setting) or ±0.1 s, |
| SIMV Rate             | bpm <sup>(4)</sup> | 4 (A)<br>6 (P)  | 40 (A)<br>60 (P)              | 12 (A)<br>20 (P)   | 1                        | ±2%                                  |
| SIMV Support Pressure | cmH <sub>2</sub> O | Current PEEP +2   | 50                            | 15                 | 0.5 < 10<br>1.0 ≥ 10     | ±0.5 cmH <sub>2</sub> O or ±5%       |
| Auto-EPAP             | -                  | Off   | On                            | Off                | -                        | -                                    |
| EPAP Min              | cmH <sub>2</sub> O | 2   | 20 or Current <i>EPAP Max</i> | 5                  | 0.5 < 10<br>1.0 ≥ 10     | ±0.5 cmH <sub>2</sub> O or ±5%       |
| EPAP Max              | cmH <sub>2</sub> O | 2 or Current <i>EPAP Min</i>                                      | 20 or <i>Pressure Limit-2</i> | 5                  | 0.5 < 10<br>1.0 ≥ 10     | ±0.5 cmH <sub>2</sub> O or ±5%       |
| EPAP Step             | cmH <sub>2</sub> O | 0.5   | 2                             | 1                  | 0.5                      | ±0.5 cmH <sub>2</sub> O or ±5%       |
| PS                    | cmH <sub>2</sub> O | 2   | 50 - Current <i>EPAP Max</i>  | 10                 | 0.5 < 10<br>1.0 ≥ 10 (6) | ±0.5 cmH <sub>2</sub> O or ±5%       |

(A)= Adult mode, (P)= Paediatric mode, (1)= 0.5 cmH<sub>2</sub>O < 10 cmH<sub>2</sub>O, 1.0 cmH<sub>2</sub>O ≥ 10 cmH<sub>2</sub>O, (2)= Only with active exhalation valve circuit, (3)= For pressure modes, also limited by Insp. Pressure- 2 cmH<sub>2</sub>O, (4)= breath per minute, (5)= *Off* is only available in Control/Assist mode, (6) Values >10 might initially be set with a 0.5 step when turning on Auto-EPAP. When changing the value, whole numbers will be used.

| Setting             | Unit               | Min   | Max                             | Default                           | Resolution                          | Tolerance                            |
|---------------------|--------------------|---|---------------------------------|-----------------------------------|-------------------------------------|--------------------------------------|
| Min PS              | cmH <sub>2</sub> O | 2   | 50 - Current<br><i>EPAP Max</i> | Variable <sup>(1)</sup>           | 0.5 < 10<br>1.0 ≥ 10 <sup>(6)</sup> | ±0.5 cmH <sub>2</sub> O or ±5%       |
| Max PS              | cmH <sub>2</sub> O | 2   | 50 - Current<br><i>EPAP Max</i> | 10                                | 0.5 < 10<br>1.0 ≥ 10 <sup>(6)</sup> | ±0.5 cmH <sub>2</sub> O or ±5%       |
| Pressure Limit      | cmH <sub>2</sub> O | Current <i>EPAP</i><br><i>Max</i> +2                                | 50                              | <i>High Pressure</i><br>Alarm - 2 | 0.5 < 10<br>1.0 ≥ 10                | ±0.5 cmH <sub>2</sub> O or ±5%       |
| Relax Time          | Minute             | 2   | 12<br>Off                       | 5                                 | 1                                   | ± (20 ms + 5% of setting) or ±0.1 s, |
| Humidifier Setting  | Step               | 1   | 5                               | 3                                 | 1                                   | -                                    |
| Heated Circuit Temp | °C/°F              | 16/61   | 30/86                           | 27/81                             | 0.5                                 |                                      |
| Flow Pattern        | -                  | Square (constant flow)<br>Decelerating (Flow is decreased linearly) |                                 |                                   | Square                              | -                                    |
| Audible Alarm Level | Step               | 1   | 5                               | 3                                 | 1                                   | -                                    |

(A)= Adult mode, (P)= Paediatric mode, (1)= 0.5 cmH<sub>2</sub>O < 10 cmH<sub>2</sub>O, 1.0 cmH<sub>2</sub>O ≥ 10 cmH<sub>2</sub>O, (2)= Only with active exhalation valve circuit, (3)= For pressure modes, also limited by Insp. Pressure- 2 cmH<sub>2</sub>O, (4)= breath per minute, (5)= *Off* is only available in Control/Assist mode, 6) Values >10 might initially be set with a 0.5 step when turning on Auto-EPAP. When changing the value, whole numbers will be used.



## 8.2.4 Monitored Values Specifications

This section describes the ranges and tolerances for monitored values on the Vivo 45 LS. All stated tolerances include measurement uncertainty. The accuracies have been tested with all allowed configurations. Stated tolerances only disclose the maximum tolerance.

### **P<sub>peak</sub>**

**Range/Performance:** 4 to 99 cmH<sub>2</sub>O.

**Tolerance:**  $\pm 0.5$  cmH<sub>2</sub>O or  $\pm 10\%$ , whichever is greatest

### **PEEP**

**Range/Performance:** 0 to 99 cmH<sub>2</sub>O.

**Tolerance:**  $\pm 0.5$  cmH<sub>2</sub>O or  $\pm 10\%$ , whichever is greatest

### **P<sub>mean</sub>**

**Range/Performance:** 0 to 99 cmH<sub>2</sub>O.

**Tolerance:**  $\pm 0.5$  cmH<sub>2</sub>O or  $\pm 10\%$ , whichever is greatest

### **CPAP<sub>pressure</sub>**

**Range/Performance:** 0 to 99 cmH<sub>2</sub>O.

**Tolerance:**  $\pm$  (4% CPAP set pressure + 0.8 cmH<sub>2</sub>O)

### **Leakage**

**Range/Performance:** 0 to 99.9 l/min (BTPS\*).

**Tolerance:**  $\pm 10\%$

### **MV<sub>i</sub>**

**Range/Performance:** 0 to 99.9 l (BTPS\*).

**Tolerance:**  $\pm 10\%$  or  $(\pm 10 \text{ ml} \times \text{bpm})$ , whichever is greatest

### **MV<sub>e</sub>**

**Range/Performance:** 0 to 99.9 l (BTPS\*).

**Tolerance:**  $\pm 10\%$  or  $(\pm 10 \text{ ml} \times \text{bpm})$ , whichever is greatest

### **Vt<sub>i</sub>**

**Range/Performance:** 0 to 9999 ml (BTPS\*).

**Tolerance:**  $\pm 10$  ml or 10%, whichever is greatest

### **Vt<sub>e</sub>**

**Range/Performance:** 0 to 9999 ml (BTPS\*).

**Tolerance:**  $\pm 10$  ml or 10%, whichever is greatest

### **FiO<sub>2</sub>**

**Range/Performance:** 0 to 100%.

**Tolerance:**  $\pm 2\%$

### **% in TgV**

**Range/Performance:** 0 to 100%.

**Tolerance:**  $\pm 1\%$

**Total Rate**

**Range/Performance:** 0 to 99 bpm.

**Tolerance:**  $\pm 1$  bpm

**Spont Rate**

**Range/Performance:** 0 to 99 bpm.

**Tolerance:**  $\pm 1$  bpm

**% Spont**

**Range/Performance:** 0 to 100%.

**SpO<sub>2</sub>**

**Range/Performance:** 70 to 100%.

**Tolerance:**  $\pm 3$  digits. No motion and flex sensor.

**Pulse Rate**

**Range/Performance:** 25 to 240 bpm.

**Tolerance:**  $\pm 3$  digits. No motion and flex sensor.

**I:E**

**Range/Performance:** 1:10 to 10:1.

**Tolerance:**  $\pm 0.1$  unit for I:E < 9.9,  $\pm 1$  unit otherwise.

**Insp. Time**

**Range/Performance:** 0.3 to 5 s.

**Tolerance:**  $\pm 0.1$  s

**Rise Time**

**Range/Performance:** 0.1 to 5 s.

**Tolerance:**  $\pm 10\%$  or  $\pm 0.1$  s, whichever is greatest

**EtCO<sub>2</sub>**

**Range/Performance:** 0 to 25%.

**Tolerance:** 0 to 15%:  $\pm (0.2 \text{ vol\%} + 2\% \text{ of reading})$ . 15 to 25%: unspecified

**InspCO<sub>2</sub>**

**Range/Performance:** 0 to 25%.

**Tolerance:** 0 to 15%:  $\pm (0.3 \text{ vol\%} + 4\% \text{ of reading})$ . 15 to 25%: unspecified

**8.2.5 Power Supply**

**AC supply:** 100 to 240 V AC, tolerance:  $+10\%/-20\%$ , 50 to 60 Hz, 1.0 - 2.0 A.

**External DC:** 19 V DC, tolerance:  $19 \text{ V} \pm 6 \text{ V}$ . Max 90 W.

**Click-in battery:** Capacity: 65Wh. Li-ion.

**Internal battery:** Capacity: 25Wh. Li-ion. Expected service life: 500 full charging cycles.

**8.2.6 Environmental Conditions**

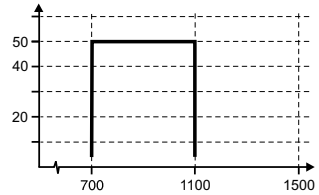
**Operating temperature range:** 5 to 40°C (41 to 104°F)

**Storage and transport temperature:** -20 to +60°C (-4 to +140°F)

**Ambient pressure range:**

700 to 1100 mbar, corresponding to ~4200 metres (13800 feet) above sea level to ~700 metres (2300 feet) below sea level, at normal atmospheric pressure.

As seen in the graph above, the ventilator is unable to deliver set max pressure at a very low ambient pressure.

**Ingress Protection:**

IP22

Solid particle protection: Hazardous parts are protected from touch by fingers and by objects greater than 12 mm.

Liquid ingress protection: The protection withstands dripping water less than 15 degrees from vertical.

The ingress protection has been tested by water drips equivalent to 3mm rain/minute for 10 minutes (2.5 minutes for each tilting direction).

**8.2.7****Other****Patient Circuit Leakage**

**Recommended leakage:** 20 to 50 l/min at 10 cmH<sub>2</sub>O (leakage circuit)

**Minimum leakage:** 12 l/min at 4 cmH<sub>2</sub>O (leakage circuit)

**Oxygen Inlet**

**Oxygen inlet port:** Maximum flow: 30 l/min (medical oxygen). Oxygen coupling is type CPC PMCD181032.

**Start-up Time**

**Start-up from unpowered state:** about 20 seconds.

**Sound Power Level**

**Sound level at 10 cmH<sub>2</sub>O in CPAP mode:** Less than 30 dB(A). Measured at 1 m.

**Alarm sound level :** Adjustable 50–80 dB(A), Measured at 1m. Tolerance:  $\pm 5$  dB(A).

**Miscellaneous**

**Maximum flow:** > 300 l/min

**Maximum flow at 20 mbar:** > 150 l/min

**Maximum limited pressure during single fault condition:** 80 cmH<sub>2</sub>O (PCV, PSV & VCV) 30 cmH<sub>2</sub>O (CPAP)

**Breathing resistance under single-fault:** <6 cmH<sub>2</sub>O at 30 l/min, <6 cmH<sub>2</sub>O at 60 l/min

**Bias-flow when using active exhalation valve:** 8 l/min

**Vivo 45 LS Dimensions**

**W × H × D:** 216 × 159 × 152 mm

**Weight:** 2.4 kg

**Patient air outlet:** 22 mm male, conical standard connector

### **EtCO<sub>2</sub> Sensor**

**W × H × D:** 38 × 37 × 34 mm

**Cable length:** 2.4 m

**Weight:** 75 g

**Warm-up time:** 10 s

**Total system response time:** 30 s

**Interference from medical gases: O<sub>2</sub>:** <-0.1% relative CO<sub>2</sub> per % O<sub>2</sub>  
(calibrated at 21% O<sub>2</sub>)

### **FiO<sub>2</sub> Sensor**

**Total system response time :** 20 s

### **Filtering/Smoothing Techniques**

**Pressure:** Low pass average time constant 16 ms

**Inspiration trigger:** Differential mass flow resolution 4 ms

**Expiration trigger:** Flow low pass filtering with level sensing

**SpO<sub>2</sub>:** No data post-processing done by the ventilator

**Effort Belt:** Low pass filter: 5Hz, High pass filter: 0.1Hz

## **8.3 Emission and Immunity Declaration**

According to IEC 60601-1-2:2014.

The performance of all functions of the ventilator is considered as essential performance for the purpose of immunity testing.

### **8.3.1 Vivo 45 LS Essential Performance**

The ventilator will deliver ventilation at the patient-connection port within its published accuracy specifications and within the alarm limits set by the operator, or generate an alarm condition for high pressure, low pressure, high PEEP, low tidal volume, low minute volume, low breath rate, high EtCO<sub>2</sub>, high and low FiO<sub>2</sub>, obstruction, low last power source, or power failure.

The ventilator will provide SpO<sub>2</sub> and pulse rate values within its published accuracy specifications and generate an alarm upon a low SpO<sub>2</sub> condition. The ventilator will provide indication when the SpO<sub>2</sub> value or pulse rate is potentially incorrect, and generate an alarm condition to indicate when the SpO<sub>2</sub> value update period has exceeded 30 seconds.

The ventilator will provide EtCO<sub>2</sub> and FiO<sub>2</sub> values within its published accuracy specifications and generate an alarm condition upon high and low EtCO<sub>2</sub> and FiO<sub>2</sub> conditions.

Under the immunity test conditions, the following allowances are acceptable:

- Error of delivered volume and PEEP of individual breaths up to 35% and error of the delivered volume and PEEP averaged over a one-minute interval up to 25%.
- Any temporary degradation of SpO<sub>2</sub>, EtCO<sub>2</sub> or FiO<sub>2</sub> performance following transient immunity test exposure shall recover from any disruption within 30 seconds.

Additionally, the following shall not be allowed:

- permanent damage or unrecoverable loss of function,
- changes in programmable parameters or settings,
- reset to default settings,
- change of operating mode,
- initiation of unintended operation.

### 8.3.2 Guidance and Manufacturer's Declaration – Electromagnetic Immunity

The ventilator is intended for use in the electromagnetic environment specified below. The customer or the user of the ventilator should assure that it is used in such an environment.

| Immunity Test   | Compliance Level   | Electromagnetic Environment - Guidance  |
|---|--|---|
| Electrostatic discharge (ESD)<br>IEC 61000-4-2  | $\pm 8$ kV contact $\pm 15$ kV air   | The relative humidity should be at least 5 %.   |
| Electrical fast transient/burst<br>IEC 61000-4-4  | $\pm 2$ kV for power supply lines<br>$\pm 1$ kV for input/output lines   | AC Power (Mains) quality should be that of a typical commercial, hospital and residential environment.  |
| Surge IEC 61000-4-5   | <b>Input power ports:</b><br>0.5 and 1 kV (Line to Line)<br>0.5, 1.0 and 2.0 kV (Line to Earth)<br><b>Signal Input/Output ports:</b><br>2.0 kV (Line to Earth) | AC Power (Mains) quality should be that of a typical commercial, hospital and residential environment.  |
| Power frequency (50/60 Hz) magnetic field IEC 61000-4-8   | 30 A/m   | Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial, hospital and residential environment. |
| Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11 | 0% $U_T$ , 0.5 cycle (multiple phase analysis);<br>0% $U_T$ , 1 cycle;<br>70% $U_T$ , 25/30 cycles (50/60 Hz);<br>0% $U_T$ , 250/300 cycles (50/60 Hz);        | Vivo 45 LS runs on internal battery during voltage dips, short interruptions and voltage variations on power supply input lines.                        |




$U_T$  is the AC Power (Mains) voltage prior to application of the test level.



#### WARNING!

Portable RF communications equipment (including peripherals such as antenna cables and external antennas) should be used no closer than 30 cm (12 inches) to any part of the ventilator, including cables specified. Otherwise, degradation of the performance of this equipment could result.

| Immunity Test              | Compliance Level  | Electromagnetic Environment - Guidance  |
|----------------------------|---|---|
| Conducted RF IEC 61000-4-6 | 3 V <sub>rms</sub> (150 kHz to 80 MHz)<br>6 V <sub>rms</sub> (inside ISM/ASR bands) | $d=0.35*\sqrt{P}$ m at 150 kHz to 80 MHz  |
| Radiated RF IEC 61000-4-3  | 10 V/m 80 MHz to 2.7 GHz  | $d= 1.2*\sqrt{P}$ m at 80 MHz to 800 MHz<br>$d= 2.3*\sqrt{P}$ m at 800 MHz to 2.5 GHz<br><b>Equation description:</b> P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and, except for portable RF communications equipment, d is the recommended separation distance in meters (m).<br>Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey <sup>a</sup> , should be less than the compliance level in each frequency range <sup>b</sup> .<br>Interference may occur in the vicinity of equipment marked with this symbol:<br> |

#### NOTE

At 80 MHz and 800 MHz, the higher frequency range applies.

These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

a) Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the ventilator is used exceeds the applicable RF compliance level above, the ventilator should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as reorienting or relocating the ventilator.

b) Over the frequency range 150 kHz to 80 MHz, field strengths should be less than 10 V/m.

### 8.3.3 Guidance and Manufacturer's Declaration – Electromagnetic Emission

The ventilator is intended for use in the electromagnetic environment specified below. The customer or the user of the ventilator should assure that it is used in such an environment.

| Emissions test   | Compliance Level | Electromagnetic Environment - Guidance  |
|--|------------------|---|
| RF emissions<br>CISPR 11                                   | Group 1          | The ventilator uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.   |
| RF emissions<br>CISPR 11                                   | Class B          | The ventilator is suitable for use in all establishments, including domestic establishments and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes. |
| Harmonic emissions<br>IEC 61000-3-2                        | Class A          |   |
| Voltage fluctuations/<br>flicker emission<br>IEC 61000-3-3 | Complies         |   |

### 8.3.4 Frequencies of portable and mobile transmitters for which the recommended separation distance is 30 cm (12 inches)

| Band (MHz)    | Service   | Immunity test level (V/m) |
|---------------|---|---------------------------|
| 380 — 390     | TETRA 400   | 27                        |
| 430 — 470     | GMRS 460, FRS 460   | 28                        |
| 704 — 787     | LTE Band 13, 17   | 9                         |
| 800 — 960     | GSM 800/900, TETRA 800, iDEN 820, CDMA 850, LTE Band 5          | 28                        |
| 1,700 — 1,990 | GSM 1800; CDMA 1900; GSM 1900; DECT; LTE Band 1, 3, 4, 25; UMTS | 28                        |
| 2,400 — 2,570 | Bluetooth, WLAN 802.11 b/g/n, RFID 2450, LTE Band 7             | 28                        |
| 5,100 — 5,800 | WLAN 802.11 a/n   | 9                         |

### 8.3.5 Recommended separation distances between external power conductors and the ventilator

| Rated maximum current in conductor (A) | Separation distance (m)                           |
|--|---|
|  | <b>50-60 Hz <math>d = I/2\pi H = I/188</math></b> |
| 1                                      | 0.005   |
| 10                                     | 0.05  |
| 30                                     | 0.16  |

For conductors rated at a maximum current not listed above, the recommended separation distance  $d$  in metres (m) can be estimated using the equation  $d = I/2\pi H$ , where  $I$  is the maximum current rating of the conductor in amperes (A) according to the transmitter manufacturer;  $H$  is the ventilator immunity compliance level to electromagnetic fields in the 50-60 Hz frequency span (30 A/m).

8.3.6 Proximity fields from RF wireless communications equipment

| Frequency Range and Level: RF wireless communication equipment  |                                     |                      |
|---|-------------------------------------|----------------------|
| Test Frequency (MHz)  | Modulation*                         | Immunity Level (V/m) |
| 385   | Pulse Modulation: 18 Hz             | 27                   |
| 450   | **FM<br>±5 Hz deviation: 1 kHz sine | 28                   |
| 710<br>745<br>780   | Pulse Modulation: 217 Hz            | 9                    |
| 810<br>870<br>930   | Pulse Modulation: 18 Hz             | 28                   |
| 1720<br>1845<br>1970  | Pulse Modulation: 217 Hz            | 28                   |
| 2450  | Pulse Modulation: 217 Hz            | 28                   |
| 5240<br>5500<br>5785  | Pulse Modulation: 217 Hz            | 9                    |
| * The carrier shall be modulated using a 50 % duty cycle square wave signal.<br>** As an alternative to FM modulation, 50 % pulse modulation at 18 Hz may be used as it would be worst case although it does not represent the actual modulation. |                                     |                      |

8.4 Delivery Settings

Delivery settings: modes and functions

- Ventilation Mode: Pressure
- Breath Mode: Support
- Patient Mode: Adult
- Device Mode: Clinical
- Home Adjust: Off
- Profile 1: Active
- Profile 2: Off
- Profile 3: Off



## **Delivery Settings, Alarms**

### **High Pressure Alarm:**

25 cmH<sub>2</sub>O (Adult)

20 cmH<sub>2</sub>O (Paediatric)

**Low Pressure Alarm:** 10 cmH<sub>2</sub>O

**High PEEP Alarm:** Off

**Low PEEP Alarm:** Off

**High V<sub>t</sub> Alarm:** 500 ml (Adult), 400 ml (Paediatric)

**High V<sub>e</sub> Alarm:** 500 ml (Adult), 400 ml (Paediatric)

**Low V<sub>t</sub> Alarm:** 300 ml (Adult), 100 ml (Paediatric)

**Low V<sub>e</sub> Alarm:** 300 ml (Adult), 100 ml (Paediatric)

**High MV<sub>i</sub> Alarm:** Off

**High MV<sub>e</sub> Alarm:** Off

**Low MV<sub>i</sub> Alarm:** Off

**Low MV<sub>e</sub> Alarm:** Off

**High Breath Rate Alarm:** Off

**Low Breath Rate Alarm:** Off

**Apnea Alarm:** Off

**Disconnection Alarm:** On

**Rebreathing Alarm:** On

**Obstruction Alarm:** Off

**High FiO<sub>2</sub> Alarm:** Off

**Low FiO<sub>2</sub> Alarm:** Off

**High SpO<sub>2</sub> Alarm:** Off

**Low SpO<sub>2</sub> Alarm:** Off

**High EtCO<sub>2</sub> Alarm:** 51 mmHg

**Low EtCO<sub>2</sub> Alarm:** Off

**High InspCO<sub>2</sub> Alarm:** Off

**Low Pulse Rate:** Off

**High Pulse Rate:** Off

## **Other**

**Patient operating time:** 0 h

**Display light:** On

**Light Intensity:** 9

**Alarm sound level:** 5

**CO<sub>2</sub> Unit:** mmHg

**Auto keypad lock:** Off

**Pre-use Test:** On

**WARNING!**

Only use accessories recommended by Breas Medical. Breas Medical cannot guarantee the performance and safety for the use of other accessories with the ventilator.

**NOTE**

Accessory equipment connected to the analogue and digital interfaces must be certified according to the respective IEC standards (e.g. IEC 60950 for data processing equipment and IEC 60601-1 for medical equipment). Furthermore all configurations must comply with the valid version of the system standard IEC 60601-1-1. Anybody who connects additional equipment to the signal input part or signal output part is configuring a medical system, and is therefore responsible for ensuring the system complies with the requirements of the valid version of the system standard IEC 60601-1-1. If in doubt, consult the technical service department or your local representative.

The following Breas accessories have been verified for use with the Vivo 45 LS.:

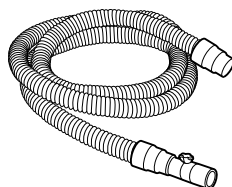
**Not all accessories are available in all markets.**

## 9.1

**Patient Circuits and Air Delivery Accessories****Circuit: 22 mm smoothbore with leak port**


**Function:** Delivers air to the patient, applied part

**Part No:** 005060

**Circuit: 1.8m x 22mm Smoothbore disposable**


**Function:** Delivers air to the patient, applied part

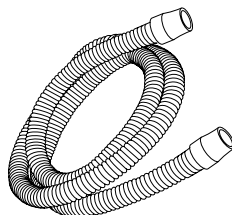**Part No:** 009118

**Circuit: Single limb 22 mm. Single patient multiple use.**

**Function:** Delivers air to the patient, applied part

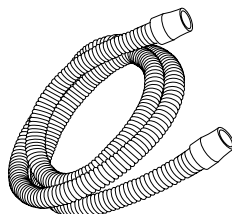
**Part No:** 008426 (30-pack of 004465)



**Circuit: Single limb heated wire 15 mm, disposable**

**Function:** Deliver heated air to the patient, non-invasively

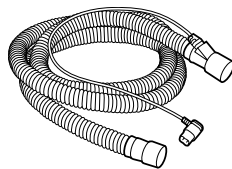
**Part No:** 006193



**Circuit: Single limb 22 mm with exhalation valve, disposable**

**Function:** Delivers air to the patient (applied part)

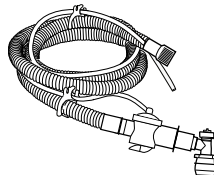
**Part No:** 007387



**Circuit: Single limb 22 mm with leakage port and bacterial filter, disposable**

**Function:** Delivers air to the patient (applied part)

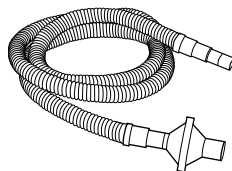
**Part No:** 007615



**Circuit: 1.8m x 15mm Smoothbore Disposable**

**Function:** Deliver air to the patient

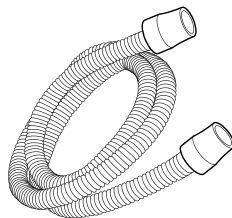
**Part No:** 009119



**Circuit: Single limb 15 mm**

**Function:** Deliver air to the patient

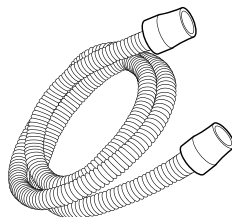
**Part No:** 008427 (30-pack of 006712)



**Circuit: Single limb 22 mm, with bacterial filter, disposable**

**Function:** Deliver air to the patient

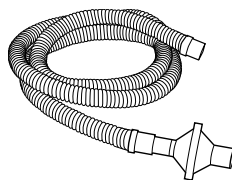
**Part No:** 007936



**Circuit: Dual limb with exhalation valve, disposable**

**Function:** Deliver air to the patient

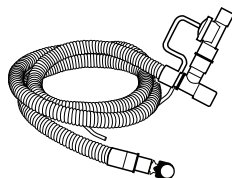
**Part No:** 007616



**Circuit: Smoothbore circuit with exhalation valve, disposable**

**Function:** Deliver air to the patient

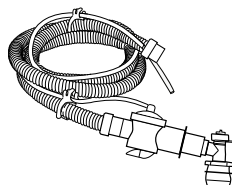**Part No:** 007857



**Pilot line connector**

Connects the exhalation valve control tube to the exhalation valve port of the Vivo 45 LS.

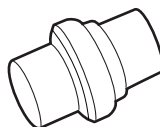
**Part No:** 007654



**Leakage Port**

**Function:** Providing a leakage for clearing exhaled gases.

**Part No:** 004426




### Low resistance bacterial filter, with CO<sub>2</sub> connector

**Function:** Filter air at ventilator outlet

#### Characteristics

- Resistance:  
0.5 cmH<sub>2</sub>O @ 30 l/m  
1.4 cmH<sub>2</sub>O @ 60 l/m  
2.76 cmH<sub>2</sub>O @ 90 l/m
- Deadspace: 33 ml
- BFE (Bacterial Filtration Efficiency): 99.9999%
- VFE (Viral Filtration Efficiency): 99.999 %

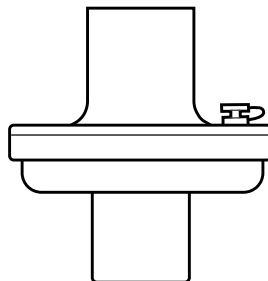
**Part No:** 007963



### Low pressure oxygen adapter

**Function:** Oxygen tube adapter with connector for the Vivo 45 LS.

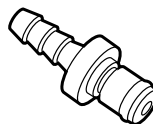
**Part No:** 005032



### Circuit: Single limb for Mouthpiece ventilation (MPV)

**Function:** Deliver air to the patient

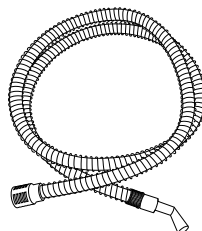**Part No:** 006093



### Mouthpiece

**Function:** Patient interface for Mouthpiece ventilation (MPV)

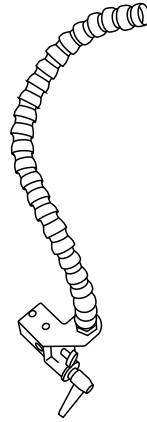
**Part No:** 006094



### MPV arm

**Function:** Hold an MPV circuit so Mouthpiece can be mounted close to the patient

**Part No:** 006095




## 9.2 Power Accessories

### Car Adapter Cable

**Function:** 12–24 VDC car adapter cable.

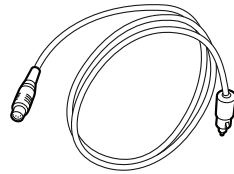
**Part No:** 007653



### Power Supply

**Function:** Deliver power to the ventilator

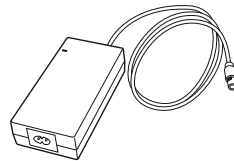**Part No:** 006396



### Power cord

**Function:** Deliver power to the AC power supply

**Part No:**

GB: 003521

CN: 005304

EU: 003520

JP: 004834

US: 003522



### XPAC - External battery with charger

**Function:** Extends usage time of supported Breas products.

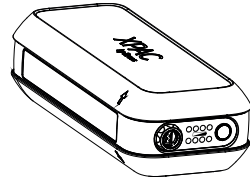
**Part No Cable for connection to device:** 007671

**Part No Charger with cable:**

*Single: Charger with one battery*

*Dual: Charger with two batteries*


Single: 007995, Dual: 007999



### Click-in battery

**Function:** Power source for transportation

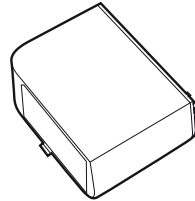**Part No:** 006265




### Click-in Battery Charger

**Function:** External charger for click-in batteries, available with bank for 2 or 4 batteries)

**Part no:**

07730 (2 batteries charger)

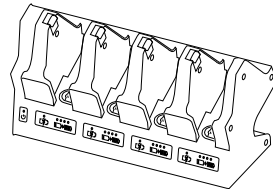
007731 (4 batteries charger)



### Cable, external DC

**Function:** External DC cable.

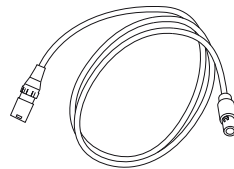
**Part No:**006709



### Cable, external DC to Ventilator Adapter

**Function:** Connect the ventilator to external DC

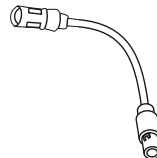
**Part No:** 006710



### **Cable, Y-adapter, Mains AC and External DC to Ventilator**

**Function:** Connects the ventilator to both mains and external DC at the same time. If the Mains power source is available, it will have precedence over the DC power source.

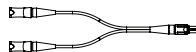**Part No:** 006711




## **9.3 Monitoring Accessories**

### **Breas PC software**

**Function:** Data monitoring software

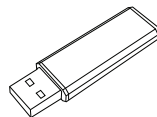**Part No:** 006718



### **USB cable**

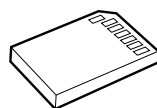
**Function:** Data cable: PC to Vivo 45 LS (USB to USB)

**Part No:** 005757

### **Memory card**

**Function:** Storage and transfer of settings, patient data and usage data


**Part No:** 006705



### **Memory card reader/writer**

**Function:** Read/write memory card

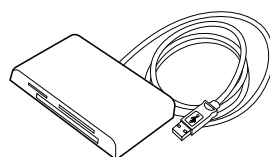
**Part No:** 002185



### **Remote alarm with cable**

**Function:** Monitor Vivo 45 LS alarms remotely

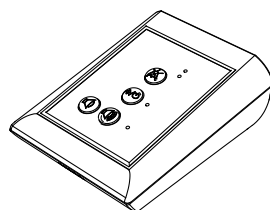**Part No:** 10 m: 006348, 25 m: 006349






### Remote alarm cable

**Function:**

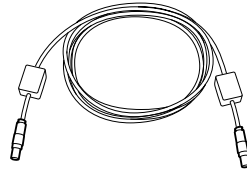**Part No:** 10 m: 006359, 25 m: 006360, 50 m: 006361



### Nurse call cable

**Function:** Connect the ventilator to a hospital nurse call system

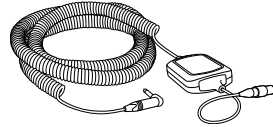
**Part No:**

NO: 006365

NC: 006364

10 k $\Omega$ , NO: 006363


10 k $\Omega$ , NC: 006362



### Remote start/stop

**Function:** Start and stop the ventilator remotely. Also, pause audio remotely.

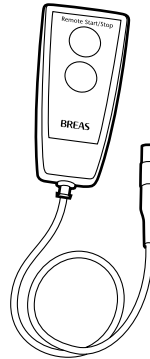**Part No:** 006649



### FiO<sub>2</sub> sensor

**Function:** Measure FiO<sub>2</sub> to the patient.

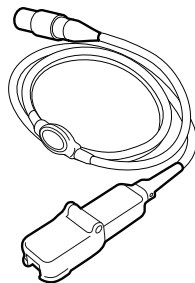
**Part No:** 006172



### **SpO<sub>2</sub> module**

**Function:** Connection interface

**Part No:** 006369




### **SpO<sub>2</sub> sensor**

**Function:** Finger Clip SpO<sub>2</sub> sensor

**Part No:**

Adult: 006589

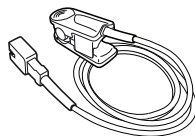
Paediatric: 006590



### **SpO<sub>2</sub> sensor**

**Function:** Multisite SpO<sub>2</sub> sensor

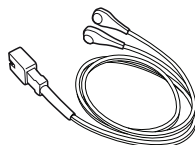
**Part No:** 006591



### **EtCO<sub>2</sub> sensor**

**Function:** Measure CO<sub>2</sub> in the airflow

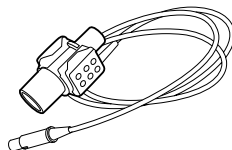
**Part No:** 006346



### **Airway adapter**

**Function:** Connects the EtCO<sub>2</sub> sensor to the patient circuit

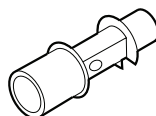
**Part No:** 005263 (25 pcs)



### PtcCO<sub>2</sub> Cable, Sentec

**Function:** Connects the ventilator to a Sentec PtcCO<sub>2</sub> monitor.

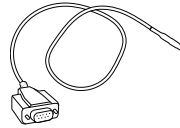**Part No:** 006179



### PtcCO<sub>2</sub> Cable, Radiometer

**Function:** Connects the ventilator to a Radiometer TCM5 PtcCO<sub>2</sub> monitor.

**Part No:** 008392



### Effort belt communication box

**Function:** Connects the ventilator to one or two effort belts.

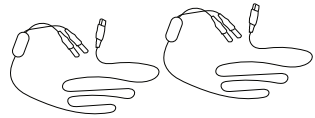
**Part No:** 006182



### Effort belt wireset

**Function:** Connects an effort belt to the communication box

**Part No:** 007083



### Effort belt

**Function:** Measures respiratory effort

**Part No:**

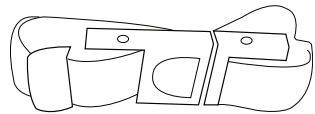
Adult:

24"-74" (107-188 cm): 007085

45"-123" (114-312 cm): 007091

Paediatric:

16"-42" (41-107 cm): 007084



## 9.4 Ventilator Filters and Detachable Parts

### Patient air inlet filter, fine, white, disposable

**Function:** Fine inlet air filtration.

**Material:** AS 100

**NaCl Penetration:**  $(0.65 \mu\text{m NaCl @ } 95 \text{ l/min}) = <7.35\%$

**Part No:** 007103 (5pcs)



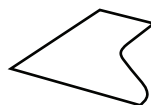
### Patient air inlet filter, coarse, grey, washable

**Function:** Coarse inlet air filtration

**Material:** Bulpren S 28133

**Filter diameter:** 1080-1580 Microns


**Part No:** 007104 (5pcs)



### Cooling air filter

**Function:** Device air inlet filtration, 5 pieces

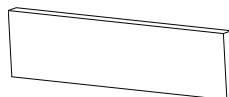
**Part No:** 007105



### Air bypass unit

**Function:** Direct the air flow within the ventilator

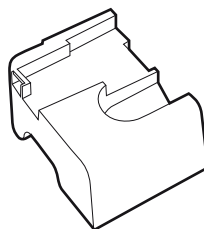
**Part No:** 007064



### Click-in water chamber

**Function:** Humidify the patient air

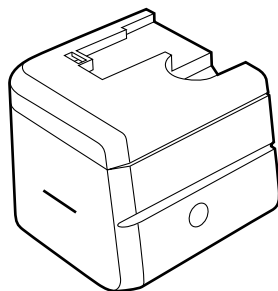**Part No:** 006490




### Side panels

**Function:** Protect the internal ventilator components.

**Part No:**

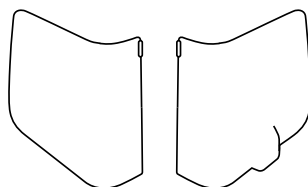Grey: 007065, Blue: 007066, Light blue: 007518




## 9.5

### Other Accessories

#### Trolley

**Function:** Mobile use, transportation

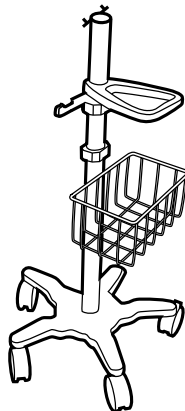
**Part No:** 007384



### Patient circuit arm

**Function:** To support a patient circuit.

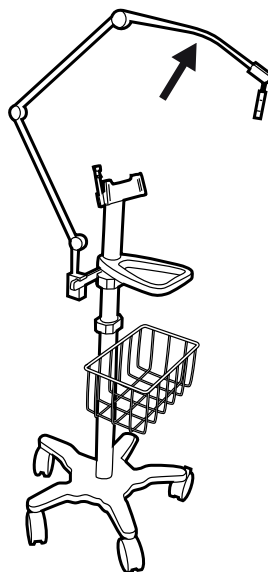
**Part No:** 007917



### Universal rail clamp

**Function:** Attach a humidifier to a trolley. This accessory is part of the trolley system.


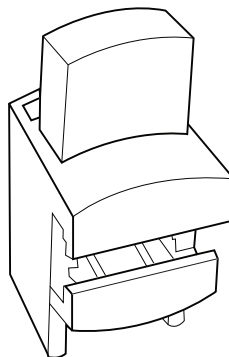
**Part No:** 007858



### E-cylinder holder

**Function:** Attach an E-cylinder to a trolley. This accessory is part of the trolley system.

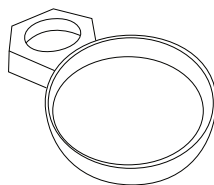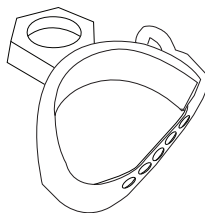
**Part No:** 005128



### IV-pole

**Function:** Pole with hooks to hang IV fluid bags.

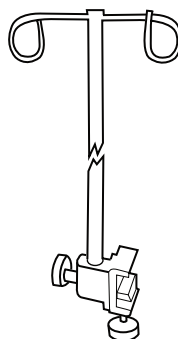
**Part No:** 007859



### Mounting bracket

**Function:** Mount the ventilator to a stand / trolley / rail system.

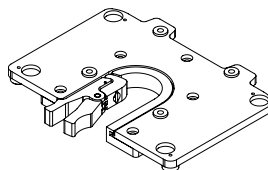**Part No:** 006761



**Protective cover**

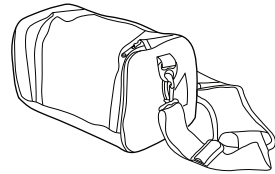
**Function:** Shock protection

**Part No:** 006067

**Lightweight Mobility Bag**

**Function:** Mobile use

**Part No:** 007555

**Carry bag**

**Function:** Storage for transportation

**Part No:** 006014





10 Patient Settings

This section can be copied and used for noting the patient’s settings.

Patient Settings - Vivo 45 LS

Patient .....

Date .....

Clinic .....

Set by .....

Ventilation mode:.....

|                  |                         |
|------------------|-------------------------|
| Patient Circuit  |                         |
| Insp. Pressure   | Inspiratory Trigger     |
| PEEP             | Expiratory Trigger      |
| Breath Rate      | Min Inspiratory Time    |
| Inspiratory Time | Max Inspiratory Time    |
| Backup Rate      | Backup Inspiratory Time |
| Target Volume    | Min Pressure            |
| Max Pressure     | CPAP                    |
| SIMV Rate        | SIMV Support Pressure   |
| Auto-EPAP        | EPAP step               |
| EPAP Min         | EPAP Max                |
| PS               | Pressure Limit )        |
| PS Min           | PS Max                  |
| Relax Time       |                         |

Notes

.....

## FAA Compliance

To whom it may concern:

In line with the FAA Advisory Circular AC 91.21-1D October 27, 2017, this kind of respiratory assistive device may be used onboard an aircraft, without further testing by the carrier, provided they have been tested for Electromagnetic Compatibility (EMC) in accordance with the current version of RTCA/DO-160, Section 21, Category M.

Breas Medical has successfully completed testing for the ventilator System. The ventilator System complies with RTCA/DO-160, Section 21, Category M and can be considered FAA compliant.

Some airlines may require advance notification before travel, and devices may need to be operated by battery. Breas Medical recommends that customers check with their airline.

### FAA Compliance (English text)

To whom it may concern:

In line with the FAA Advisory Circular AC 91.21-1D October 27, 2017, this kind of respiratory assistive device may be used onboard an aircraft, without further testing by the carrier, provided they have been tested for Electromagnetic Compatibility (EMC) in accordance with the current version of RTCA/DO-160, Section 21, Category M.

Breas Medical has successfully completed testing for the ventilator System. The ventilator System complies with RTCA/DO-160, Section 21, Category M and can be considered FAA compliant.

Some airlines may require advance notification before travel, and devices may need to be operated by battery. Breas Medical recommends that customers check with their airline.

# Index

## A

|   |     |
|---|-----|
| AC power  |     |
| connect .....   | 34  |
| Accessories   |     |
| usage .....   | 74  |
| Active exhalation valve patient circuit   |     |
| connect .....   | 35  |
| Adjust  |     |
| patient settings .....  | 37  |
| Adult mode .....  | 54  |
| Air filters   |     |
| clean and replace .....   | 126 |
| Air inlet, position .....   | 26  |
| Air pathway   |     |
| disinfection .....  | 125 |
| Air Temp. Sensor Fail alarm .....   | 118 |
| Alarm   |     |
| Air Temp. Sensor Fail .....   | 118 |
| Alarm Battery Error .....   | 115 |
| Ambient Pressure Compensation   |     |
| Lost .....  | 114 |
| Ambient Temperature Compensation  |     |
| Lost .....  | 114 |
| Apnea .....   | 101 |
| Check CO <sub>2</sub> Adapter .....   | 113 |
| Clock failure .....   | 119 |
| CO <sub>2</sub> Disconnected (CO <sub>2</sub> Sensor Failure/Disconnection Alarm) .....   | 112 |
| CO <sub>2</sub> Sensor Error .....  | 113 |
| CO <sub>2</sub> Sensor Error (CO <sub>2</sub> Accuracy Error Alarm) .....                 | 112 |
| Cooling Fan Error .....   | 119 |
| Crit. Low Last Power Source .....   | 110 |
| Database Integrity Fail .....   | 119 |
| Disconnection .....   | 101 |
| Exhalation Valve Control  |     |
| Error .....   | 111 |
| FiO <sub>2</sub> Disconnected (FiO <sub>2</sub> Sensor Failure/Disconnection Alarm) ..... | 113 |
| Heated Circuit Fault .....  | 117 |
| Heated Circuit Temp. .....  | 116 |
| High breath rate .....  | 100 |
| High EtCO <sub>2</sub> .....  | 106 |
| High FiO <sub>2</sub> .....   | 104 |
| High Humidifier Temp. .....   | 116 |
| High InspCO <sub>2</sub> .....  | 107 |

|   |     |
|---|-----|
| High MVe .....  | 99  |
| High Mvi .....  | 97  |
| High Patient Air Temp (High Patient Air Temperature) .....                                | 109 |
| High PEEP .....   | 95  |
| High pressure .....   | 93  |
| High pulse rate .....   | 107 |
| High SpO <sub>2</sub> .....   | 105 |
| High Vte .....  | 98  |
| High Vti .....  | 96  |
| Humidifier Fault .....  | 117 |
| Humidifier/bypass loose .....   | 120 |
| Humidity Compensation Lost .....  | 114 |
| Internal Error .....  | 118 |
| Internal Function Failure .....   | 118 |
| Internal temp high .....  | 119 |
| Internal/Click-In Battery Hot .....   | 116 |
| LED Failure .....   | 115 |
| Lost AC Power (Mains) .....   | 110 |
| Low Alarm Battery .....   | 115 |
| Low Battery .....   | 110 |
| Low breath rate .....   | 100 |
| Low EtCO <sub>2</sub> .....   | 106 |
| Low FiO <sub>2</sub> .....  | 104 |
| Low MVe .....   | 99  |
| Low MVi .....   | 97  |
| Low Patient Air Temp. (Low Patient Air Temperature Alarm) .....                           | 109 |
| Low PEEP .....  | 95  |
| Low pressure .....  | 94  |
| Low pulse rate .....  | 108 |
| Low SpO <sub>2</sub> .....  | 105 |
| Low Vte .....   | 98  |
| Low Vti .....   | 96  |
| Obstruction .....   | 103 |
| Poor SpO <sub>2</sub> Signal .....  | 112 |
| Power Failure Alarm .....   | 109 |
| PtcCO <sub>2</sub> Alarm .....  | 108 |
| Rebreathing .....   | 103 |
| SpO <sub>2</sub> Disconnected (SpO <sub>2</sub> Sensor Failure/Disconnection Alarm) ..... | 111 |
| SpO <sub>2</sub> Signal Lost .....  | 111 |
| Alarm Battery Error alarm .....   | 115 |
| Alarms  |     |
| technical .....   | 109 |
| Ambient Pressure Compensation Lost alarm .....  | 114 |

|   |     |
|---|-----|
| Ambient Temperature Compensation Lost alarm .....           | 114 |
| Apnea alarm .....   | 101 |
| Artificial nose .....                                       | 19  |
| Assisted Pressure Controlled Ventilation                    |     |
| definition .....  | 56  |
| Assisted Pressure Controlled Ventilation with Target Volume |     |
| definition .....  | 56  |
| Audible range   |     |
| operator's position .....                                   | 92  |
| Audio pause and reset .....                                 | 92  |

## B

|                                |    |
|--------------------------------|----|
| Backup insp. time .....        | 63 |
| Bacterial filter .....         | 18 |
| Battery                        |    |
| charge .....                   | 69 |
| click-in .....                 | 69 |
| connect click-in battery ..... | 69 |
| Breath mode .....              | 54 |
| Breath rate .....              | 62 |

## C

|                               |     |
|-------------------------------|-----|
| Cables                        |     |
| inspect .....                 | 36  |
| Calibration                   |     |
| CO2 .....                     | 77  |
| Caution                       |     |
| Icon .....                    | 10  |
| Change patients .....         | 127 |
| Charge                        |     |
| batteries .....               | 72  |
| Check                         |     |
| before first use .....        | 32  |
| Check CO2 Adapter alarm ..... | 113 |
| Circuit                       |     |
| connect .....                 | 35  |
| Cleaning .....                | 124 |
| CO2 sensor .....              | 77  |
| main unit .....               | 124 |
| patient .....                 | 125 |
| Cleaning and maintenance      |     |
| safety information .....      | 20  |
| Click-in battery .....        | 71  |
| charger .....                 | 72  |

|   |     |
|---|-----|
| connect .....   | 72  |
| Click-in humidifier   |     |
| using .....   | 83  |
| Clock failure alarm .....   | 119 |
| CO2 Accuracy Error Alarm (CO2 Sensor Error) alarm .....               | 112 |
| CO2 sensor  |     |
| clean .....   | 77  |
| connect .....   | 77  |
| safety information .....  | 77  |
| specifications .....  | 135 |
| usage .....   | 77  |
| CO2 Sensor Error alarm .....  | 113 |
| CO2 Sensor Failure/Disconnection Alarm (CO2 Disconnected) alarm ..... | 112 |
| CO2 unit .....  | 49  |
| CO2 zeroing .....   | 78  |
| Compliance data   |     |
| read .....  | 49  |
| Connect   |     |
| click-in battery .....  | 71  |
| CO2 sensor .....  | 77  |
| nurse call .....  | 75  |
| to mains .....  | 34  |
| Continuous Positive Airway Pressure                                   |     |
| definition .....  | 58  |
| Cooling air filter  |     |
| replace .....   | 126 |
| Cooling Fan alarm .....   | 119 |
| CPAP .....  | 67  |
| definition .....  | 58  |
| Crit. Low Last Power Source alarm .....                               | 110 |
| Curves screen .....   | 49  |

## D

|                            |     |
|----------------------------|-----|
| Date                       |     |
| set .....                  | 49  |
| Date format .....          | 49  |
| DC power                   |     |
| click-in battery LED ..... | 25  |
| external DC .....          | 25  |
| external DC LED .....      | 25  |
| internal battery LED ..... | 25  |
| Delivery settings .....    | 142 |
| Device information .....   | 49  |
| Device mode .....          | 54  |
| Device settings .....      | 49  |

|                          |     |
|--------------------------|-----|
| Dimensions .....         | 135 |
| Disconnection alarm..... | 101 |
| Disinfection             |     |
| air pathway .....        | 125 |
| Display                  |     |
| home mode overview ..... | 46  |
| light.....               | 49  |
| navigate.....            | 41  |
| overview.....            | 45  |
| symbols.....             | 42  |

## E

|                                |     |
|--------------------------------|-----|
| Electrical safety              |     |
| user precautions.....          | 15  |
| Enter                          |     |
| standby mode .....             | 40  |
| Environmental conditions       |     |
| safety information.....        | 15  |
| specifications.....            | 130 |
| Equipment designation and      |     |
| safety label.....              | 29  |
| Exhalation Valve Control Error |     |
| alarm.....                     | 111 |
| Exp. Trigger.....              | 65  |
| External DC .....              | 72  |

## F

|   |     |
|---|-----|
| Filter  |     |
| bacterial .....   | 18  |
| disposable .....  | 23  |
| safety information.....   | 18  |
| Filtering/smoothing techniques  |     |
| specifications.....   | 130 |
| Filters   |     |
| clean and replace.....  | 126 |
| FiO <sub>2</sub> Sensor Failure/Disconnection Alarm (FiO <sub>2</sub> Disconnected) alarm ..... | 113 |
| First use check .....   | 32  |
| Font panel  |     |
| main unit .....   | 25  |
| Function  |     |
| navigation buttons .....  | 25  |

## G

|                                |    |
|--------------------------------|----|
| General user precautions ..... | 12 |
|--------------------------------|----|

## H

|  |        |
|--|--------|
| HCH .....                              | 19     |
| Heat and moisture exchanger .....      | 19     |
| Heated circuit                         |        |
| settings .....                         | 49     |
| Heated circuit (on/off. setting).....  | 67     |
| Heated circuit (temp. setting).....    | 67     |
| Heated Circuit Fault alarm .....       | 117    |
| Heated Circuit Temp alarm.....         | 116    |
| Heated patient circuit                 |        |
| connect.....                           | 35     |
| using.....                             | 88     |
| High breath rate alarm.....            | 100    |
| High EtCO <sub>2</sub> alarm.....      | 106    |
| High FiO <sub>2</sub> alarm .....      | 104    |
| High Humidifier Temp. alarm.....       | 116    |
| High InspCO <sub>2</sub> alarm.....    | 107    |
| High MVe alarm .....                   | 99     |
| High Mvi alarm .....                   | 97     |
| High Patient Air Temperature           |        |
| Alarm (High Patient Air                |        |
| Temp) alarm .....                      | 109    |
| High PEEP alarm.....                   | 95     |
| High pressure alarm .....              | 93     |
| High pulse rate alarm .....            | 107    |
| High SpO <sub>2</sub> alarm.....       | 105    |
| High Vte alarm .....                   | 98     |
| High Vti alarm .....                   | 96     |
| HME .....                              | 19     |
| Home Adjust.....                       | 53     |
| Home mode.....                         | 49, 53 |
| activate .....                         | 53     |
| overview.....                          | 44     |
| Humidification                         |        |
| safety information.....                | 19     |
| Humidifier                             |        |
| detaching .....                        | 87     |
| filling .....                          | 84     |
| installing .....                       | 85     |
| Humidifier (on/off setting).....       | 67     |
| Humidifier (value setting) .....       | 67     |
| Humidifier Fault alarm .....           | 117    |
| Humidifier/bypass loose alarm .....    | 120    |
| Humidity Compensation Lost             |        |
| alarm .....                            | 114    |
| Hygroscopic condenser humidifier... 19 |        |
| Hypoventilation.....                   | 100    |

## I

|                                       |     |
|---------------------------------------|-----|
| Icon .....                            | 70  |
| Caution .....                         | 10  |
| Note .....                            | 10  |
| Reference .....                       | 10  |
| Warning .....                         | 10  |
| Information button .....              | 25  |
| Insp. Pressure .....                  | 62  |
| Insp. Time .....                      | 63  |
| Insp. trigger .....                   | 64  |
| Inspect                               |     |
| before use .....                      | 36  |
| cables .....                          | 36  |
| placement .....                       | 36  |
| Inspiratory time                      |     |
| sigh .....                            | 63  |
| Intended use .....                    | 9   |
| Internal battery .....                | 70  |
| Internal Error alarm .....            | 118 |
| Internal Function Failure alarm ..... | 118 |
| Internal temp high alarm .....        | 119 |
| Internal/Click-In Battery Hot         |     |
| alarm .....                           | 116 |

## L

|                                   |     |
|-----------------------------------|-----|
| Label .....                       | 29  |
| Language                          |     |
| selecting for clinical mode ..... | 49  |
| selecting for home mode .....     | 49  |
| Leakage circuit (single limb)     |     |
| pneumatic diagram .....           | 130 |
| LED                               |     |
| click-in battery .....            | 25  |
| external DC .....                 | 25  |
| front panel .....                 | 25  |
| internal battery .....            | 25  |
| LED Failure alarm .....           | 115 |
| Light intensity .....             | 49  |
| Lost AC Power (Mains) alarm ..... | 110 |
| Low Alarm Battery alarm .....     | 115 |
| Low Battery alarm .....           | 110 |
| Low breath rate alarm .....       | 100 |
| Low EtCO <sub>2</sub> alarm ..... | 106 |
| Low FiO <sub>2</sub> alarm .....  | 104 |
| Low MVe alarm .....               | 99  |
| Low MV <sub>i</sub> alarm .....   | 97  |

## Low Patient Air Temperature

|                                  |     |
|----------------------------------|-----|
| Alarm (Low Patient Air Temp)     |     |
| alarm .....                      | 109 |
| Low PEEP alarm .....             | 95  |
| Low pressure alarm .....         | 94  |
| Low pulse rate .....             | 108 |
| Low SpO <sub>2</sub> alarm ..... | 105 |
| Low Vte alarm .....              | 98  |
| Low Vti alarm .....              | 96  |

## M

|  |     |
|--|-----|
| Main components .....                      | 23  |
| Main unit                                  |     |
| clean .....                                | 124 |
| Mains                                      |     |
| connect .....                              | 34  |
| Maintenance                                |     |
| filters .....                              | 126 |
| service information .....                  | 127 |
| Maximum inspiratory time .....             | 65  |
| Memory card slot                           |     |
| position .....                             | 26  |
| Menu                                       |     |
| navigate .....                             | 41  |
| symbols .....                              | 42  |
| Minimum inspiratory time .....             | 65  |
| Modes                                      |     |
| definitions .....                          | 53  |
| Monitor section .....                      | 49  |
| Monitored value                            |     |
| % in TgV .....                             | 52  |
| % Spont .....                              | 52  |
| EtCO <sub>2</sub> .....                    | 53  |
| FiO <sub>2</sub> .....                     | 52  |
| I:E .....                                  | 53  |
| Insp. time .....                           | 53  |
| InspCO <sub>2</sub> .....                  | 53  |
| Leakage .....                              | 51  |
| MVe .....                                  | 52  |
| MV <sub>i</sub> .....                      | 52  |
| PEEP .....                                 | 51  |
| Pmean .....                                | 51  |
| Ppeak .....                                | 51  |
| PtcCO <sub>2</sub> .....                   | 53  |
| Pulse Rate .....                           | 53  |
| Rise time .....                            | 53  |
| SpO <sub>2</sub> (oxygen saturation) ..... | 53  |
| Spont Rate .....                           | 52  |
| Total Rate .....                           | 52  |

|                            |     |  |     |
|----------------------------|-----|--|-----|
| Vte.....                   | 52  | safety information.....                  | 16  |
| Vti .....                  | 52  | Patient mode .....                       | 54  |
| Monitored values .....     | 50  | Patient operating time .....             | 49  |
| Mouthpiece                 |     | Patient profiles .....                   | 49  |
| pressure mode.....         | 57  | Patient settings                         |     |
| volume mode.....           | 58  | adjust .....                             | 37  |
| MPV                        |     | Patients                                 |     |
| pressure mode.....         | 57  | change .....                             | 127 |
| volume mode.....           | 58  | PCV                                      |     |
| <b>N</b>                   |     | definition .....                         | 56  |
| Navigate                   |     | mouthpiece.....                          | 57  |
| menu .....                 | 41  | PCV-SIMV                                 |     |
| Note                       |     | definition .....                         | 56  |
| Icon.....                  | 10  | PCV(A)                                   |     |
| Nurse call                 |     | definition .....                         | 56  |
| connect.....               | 75  | PCV(A+TgV)                               |     |
| electrical safety .....    | 15  | definition .....                         | 56  |
| <b>O</b>                   |     | PCV(TgV)                                 |     |
| Obstruction alarm .....    | 103 | definition .....                         | 56  |
| On/Off                     |     | PEEP.....                                | 62  |
| switch.....                | 40  | Perform pre-use test .....               | 37  |
| Operating conditions       |     | Placement .....                          | 32  |
| specifications .....       | 130 | inspect .....                            | 36  |
| Operating mode .....       | 54  | Pneumatic diagram                        |     |
| Operator's position        |     | leakage circuit (single limb) .....      | 130 |
| audible range .....        | 92  | single limb circuit .....                | 130 |
| Others section .....       | 49  | Poor SpO <sub>2</sub> Signal alarm ..... | 112 |
| Oxygen                     |     | Position                                 |     |
| safety information.....    | 21  | operator, audible range .....            | 92  |
| usage.....                 | 21  | Power Failure Alarm alarm .....          | 109 |
| Oxygen inlet               |     | Power source priority .....              | 68  |
| position.....              | 29  | Power supply specification .....         | 130 |
| specification .....        | 130 | Pre-use test.....                        | 49  |
| <b>P</b>                   |     | Pre-Use Test                             |     |
| Paediatric mode.....       | 54  | Failure .....                            | 38  |
| Part numbers               |     | Pressure controlled ventilation          |     |
| main components.....       | 23  | mouthpiece.....                          | 57  |
| Pathway                    |     | Pressure Controlled Ventilation          |     |
| internal disinfection..... | 125 | definition .....                         | 56  |
| Patient air filters        |     | Pressure Controlled Ventilation with     |     |
| clean and replace.....     | 126 | Target Volume                            |     |
| Patient air outlet         |     | definition .....                         | 56  |
| position.....              | 26  | Pressure Support Ventilation             |     |
| Patient circuit            |     | definition .....                         | 55  |
| connect.....               | 35  | Pressure Support Ventilation with Tar-   |     |
|                            |     | get Volume                               |     |
|                            |     | definition .....                         | 55  |
|                            |     | Pressure unit .....                      | 49  |
|                            |     | Priority                                 |     |

|                               |     |   |     |
|-------------------------------|-----|---|-----|
| power source.....             | 68  | Setup .....                             | 47  |
| Product label.....            | 29  | Serial number .....                     | 49  |
| Profile                       |     | Service.....                            | 127 |
| select.....                   | 47  | Settings                                |     |
| Profiles                      |     | at delivery .....                       | 142 |
| delivery settings .....       | 142 | specification .....                     | 130 |
| Protective cover.....         | 81  | Side panel .....                        | 26  |
| PSV                           |     | Sigh                                    |     |
| definition .....              | 55  | settings .....                          | 63  |
| PSV(TgV)                      |     | SIMV                                    |     |
| definition .....              | 55  | in PCV mode.....                        | 56  |
| PtcCO <sub>2</sub> Alarm..... | 108 | in VCV mode .....                       | 57  |
| <b>R</b>                      |     | SIMV rate .....                         | 62  |
| Rebreathing alarm .....       | 103 | SIMV support pressure.....              | 65  |
| Reference                     |     | Single limb circuit                     |     |
| Icon.....                     | 10  | pneumatic diagram.....                  | 130 |
| Regular maintenance.....      | 127 | Sound level                             |     |
| Remote alarm                  |     | specification .....                     | 130 |
| using.....                    | 76  | Specifications                          |     |
| Remote Start/Stop.....        | 81  | CO <sub>2</sub> sensor.....             | 130 |
| Repair .....                  | 127 | dimension .....                         | 130 |
| Reset alarm.....              | 92  | environmental conditions.....           | 130 |
| Rise time .....               | 64  | Filtering/smoothing                     |     |
| <b>S</b>                      |     | techniques.....                         | 130 |
| Safety information            |     | miscellaneous.....                      | 130 |
| cleaning and maintenance..... | 20  | monitored values.....                   | 130 |
| CO <sub>2</sub> sensor.....   | 77  | operating conditions .....              | 130 |
| environmental conditions..... | 15  | oxygen inlet .....                      | 130 |
| filters .....                 | 15  | power supply .....                      | 130 |
| humidification.....           | 19  | settings .....                          | 130 |
| oxygen.....                   | 21  | sound level .....                       | 130 |
| patient circuit.....          | 16  | SpO <sub>2</sub>                        |     |
| Safety label .....            | 29  | usage.....                              | 81  |
| Screen                        |     | SpO <sub>2</sub> Sensor Failure/Discon- |     |
| alarm .....                   | 48  | nection Alarm (SpO <sub>2</sub> Discon- |     |
| alarm/event history.....      | 48  | nected) alarm .....                     | 111 |
| curves .....                  | 49  | Standby mode                            |     |
| home mode main screen         |     | enter .....                             | 54  |
| overview .....                | 46  | Start                                   |     |
| Monitor.....                  | 49  | switch on and enter standby             |     |
| Others .....                  | 49  | mode .....                              | 40  |
| Trends .....                  | 49  | treatment .....                         | 40  |
| Section                       |     | Start-Stop                              |     |
| Alarm.....                    | 48  | using remote .....                      | 81  |
| Monitor.....                  | 49  | Start-up time                           |     |
| Others .....                  | 49  | specification .....                     | 130 |
|                               |     | Stop                                    |     |
|                               |     | treatment .....                         | 40  |
|                               |     | storage                                 |     |



|                                  |     |
|----------------------------------|-----|
| battery .....                    | 72  |
| Storage .....                    | 128 |
| Sup. pressure .....              | 65  |
| Switch off .....                 | 40  |
| Switch on/off .....              | 40  |
| Symbols                          |     |
| equipment designation and safety |     |
| label .....                      | 29  |
| menu .....                       | 42  |

## T

|                              |     |
|------------------------------|-----|
| Target volume .....          | 66  |
| in PCV + A, definition ..... | 56  |
| in PCV, definition .....     | 56  |
| in PSV, definition .....     | 55  |
| Technical alarms .....       | 109 |
| Time                         |     |
| set .....                    | 49  |
| Time format .....            | 49  |
| Treatment                    |     |
| select preset profile .....  | 47  |
| start .....                  | 40  |
| stop .....                   | 40  |
| Trends screen .....          | 49  |

## U

|                                |    |
|--------------------------------|----|
| Up/Down buttons .....          | 25 |
| User interface                 |    |
| navigate .....                 | 41 |
| symbols .....                  | 42 |
| User precautions               |    |
| electrical safety .....        | 15 |
| environmental conditions ..... | 15 |
| general .....                  | 12 |
| patient circuit .....          | 16 |

## V

|                               |    |
|-------------------------------|----|
| Values                        |    |
| monitored .....               | 50 |
| VCV                           |    |
| definition .....              | 57 |
| mouthpiece .....              | 58 |
| VCV-SIMV                      |    |
| definition .....              | 57 |
| VCV(A)                        |    |
| definition .....              | 57 |
| Ventilation mode .....        | 54 |
| Volume controlled ventilation |    |

|                  |    |
|------------------|----|
| mouthpiece ..... | 58 |
|------------------|----|

## W

|                        |     |
|------------------------|-----|
| Warning                |     |
| Icon .....             | 10  |
| Water chamber          |     |
| detaching .....        | 87  |
| filling .....          | 84  |
| installing .....       | 85  |
| open lid .....         | 88  |
| Water trap .....       | 19  |
| Weight .....           | 130 |
| Wheelchair             |     |
| protective cover ..... | 81  |